



Dec 07, 2021

Nextera XT protocol for MiSeq HIV^{PR-RT} sequencing

Paula Aulicino¹¹Hospital de Pediatría "Prof. Dr. Juan P. Garrahan"

1

dx.doi.org/10.17504/protocols.io.b2pjgdkn**HIV GENOTYPIC STUDY BY NGS (ILLUMINA)** Paula Aulicino

The protocol allows amplification and NGS sequencing of a PR-RT HIV-1 sequence using tagmentation strategy and Illumina MiSeq equipment, followed by analysis of drug-resistance associated mutations by HyDRA

DOI

dx.doi.org/10.17504/protocols.io.b2pjgdkn

Paula Aulicino 2021. Nextera XT protocol for MiSeq HIVPR-RT sequencing.

protocols.io<https://dx.doi.org/10.17504/protocols.io.b2pjgdkn> protocol ,

Dec 07, 2021

Dec 07, 2021


55755

- 1 Prepare the sample in BSL2 laboratory. This protocol can be used to amplify and sequence HIV-1 from plasma or culture supernatants, as well as from cell-associated HIV-1. In the latter, skip the step of reverse transcription and go straight to PCR amplification.
- 2 Perform reverse transcription of HIV-1 and PCR amplification of PR-RT 1168bp fragment using primers and conditions published in Aulicino et al, Antivir Ther 2010;15(4):641-50. PCR amplifications are performed in 25-cycle reactions to avoid introducing artificial mutations.

Aulicino PC, Rocco CA, Mecikovsky D, Bologna R, Mangano A, Sen L (2010). HIV type-1 genotypic resistance profiles in vertically infected patients from Argentina reveal an association between K103N+L100I and L74V mutations.. Antiviral therapy.
<https://doi.org/10.3851/IMP1571>

- 3 Check second PCR products with agarose gel electrophoresis (1%). A single 1168bp band should be visible.
- 4 Perform purification of the PCR product with magnetic beads (Agencourt AMPure XP kit, Beckman Coulter, USA)
- 5 Perform tagmentation and indexation with Nextera XT DNA Library Preparation Kit (Illumina, San Diego, CA, USA), following manufacturer's instructions.
- 6 Check electrophoresis tagmentation patterns with QIAxel or another similar equipment. If none is available, a polyacrylamide gel can be used.
- 7 Quantify DNA using the fluorometric Qubit dsDNA high-sensitivity (HS) assay (Qubit 1X dsDNA HS Assay Kit, ThermoFisher Scientific, USA)
- 8 Calculate PCR samples concentration and volume to normalize all samples to 2nM, using PCR water.
- 9 Pool all samples using 5ul of each 2nM dilution of each sample.
- 10 Prepare PhiX at 12.5 pM. Dilute the library to 8 picomolar in HT1 buffer, and add 20% of PhiX. The purpose of PhiX is to add diversity to the samples.

11 Run in Illumina MiSeq

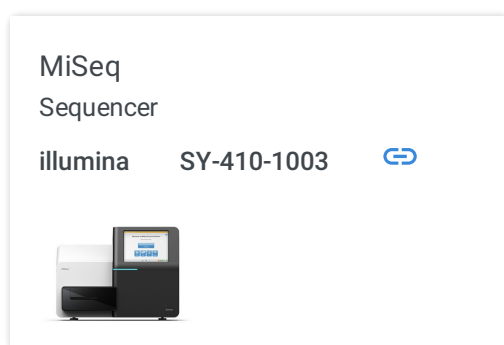

Controls and Error Rate calculation

- 12 pNL4-3 HIV-1 was used as control. Using BWA, align fastq reads belonging to the control sample to the reference pNL4-3 HIV-1 genome (GenBank Accession Number AF324493). Calculate the average number of mismatches/ total bases. Error Rate should not exceed 0,5%.

Sequence analysis using HyDRA

- 13 HyDRA (<https://hydra.canada.ca/pages/home?lang=en-CA>) is an annotated reference-based bioinformatics pipeline scripted in Perl, which analyses next generation sequencing (NGS) data for genotyping HIV-1 drug resistance (HIVDR) mutations. It utilizes an annotated HXB2 sequence (GenBank Accession number: K03455) for reference mapping by Bowtie2 2 , and stringent data quality assurance and variant calling criteria to identify HIVDR associated mutations based on the Stanford HIV Drug Resistance Database (<http://hivdb.stanford.edu/>) and 2009 WHO list for Surveillance of Transmitted HIVDR. All HIVDR mutations found in the pol genes; protease (PR), reverse transcriptase (RT), and integrase (IN), are reported according to classifications outlined in the Stanford Surveillance Drug Resistance Mutation list .

Taylor T, Lee ER, Nykoluk M, Enns E, Liang B, Capina R, Gauthier MK, Domselaar GV, Sandstrom P, Brooks J, Ji H (2019). A MiSeq-HyDRA platform for enhanced HIV drug resistance genotyping and surveillance.. Scientific reports.
<https://doi.org/10.1038/s41598-019-45328-3>

- 14 Submit sequences to HyDRA. The following Advanced options for Nextera seqs can be used to inform amino acids in 1% frequency or higher, with a minimum read depth of 1000 x:
Consensus percent: 20
Target coverage: -1
Length cutoff: 100
Score cutoff: 2

Error rate: include error rate calculated with the control (e.g. 0.005)
Minimum variant quality: 20
Minimum read depth: 1000
Minimum allele count: 7
Minimum amino acid frequency: 0.01