

May 08, 2024

WATER PRODUCTION FOR AWARE (Organic Contaminants)

DOI

dx.doi.org/10.17504/protocols.io.14egn67jql5d/v1



Celia Manaia¹

¹Universidade Católica Portuguesa



AWARE Project

Horizon Europe 101084245

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.14egn67jql5d/v1

Protocol Citation: Celia Manaia 2024. WATER PRODUCTION FOR AWARE (Organic Contaminants). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.14egn67jql5d/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: May 08, 2024

Last Modified: May 08, 2024

Protocol Integer ID: 99427

Keywords: water sampling, water processing, water analysis, waste water treatment, advanced tertiary treatment SOP

Funders Acknowledgement:

Horizon Europe

Grant ID: 101084245

Abstract

The protocol summarises the procedures used for analytical control. The protocol describes the Standard Operating Procedure (SOP) for the optimization of advanced tertiary treatment of water, based on a comprehensive quality and risk assessment.

Guidelines

Directive 2008/105/EC of the European Parliament and of the Council of 16 December 2008 on environmental quality standards in the field of water policy.


Materials

A	B	C	D	E	F	G	H
Parameter	V (mL) x R	S	Processing	Analytical method	Result	LOD / LOQ	Goal value
Organic contaminants	500 x 2	On ice	Solid-phase extraction	HPLC/MS-TOF	Presence/absence	-	Not present

Materials:

A	B	C	D
	Becker		cartridges (OASIS HLB 200 mg)
Pre-filtration	Buchner flask	SPE	vacuum pump
	PVDF 0.45 µm filters		graduated cylinder
	trap flask		becker
	vacuum pump		trap flask
			SPE vacuum manifold system

Safety warnings

 ---

Organic Contaminants:

- 1 The water production for AWARE main activities includes three stages – disinfection by ultraviolet C radiation (UVC), storage for 🕒 12:00:00 - 🕒 24:00:00 (according to water load and season) and ozonation. The water quality is monitored at these three stages, for the parameters indicated in Figure 1 below.

1d 12h

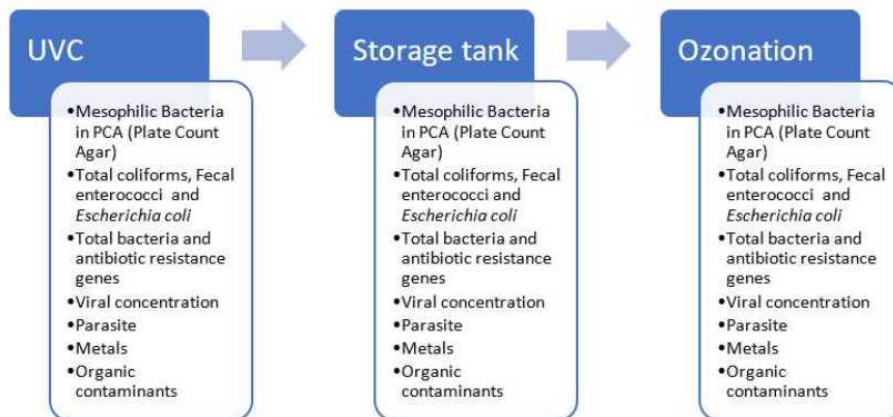


Figure 1. Treatment and storage of municipal treated wastewater used for integrated aquaponics and an indication of the comprehensive quality and risk assessment.

1.1 Sampling, Processing, and Analyses

9h

Water samples are collected (see Figure 2) and processed within a 🕒 06:00:00 interval, before being shipped for the partner responsible for the analyses (Table 1). In case no processing is needed, samples are frozen and stored at 🌡️ -80 °C within 🕒 03:00:00 .

For each sampling event, the date, day of the week and hour; the temperature and rain.

Sampling points, indicated in Figure 2 were designated from A to I:

- Influent of primary treatment (A)
- Influent of biological treatment (activated sludge) (B)
- Treated secondary effluent (C)
- Sand filter effluent (D)
- UVC effluent (E)
- Storage for reuse tank effluent (F)
- Ozonation effluent (1 dose, e.g., 🧪 5 mg O₃) - MITO3X technology - (G)

- [illegible]

1

1d 0h 30m

2

Analysis: Screening of Organic Contaminants in Water.

2.1

2.1 Sample filtration (**0.45um** PVDF).





2.2

2.3




3.1. Pass  5 mL of MeOH

3.2. Pass 5 mL of ultrapure water










- 2.4 Pass the  200 mL of sample
- 2.5 Rinse the volumetric flask or beaker that contained the sample with 2x  10 mL of ultrapure water, which are then passed through the cartridge
- 2.6 Drying the cartridge resin (e.g. N2 flow for  00:30:00) 30m
- 2.7 Store the cartridges in a freezer ( -20 °C) 1d

All samples are to be processed in duplicate.

For the cartridge blanks (n=2): start with  [go to step #2.3](#) , and then spike directly on the cartridge the same amount of internal standard as for samples ( 50 µL mix de  0 mL).

Then go directly to  [go to step #2.6](#) .

For the ultrapure water blanks (n=2), the same type of sampling flasks are filled with ultrapure water  24:00:00 before sampling. Then they are treated using the same protocols as for samples.

Elution of the loaded cartridges was carried out by gravity using  10 mL of methanol followed by  10 mL of dichloromethane. The eluates were concentrated approximately to ca.  0.5 mL using a Turbovap II concentrator (Zymark, Hopkinton, MA, USA), then to dryness under a gentle nitrogen stream and finally reconstituted in  500 µL of methanol and filtered with a  1 µg mL GHP%₀₀₀  13 mm syringe filter (Pall Corporation). The extracts were injected in an Agilent 1290 LC coupled with and Agilent 550 QTOF system.

Parameters framed by Legal and Regulatory Requirements:

- 3 Using the EU Drinking Water Directive:

Organic contaminants - DIRECTIVE 2008/105/EC OF THE EUROPEAN PARLIAMENT AND THE COUNCIL of 16 December 2008 on environmental quality standards in the field of water policy.



Protocol references

Montes, R.; Méndez, S.; Carro, N.; Cobas, J.; Alves, N.; Neuparth, T.; Santos, M.M.; Quintana, J.B.; Rodil, R. Screening of Contaminants of Emerging Concern in SurfaceWater and Wastewater. Effluents, Assisted by the Persistency-Mobility-Toxicity Criteria. *Molecules* 2022, 27, 3915. <https://doi.org/10.3390/molecules27123915>