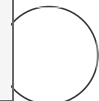




BARseq - high-throughput cell typing with in situ sequencing /v2

Xiaoyin Chen¹, Anthony M. Zador², Mararue¹
¹Allen Institute for Brain Science; ²Cold Spring Harbor Laboratory

 Xiaoyin Chen
Allen Institute for Brain Science

VERSION 2
AUG 16, 2023

ABSTRACT

This protocol describes the application of BARseq as a standalone in situ sequencing method to achieve multiplexed interrogation of endogenous genes. In this variation, BARseq is similar to in situ sequencing (ISS), but uses Illumina SBS for sequencing readout.

GUIDELINES

Standard precautions with RNA samples should be taken to reduce RNA degradation during tissue processing and library preparation. Pipetting and suctioning should be gentle throughout the whole procedure, and sample should not be left dried.

MATERIALS

MATERIALS

DNA oligos:

A	B	C
XC2757	/5AmMC12/NNNNNNNNNNNNNNNNNNNNNNNN	N20 for random priming
YS220	GATCGTCGGACTGTAGAACTCTGAACCTGTCTG	sequencing primer
YS221	/5Alex594N/GATCGTCGGACTGTAGAACTCTGAACCTGTCTG	Hybridization probe, for RPI gene detection
XC2758	/5Alex488N/AGTCAGCGTCGAGCACGCGGCACTTATTGCA	Hybridization probe, Slc17a7
XC2759	/5Alex532N/TGAGTAGAGTTGACTAAGAGCCGTTAGTGCC	Hybridization probe, Gad1
XC2760	/5Alex647N/TCGCTGTACTAATAGTTGTCTGACAGATCGTCA	Hybridization probe, Slc30

Reagents:

OPEN ACCESS



DOI:
[dx.doi.org/10.17504/protocols.io.81wgbp4j3vpk/v2](https://doi.org/10.17504/protocols.io.81wgbp4j3vpk/v2)

Protocol Citation: Xiaoyin Chen, Anthony M. Zador, Mararue 2023. BARseq - high-throughput cell typing with in situ sequencing. [protocols.io](https://doi.org/10.17504/protocols.io.81wgbp4j3vpk/v2) <https://doi.org/10.17504/protocols.io.81wgbp4j3vpk/v2> Version created by Xiaoyin Chen

MANUSCRIPT CITATION: <https://doi.org/10.1101/2022.11.06.515380>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Sep 15, 2022

Last Modified: Aug 16, 2023

PROTOCOL integer ID:
70065

Keywords: BARseq, in situ sequencing, spatial transcriptomics

Phusion high-fidelity PCR kit Thermo Scientific Catalog #F553S

Tween-20 Sigma-aldrich Catalog #P-7949

BS(PEG)9, 100 mg (Note: BS(PEG)9 loses its effectiveness 1 month after reconstitution in DMSO. Prepare a fresh batch every month, especially if it has been frozen and thawed repeatedly. Thermo Scientific Catalog #21582

Formamide Thermo Fisher Scientific Catalog #AM9342

10x PBS Thermo Fisher Scientific Catalog #AM9624

RNase-free water Contributed by users

dNTP Mix (dATP, dCTP, dGTP, and dTTP, each at 10mM) Thermo Fisher Scientific Catalog #R0192

phi29 DNA Polymerase (10 U/ μ L) Thermo Scientific Catalog #EP0091

RNase H Enzymatics Catalog #Y9220L

Glycerol Sigma Aldrich Catalog #G5516

Ethanol Merck Millipore Catalog #100983

Pierce™ MMTS (methyl methanethiosulfonate) Thermo Fisher Catalog #23011

SSC (20X), RNase-free Thermo Fisher Catalog #AM9770

RiboLock RNase Inhibitor (40 U/ μ L) Thermo Fisher Catalog #EO0381

RevertAid H Minus Reverse Transcriptase (200 U/ μ L) Thermo Fisher Catalog #EP0452

Paraformaldehyde 20% Electron Microscopy Sciences Catalog #15713

BSA Molecular grade New England Biolabs Catalog #B9000S

Ampligase DNA Ligase Kit Lucigen Catalog #A8101

KCl (2 M) RNase-free Thermo Fisher Scientific Catalog #AM9640G

Aminoallyl-dUTP Solution (50 mM) Thermo Fisher Scientific Catalog #R1101

Tris (1 M) pH 8.0 RNase-free Thermo Fisher Scientific Catalog #AM9855G

HiSeq SBS Kit v4 illumina Catalog #FC-401-4003

Grace Bio-Labs HybriWell-FL™ sealing system Fluor-friendly adhesive chamber Sigma Aldrich Catalog #GBL612204

Other equipment required include incubators set at 37 °C, 45 °C, and 60 °C. All tubes should be RNase-free. RNase-free filter tips should be used. A Crest Xlight v3 spinning disk confocal on an Nikon Ti2E with Photometrics Kinetix, and Lumencor Celesta was used for imaging the sequencing steps. The filters and lasers used are indicated in Table 1.

A	B	C	D
Channels	Laser	Dichroic	Emission filter
G/YFP	514	Zt405/514/63 5rpc	FF01-565/24
T/RFP	561	FF421/491/5 67/659/776-Di01	FF01-441/511/593/684/817
A	640	Zt405/514/63 5rpc	FF01-676/29
C	640	Zt405/514/63 5rpc	FF01-775/140
GFP	488	FF421/491/5 72-Di01	69401m
DAPI	405	FF421/491/5 72-Di01	69401m
TexasRed	561	FF421/491/5 72-Di01	69401m
Cy5	640	Zt405/514/63 5rpc	ZET532/640 m

Table 1. Laser and filter settings for sequencing imaging.

SAFETY WARNINGS

Use caution when handling liquids containing formaldehyde and formamide.

Library preparation

- 1 Tissues with barcoded neurons should be cryo-sectioned to 20 µm and mounted on slides.

Slides can be stored at -80 °C for up to a month.

2 DAY 1

Take slide(s) out of -80 °C and immerse immediately in 4% paraformaldehyde in 1x PBS (2 slides per 50mL falcon tube, back-to-back)

3 Incubate for 1 hour at room temperature on slow shaker

4 Wash the slides by immersing in 1x PBS (2 slides per 50ml falcon tube, back to back)

5 Wipe excess PBS off the surface of the chamber, then stick on the Hybriwell-FL chambers. Note that the ports on the chamber should be placed as far away from the tissue slices as possible.

6 Wash twice in PBST (1x PBS + 0.5% Tween-20)

7 Wash in 70% Ethanol for 5 mins

8 Wash in 85% Ethanol for 5 mins

9 Wash in 100% Ethanol for 5 mins

- 10 Replace with new 100% Ethanol, drop extra 100% Ethanol on top of slides and cover with ParaFilm to avoid evaporation. Incubate for at least 1.5 hrs at 4 °C (up to 3 hours)
- 11 Wash in PBST for 4-6 times, until all bubbles are cleared in the chamber and PBST flows into and out of the chamber smoothly.
- 12 Make reverse transcription mix: 50 µM N20 primer (XC2757), 20 U/µL RevertAid H Minus M-MuLV reverse transcriptase, 500 µM dNTP, 0.2 µg/µL BSA, 1 U/µL RiboLock RNase Inhibitor, 1x RevertAid RT buffer.
- 13 Incubate in reverse transcription mix overnight at 37 °C. Create a humidity chamber to avoid the slides drying out using kim-wipes and DI water.
- 14 **DAY 2:**

Wash with PBST once
- 15 Incubate in a mixture of 1µL BS(PEG)9 per 4 µL PBST (e.g. 200ul BS(PEG)9 and 800ul PBST) for one hour at room temperature
- 16 Wash with 1M Tris pH 8.0, then incubate in new 1M Tris pH 8.0 for 30 mins
- 17 Wash twice in PBST

- 18** Make ligation mix: 1x Ampligase buffer, 100 nM padlock probe each, 0.5 U/μL Ampligase, 0.4 U/μL RNase H, 1 U/μL RiboLock RNase Inhibitor, 50 mM KCl (extra of those already provided by the ampligase buffer), 20% formamide.
- 19** Incubate in ligation mix for at least 30 mins at 37 °C (can go longer but not shorter), then at least 45 mins at 45 °C (can go longer but not shorter).
- 20** Wash twice in PBST
- 21** Make RCA mix: 1 U/μL phi29 DNA polymerase, 1x phi29 polymerase buffer, 0.25 mM dNTP, 0.2 μg/μL BSA, 5% glycerol (extra of those from the enzymes), 125 μM aminoallyl dUTP
- 22** Incubate in RCA mix overnight at room temperature
- 23** **DAY 3:**

Wash with PBST once
- 24** Incubate in a mixture of 1 μL BS(PEG)9 per 4 μL PBST (e.g. 200ul BS(PEG)9 and 800ul PBST) for one hour at room temperature
- 25** Wash with 1M Tris pH 8.0, then incubate in new 1M Tris pH 8.0 for 30 mins
- 26** Wash twice in PBST

Sequencing

27 Hybridization of Gene sequencing primer:

Wash with FISH wash (2x SSC with 10% formamide)

28 Hybridize sequencing primer (YS220) with a primer concentration of 1 μ M in FISH wash for 10 mins at room temperature

29 Wash with FISH wash three times, 2 mins each

30 Wash with PBST twice

31 Sequence first cycle:

Do the following incubations. Unless noted with incubation temperature, each step is performed at room temperature. For steps without incubation time, treat these as quick washes. large flat metal blocks can be used to place sample slides to quickly cycle through high and low temperatures. This version uses MiSeq Nano v2 kit:

<https://www.illumina.com/products/by-type/sequencing-kits/cluster-gen-sequencing-reagents/miseq-reagent-kit-v2.html>

31.1 Incorporation Buffer 60 °C 3 mins x1

31.2 2% PBST x1

31.3 Idoacetamide blocker: For 9.3mg vial dilute pellet in between 2.5-3.5mL 2% PBST. Make fresh tube daily and store out of light.

Idoacetamide blocker 60 °C 3 mins x1

31.4 2% PBST x1

31.5 Incorporation Buffer x2

31.6 IRM 60 °C 3 mins x2

31.7 2% PBST x1

31.8 2% PBST 60 °C 3 mins x4

31.9 Replace 2% PBST with USM and Image **if slides are dirty, clean with 70% Ethanol before adding USM**

32 **Sequence subsequent cycles:** Do the following incubations. Unless noted with incubation temperature, each step is performed at room temperature. For steps without incubation time, treat these as quick washes. large flat metal blocks can be used to place sample slides to quickly cycle through high and low temperatures.

32.1 Incorporation buffer x2

- 32.2** CRM 60 °C 3 mins x2
- 32.3** Incorporation buffer x1 - Wipe ports after adding the incorporation buffer, to ensure that no CRM is left on the slide's surface
- 32.4** 2% PBST x1
- 32.5** Idoacetamide blocker: For 9.3mg vial dilute pellet in between 2.5-3.5mL 2% PBST. Make fresh tube daily and store out of light.
- Idoacetamide blocker 60 °C 3 mins x1
- 32.6** 2% PBST x1
- 32.7** Incorporation buffer x2
- 32.8** IRM 60 °C 3 mins x2
- 32.9** 2% PBST x1

32.10 2% PBST 60 °C 3 mins x4

32.11 Replace 2% PBST with USM and image **if slides are dirty, clean with 70% Ethanol before adding USM**

33 Hybridization cycle

33.1 Hybridize probes:

Make strip buffer: 60% formamide 2xSSC 0.01% Tween20

Strip buffer 60 °C 5 mins x3

Cool down quickly on metal plates between washes, place on metal plates in 60 °C oven to heat up quickly.

33.2 FISH wash (2x SSC with 10% formamide) 1x

33.3 Hybridize probes (YS221, XC2758, XC2759, XC2760) with a primer concentration of 1 µM in FISH wash at 60 °C for 2 minutes, then for 10 mins at room temperature. Rotate plates in holder to ensure they cool down slowly.

33.4 FISH wash x1

33.5 0.002 mg/ML DAPI in 2% PBST, room temperature for 5 mins

33.6 Replace PBST with USM and image ****if slides are dirty, clean with 70% Ethanol before adding USM****

33.7

33.8