



VERSION 1

DEC 29, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.n92ldzox7v5b/v1

Protocol Citation: freeman.lan 2023. DoTA-seq V3. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n92ldzox7v5b/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Aug 22, 2022


DoTA-seq V3 V.1

freeman.lan¹

¹UW Madison

DoTA-seq

Droplet Microfluidics

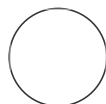

freeman.lan

ABSTRACT

This protocol describes the process of DoTA-seq generating a single cell sequencing library from a cell suspension. This workflow can be performed in two days, with the PCR step happening overnight. Before beginning this workflow make sure to have:

1. The necessary microfluidics devices prepared and ready to go
2. The multiplex DoTA-seq target primers validated to work together without generating large molecular weight primer dimers.

Please read the publication for further details.

MATERIALS



ddPCR Supermix for probes (no dUTP) BioRad
 Sciences Catalog #1863024



MetaPolyzyme Sigma
 Aldrich Catalog #MAC4L-5MG



Lysozyme from chicken egg white Sigma
 Aldrich Catalog #L6876



HFE 7500 Perfluorinated Oil



Perfluorooctanol Sigma
 Aldrich Catalog #370533



TCEP-HCl Gold
 Biotechnology Catalog #TCEP



NN'-Bis(acryloyl)cystamine Santa Cruz
 Biotechnology Catalog #sc-215506



Ammonium persulfate Catalog #A3678



Acrylamide P212121

Last Modified: Dec 29, 2023

PROTOCOL integer ID: 69011



TEMED (Tetramethyl-ethylenediamine) Sigma-
aldrich Catalog #T9281



Biorad Evagreen Droplet Oil BioRad
Sciences Catalog ##1864005



DNA Clean & Concentrator™-5 Zymo
Research Catalog #D4003



DNA Clean & Concentrator™-5 Zymo
Research Catalog #D4003



Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR
Tubes Corning Catalog #PCR-02-L-C



NEBNext Library Quant Kit for Illumina - 100 rxns New England
Biolabs Catalog #E7630S



SYBR Green Thermo Fisher
Scientific



Proteinase K solution, 20 mg ml -
1 Ambion Catalog #AM2546



Lysozyme from chicken egg white Sigma
Aldrich Catalog #L6876



Pre-injection buffer 10mM HEPES pH 7.5 Pluronic
0.1%

Safety information

Unpolymerized Acrylamide is toxic, handle with care and dispose according to regulations

PROTOCOL MATERIALS



Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR
Tubes Corning Catalog #PCR-02-L-C

Materials, Step 32





















EDTA (0.5 M), pH 8.0 Life
Technologies Catalog #AM9260G

Step 35



High Sensitivity D1000 ScreenTape Agilent
Technologies Catalog #5067-5584

Step 41

 PBS 2% Tween 20	Step 25
 DTT Merck MilliporeSigma (Sigma-Aldrich) Catalog #D0632	Step 29
 SYBR Green Thermo Fisher Scientific	Materials, Step 3
 TCEP-HCl Gold Biotechnology Catalog #TCEP	Materials, Step 38
 Biorad Evagreen Droplet Oil Bio-Rad Laboratories Catalog ##1864005	
In Materials and 3 steps	
 PBS EDTA 1mM 0.1% w/v Tween 20	Step 18
 SPRiselect reagent kit Beckman Coulter Catalog #B23317	Step 40
 Ammonium persulfate Catalog #A3678	Materials, Step 6
 Acrylamide P212121	Materials, Step 6
 TEMED (Tetramethyl-ethylenediamine) Merck MilliporeSigma (Sigma-Aldrich) Catalog #T9281	
Materials, Step 10	
 MetaPolyzyme Merck MilliporeSigma (Sigma-Aldrich) Catalog #MAC4L-5MG	
Materials, Step 20	
 PBS 0.1% Tween20	In 2 steps
 PBS 0.1% Tween 20	Step 2
 Cellbrite Fix 555 Biotium Catalog ##30088	Step 5
 Acetone	Step 16
 10% SDS Bio-Rad Laboratories Catalog #161-0146	Step 23
 DNA Clean & Concentrator™-5 Zymo Research Catalog #D4003	
In Materials, Materials, Step 39	
 Proteinase K solution, 20 mg ml – 1 Ambion Catalog #AM2546	
Materials, Step 23	



HFE 7500 Perfluorinated
Oil

Materials, Step 36



NEBNext Library Quant Kit for Illumina - 100 rxns New England
Biolabs Catalog #E7630S

Materials, Step 42



Perfluorooctanol Merck MilliporeSigma (Sigma-
Aldrich) Catalog #370533

In Materials and [2 steps](#)



Isopropanol Step 17



NN'-Bis(acryloyl)cystamine Santa Cruz
Biotechnology Catalog #sc-215506

Materials, Step 6



ddPCR Supermix for probes (no dUTP) Bio-Rad
Laboratories Catalog #1863024

Materials, Step 29



Lysozyme from chicken egg white Merck MilliporeSigma (Sigma-
Aldrich) Catalog #L6876

In Materials, Materials, Step 20





Pre-injection buffer 10mM HEPES pH 7.5 Pluronic
0.1%

In Materials and [2 steps](#)

Preparing Cells

10m

1

Prepare a cell suspension by washing twice in 1mL of  PBS 0.1% Tween20 Sigma
Aldrich by
spinning down at  5000 x g, 00:01:00

2


Resuspend cells in  100 μ L  PBS 0.1% Tween 20 Sigma
Aldrich

3

Add  1 μ L  SYBR Green Sigma
Aldrich 10,000X dye to the cells to stain them




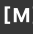







4

Count cells using a hemacytometer using the SYBR signal, calculate concentration of the cell
suspension.

- 5 (Optional) Stain with  Cellbrite Fix 555 Sigma Aldrich Catalog ##30088 to get a cell membrane/wall stain

Preparing Gel



30m

- 6 Make  200 μL Hydrogel Precursor Solution - Mix together in a tube:
-  100 μL  Acrylamide Sigma Aldrich monomer in water  25 Mass / % volume
-  15 μL  NN'-Bis(acryloyl)cystamine Sigma Aldrich Catalog #sc-215506 in Methanol  5 Mass / % volume
-  10 μL  Ammonium persulfate Sigma Aldrich Catalog #A3678  10 Mass / % volume
-  75 μL Cell suspension diluted in PBS (a total of **7e6 cells** to achieve a final concentration of **3.5e7 cells/mL**)

Vortex Vigorously to Mix

Generate Gel Droplets

10m

- 7 Prepare and Load the Syringes with the gel sample and  600 μL  Biorad Evagreen Droplet Oil Sigma Aldrich Catalog ##1864005 and connect to the microfluidic devices by following this protocol

Protocol



NAME

Loading Syringes to Inject into Microfluidics Device

CREATED BY

freeman.lan

PREVIEW

8 Run the syringe pumps at **600uL/hr** for the gel, and **1000uL/hr** for the oil syringe.

9 Collect gel droplets for  **00:20:00** in a 1mL tube.

20m

Note


Sometimes the initial droplet formation produces polydisperse droplets. In this case, wait 2 min for the bad emulsion to leave the outlet tubing into a waste tube, then begin collecting in the collection tube.



10 Make  **200 µL** Gel Polymerization Oil - mix together in a tube:

 **195 µL**

 Biorad Evagreen Droplet Oil Sigma
Aldrich Catalog ##1864005

 **5 µL**


 TEMED (Tetramethyl-ethulenediamine) Sigma
Aldrich Catalog #T9281

11 Add the Gel polymerization oil to the collected droplets and Incubate the tube containing droplets at  **37 °C** for  **00:10:00** to complete polymerization of the gel matrix.

10m

Note

You can now look at the emulsion under the microscope using



 Countess slides Thermo Fisher Scientific Catalog #C10228 to determine the encapsulation ratio of your cells. SYBRGreen and CF555 signal should be concordant and correspond to cells.


Breaking out gels from emulsion

12 Pulse spin the emulsion in a centrifuge to close pack the emulsion and drain the oil to the bottom of the tube.

13 Use a pipette to remove the oil at the bottom of the tube, leaving just the emulsion

14

Add  200 μL of  Perfluorooctanol Sigma Aldrich Catalog #370533 to break the emulsion



Vortex, then Wait  00:01:00 for the emulsion to break.


1m

15

Pulse spin again and remove the oil in the bottom of the tube with a pipette.

16



Add  1000 μL of  Acetone Sigma Aldrich to the tube.


Wait  00:00:05 then remove with a pipette.

The gels should begin to flocculate and dehydrate.

5s

17

Add  1000 μL of  Isopropanol Sigma Aldrich



Wait  00:00:05 then remove with a pipette.


The gels should dehydrate and become hard.

5s

Note: Do not wait too long as it could cause the gels to irreversibly aggregate into clumps.


18

Resuspend in  1000 μL of  PBS EDTA 1mM 0.1% w/v Tween 20 Sigma Aldrich

The gels can be stored at  4 $^{\circ}\text{C}$ for several days without changing DoTA-seq results.






Note

You can now look at the gels under the microscope using

 Countess slides Thermo Fisher Scientific Catalog #C10228 to determine the encapsulation ratio of your cells.

You should see some loss in CF555 signal as the acetone and alcohol wash removes some bacterial membranes.


Lysing Bacteria

- 19 Wash gels 3 times in  1000 μL  1X PBS (Phosphate-buffered saline) Sigma Aldrich (No Tween) by centrifugation at  500 x g, 00:00:15 each time 15s
- 20 Make a Enzymatic Lysis Solution by adding:
-  20 mg  Lysozyme from chicken egg white Sigma Aldrich Catalog #L6876
-  100 μL [M] 1 mg/mL  MetaPolyzyme Sigma Aldrich Catalog #MAC4L-5MG
-  900 μL  1X PBS (Phosphate-buffered saline) Sigma Aldrich (No Tween)
- 21 Resuspend the gels in this lysis solution. Incubate at  37 $^{\circ}\text{C}$ for  01:00:00 1h
- 22 Wash the gels 3 times in  1000 μL  PBS 0.1% Tween20 Sigma Aldrich
- 23 Make a SDS Lysis solution by adding:
-  20 μL [M] 20 mg/mL  Proteinase K solution, 20 mg ml – 1 Sigma Aldrich Catalog #AM2546
-  100 μL  10% SDS Sigma Aldrich Catalog #161-0146
-  880 μL  1X PBS (Phosphate-buffered saline) Sigma Aldrich
- 24 Resuspend the gels in this SDS Lysis solution, incubate at  55 $^{\circ}\text{C}$ for  01:00:00 1h
- 25 Wash the gels three times in  1000 μL  PBS 2% Tween 20 Sigma Aldrich
- Note: Use PBS **2% Tween 20**, not 0.1% Tween

These gels can be stored at  4 °C for several days without impacting DoTA-seq results.

Note


You can now look at the gels under the microscope using


 Countess slides Thermo Fisher Scientific Catalog #C10228 to determine the encapsulation ratio and lysis efficiency of your cells.


It is advised to restrain with SYBR and CF555 to get best signal. Lysed cells should exhibit SYBR signal but no CF555 Signal.


Barcoding the Cells

7m

26 Wash the gels three times in  1000 µL

 Pre-injection buffer 10mM HEPES pH 7.5 Pluronic 0.1% Sigma
Aldrich

27 Resuspend gels in  100 µL

 Pre-injection buffer 10mM HEPES pH 7.5 Pluronic 0.1% Sigma
Aldrich

28 Load the gels into a syringe following the protocol described in this excellent visual protocol.

CITATION



Demaree B, Weisgerber D, Lan F, Abate AR (2018). An Ultrahigh-throughput Microfluidic Platform for Single-cell Genome Sequencing.. Journal of visualized experiments : JoVE.

LINK

<https://doi.org/10.3791/57598>

29 Generate a PCR Master Mix (This mix gives about ~10,000 cells per library - Scale up as required)

20m

 25 µL  ddPCR Supermix for probes (no dUTP) Sigma
Aldrich Catalog #1863024

 0.4 µL [M] 50 micromolar (µM) P7 Primer

 0.4 µL [M] 50 micromolar (µM) P5 Primer with appropriate I5 index

🧪 0.2 µL Variable	[M] 10 micromolar (µM)	DoTA-seq multiplex primer mix (10uM concentration per primer)
🧪 0.2 µL Variable	[M] 10 micromolar (µM)	16S DoTA-seq primers
🧪 0.5 µL Variable	[M] 1 picomolar (pM)	Freshly diluted from 500pM stock Barcode Oligo
🧪 0.25 µL	[M] 500 millimolar (mM)	Single use aliquots
		⚗ DTT Sigma Aldrich Catalog #D0632

Note

The ratio of 16S to DoTA-seq target primers mix can be varied depending on the relative amplification efficiencies. The best way to determine is to start from equal concentrations, then adjust based on the sequencing results (do most cells contain more 16S reads than target reads?)

Note

Typically, 0.5uL of 1pM barcode will give approximately 1 barcode for every 10 droplets. However, it is best to measure the barcode encapsulation rate by making PCR droplets containing the barcodes at an estimated dilution and P7 and Barrev primers targeting the barcode for amplification. Visualize the resulting PCR emulsion using SYBRgreen staining under the microscope to obtain the real encapsulation ratio.

Note

Barcode oligos should always be freshly diluted from 500pM to 1pM before use, as we have found gradual loss of barcodes over time in a 1pM solution.

30 Load the PCR mastermix into the syringe following this protocol

Protocol





NAME

Loading Syringes to Inject into Microfluidics Device

CREATED BY

freeman.lan

PREVIEW

- 31 Load  500 μL of  Biorad Evagreen Droplet Oil Sigma Aldrich Catalog ##1864005 into a syringe following this protocol

Protocol






NAME

Loading Syringes to Inject into Microfluidics Device

CREATED BY

freeman.lan

PREVIEW

- 32 Run the syringe pumps at 200uL/hr for the gel and PCR mastermix, and 900uL/hr for the oil syringe. 7m
Collect droplets in an  Axygen® 0.2 mL Maxymum Recovery® Thin Wall PCR TubesSigma Aldrich Catalog #PCR-02-L-C
for  00:07:00 for every  25 μL of PCR mastermix or until the PCR mastermix runs out.

- 33 Use a pipette to remove the oil in the PCR tube, leaving just the emulsion layer

- 34 Thermocycle the PCR emulsion as follows: 4h

 95 °C 5 min

20 cycles of:

 95 °C 30s

 72 °C 10s

 60 °C 5 min

 72 °C 30s

20 cycles of:

 95 °C 30s

🔥 72 °C 10s

🔥 60 °C 90s

🔥 72 °C 30s

Final incubation of:

🔥 72 °C 10min

🔥 12 °C Hold

All ramp times are at 🔥 1 °C per second

PCR Cleanup

1m

35 Keep the emulsion on ice to prevent polymerase activity

Add 🧪 25 µL ⚗️ EDTA (0.5 M), pH 8.0 Sigma
Aldrich Catalog #AM9260G to the emulsion

Vortex the emulsion to mix

36 Add 🧪 25 µL ⚗️ HFE 7500 Perfluorinated Oil Sigma
Aldrich to the emulsion

Add 🧪 25 µL ⚗️ Perfluorooctanol Sigma
Aldrich Catalog #370533 to the emulsion

Vortex the emulsion to mix


37 Wait ⌚ 00:01:00 , then pulse centrifuge to separate the PCR mix from the oil
Transfer the top aqueous phase to a new 1mL tube.

1m


38 Add 🧪 20 µL [M] 1 Molarity (M) ⚗️ TCEP-HCl Sigma
Aldrich Catalog #TCEP to the tube and
vortex to completely decrosslink any remaining gels.

Note


You should be unable to obtain any "jellyish" substance by centrifugation! If there is any jellyish substance left it is not fully de-crosslinked. **ADD MORE TCEP.**

39 Clean up the PCR reaction using the  DNA Clean & Concentrator™-5 Sigma
Aldrich Catalog #D4003 10m

kit.

Elute in  50 µL Elution Buffer.

40 Remove primer dimers and free barcodes using the 10m

 SPRIselect reagent kit Sigma
Aldrich Catalog #B23317 with 0.7X volume of beads.

41 Check the resulting library for primer dimers using 10m

Equipment

TapeStation

NAME

Agilent

BRAND

G2991AA

SKU


<https://www.agilent.com/en/product/tapestation-automated-electrophoresis/tapestation-instruments/4200-tapestation-system-228263>

LINK

with a  High Sensitivity D1000 ScreenTape Sigma
Aldrich Catalog #5067-5584

Other high sensitivity capillary electrophoresis methods will also work.

There should be minimal primer dimers on the trace. Below is an example of an acceptable trace.

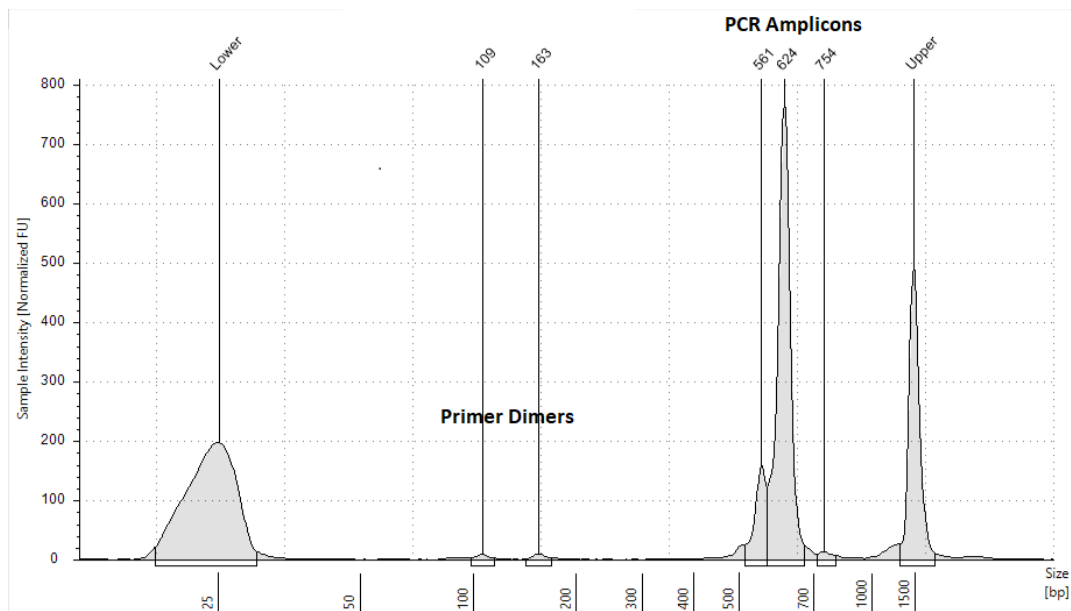

Example of an acceptable Tapestation trace.

42 Quantify the library using a qPCR library quantification kit such as

1h

NEBNext Library Quant Kit for Illumina - 100 rxns Sigma
Aldrich Catalog #E7630S

Note

Note that you must use a PCR based library quantification kit as not all amplicons contain all the adaptors for sequencing and therefore will throw off sequence non-specific forms of quantification!

43 Sequence the library on an Illumina sequencer using Custom Sequencing Primers listed here.

DoTA-seq-Oligo-Sequences.xlsx