



JAN 15, 2024

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.3byl4brwrvo5/v1](https://dx.doi.org/10.17504/protocols.io.3byl4brwrvo5/v1)

**Protocol Citation:** andrew. west, arpine.sokratian 2024. a-Synuclein protein expression and purification. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.3byl4brwrvo5/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working  
 We use this protocol and it's working

**Created:** May 21, 2022

## 🌐 a-Synuclein protein expression and purification

andrew.west<sup>1</sup>, arpine.sokratian<sup>1</sup>

<sup>1</sup>Duke University

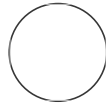

andrew.west

### ABSTRACT

Protocol for recombinant a-synuclein purification, useful as monomer template for seeded-amplification assays (like RT\_QUIC or PMCA) . Recommendations are to store the protein always on ice while not running and do not stop purification when it is started.

### MATERIALS

1. BL21-CodonPlus (DE3)-RIL Chemical Competent Cells (Agilent #230245-41)
2. Thermo Scientific™ Low Protein Binding Collection Tubes (1.5 mL) PI90411



Nalgene™ Rapid-Flow™ Sterile Single Use Vacuum Filter Units Thermo Fisher Scientific Catalog #565-0010



SnakeSkin Dialysis Tubing 3.5K MWCO 35 Thermo Fisher Scientific Catalog #88244



Endotoxin detection kit LAL Genscript Catalog #95045-024



ToxinEraser™ Endotoxin Removal Kit Genscript Catalog #89233-330



ToxinEraser™ Endotoxin Removal Resin Genscript Catalog #L00402



HiPrep Q HP anion exchange chromatography column Cytiva Catalog #29018182



MilliporeSigma™ Amicon™ Ultra-15 Centrifugal Filter Units Catalog #MilliporeSigma™ UFC901024


Last Modified: Jan 15, 2024

PROTOCOL integer ID:  
62989


Keywords: ASAPCRN

Funders  
Acknowledgement:  
Aligning Science Across  
Parkinson’s  
Grant ID: ASAP-020527


PROTOCOL MATERIALS

 Nalgene™ Rapid-Flow™ Sterile Single Use Vacuum Filter UnitsThermo Fisher Scientific Catalog #565-0010


Materials, Step 21

 SnakeSkin Dialysis Tubing 3.5K MWCO 35 Thermo Fisher Scientific Catalog #88244


Materials, Step 22

 Endotoxin detection kit LAL Genscript Catalog #95045-024


Materials, Step 43

 ToxinEraser™ Endotoxin Removal Kit Genscript Catalog #89233-330


Materials, Step 41

 ToxinEraser™ Endotoxin Removal Resin Genscript Catalog #L00402

Materials, Step 41

 HiPrep Q HP anion exchange chromatography column Cytiva Catalog #29018182

Materials, Step 24

 MilliporeSigma™ Amicon™ Ultra-15 Centrifugal Filter Units Catalog #MilliporeSigma™ UFC901024


Materials, Step 35

Transformation


1d

1


Thaw down an aliquot of plasmid construct (pRK172) encoding WT-human-a-synuclein [M] 0.3 mg/r... 15m


 On ice

2

Thaw down  On ice an aliquot of BL21 (DE3) RIL competent E Coli cells 15m

3

Add  1 µL of plasmid construct to the thawed competent cells and gently mix by flicking the bottom of the tube with a finger a few times




Safety information


do not resuspend

3.1 Incubate the reaction mix  On ice 15m



4 Perform heat-shock transformation  42 °C in water bath incubator with manually shaking at 1m



 100 rpm, 00:00:45

#### Equipment

**Precision™ General Purpose Water Bath**

NAME

Water Bath

TYPE


Thermo Scientific

BRAND

TSGP10

SKU

5 Immediately transfer the tube on ice and incubate for 1 min. 1m

6 Add  1000 µL of SOC media to a chilled reaction 10s

7 Incubate the bacteria  200 rpm, 37°C, 00:30:00 30m

### Equipment

ThermoMixer® C

NAME

ThermoMixer

TYPE

Eppendorf

BRAND

5382000023

SKU

7.1 Prepare sterile 10cm LB agar plate containing **[M] 0.1 mg/mL** of ampicillin

8 Centrifuge at **[S] 500 x g, 10°C, 00:03:00** and discard the supernatant leaving **[A] 50 µL** of media **3m**

8.1 Spread 50 ul of cell suspension onto a selection plate and incubate overnight at **[T] 37 °C** in **10m** bacterial incubator



### Equipment

Isotemp™ Microbiological Incubator, 178 L, Stainless Steel

NAME

Microbiological Incubator

TYPE

Fisherbrand



BRAND


15-103-0513

SKU


## Protein expression

12h




9 Pick one colony and transfer into  10 mL LB media with  0.1 mg/mL of ampicillin start in the morning (9:00 am)

9.1 Incubate the bacteria  250 rpm, 37°C, 05:00:00 until it reaches OD 0.2-0.3 5h

Equipment	
Natural convection incubator	NAME
Bacterial shaker	TYPE
Innova	BRAND
M1335-0000	SKU

10 Transfer a starter culture to 2X2L flasks filled with 0.5L LB media with  0.1 mg/mL of ampicillin 5 mL to each flask

11 Incubate the culture at the same conditions until it reaches OD 0.8 (use nanodrop or cuvette) (reaches optimal density at 6-7 pm) 5h

12 Induce protein expression by adding  0.05 millimolar (mM) IPTG, incubate at  18 °C for 12h  
 12:00:00 overnight




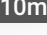






















Note

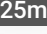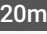To cool down the grown culture, transfer the flasks into ice-bath and incubate until it reaches desired temperature

## Cell lysis

11m 45s

- 13** Collect the pellets by centrifugation (JA14 rotor) at  5000 x g x g at  4 °C for  00:10:00 .  10m  
Used 250 ml Beckman tubes  
Usually get 10-12 g from 2L
- 14** Add to pellets 80 ml of lysis buffer (total):  10 millimolar (mM) TrisHCl  7.6 ,  
 750 millimolar (mM) NaCl,  1 millimolar (mM) EDTA,  1 millimolar (mM) PMSF (add just before using, have aliq frozen  0.1 Molarity (M) ), protease inhibitors (use MAXI version, need only one tablet);
- 15** Carefully resuspend the pellets to homogenize the solution
- 15.1** Heat up  1 L of water in a high temperature resistant glass beaker (turn heat to the max on the magnetic stirrer)
- 16** While waiting on water to get to the boiling point sonicate the lysates (use thick prob-tip) for  00:01:00 , 30%,  00:00:15 ON  00:00:30 OFF of amplitude then go to next falcon, had 3 falcons (repeat 3 times, avoid overheating)  10m
- 17** After sonication samples need to get boiled thereby put the falcon tubes into glass beaker and boil for  00:25:00 . Use tweezers to pull out the tubes  25m
- 18** Transfer boiled homogenates into new 50 mL falcon tubes; chill down suspensions at room temperature for 20 min  20m
- 19** Prepare  4 L of buffer  10 millimolar (mM) TrisHCl  7.6 ,  50 millimolar (mM) NaCl,  1 millimolar (mM) EDTA,  1 millimolar (mM) PMSF for dialysis


20

Centrifuge the homogenates at  20000 x g for  01:00:00 at  4 °C

1h


21

Filter the supernatant using 0.45 um filter unit


 Nalgene™ Rapid-Flow™ Sterile Single Use Vacuum Filter UnitsThermo Fisher Scientific Catalog #565-0010


22

Transfer filtered supernatant into dialysis bag which is: SnakeSkin Dialysis tubing, 3.5K MWCO, 35 mm dry I.D., 35 feet.



*Measure the dialysis tube taking into consideration that 5 cm length of tube holds 48 mL of the sample (plus 2.5cm at each end for closure). Clip the tube using green clips, make sure it does not leak.*


Place the dialysis bag into  4 L plastic beaker filled with dialysis buffer, incubate overnight on magnetic plate on the slow mode (Chromatography fridge)

 SnakeSkin Dialysis Tubing 3.5K MWCO 35 Thermo Fisher Scientific Catalog #88244

Equipment			
ÄKTA pure 25 L1			NAME
ÄKTA pure chromatography system			TYPE
Cytiva			BRAND
29018225			SKU

Protein purification (anion-exchange chromatography)



23

After a night of dialysis ( 4 °C slow mixing) collect the suspension into 100 mL glass bottle (filter the sample before running on the column, 0.22 um filter).

24

Column - HiPrep Q HP 16/10 column 1x20 ml (stored in 70% ethanol);



- 24.1 Wash the column 2V of miliQ degassed water
- 24.2 Wash the column with 2V of **STARTING BUFFER** [M] 10 millimolar (mM) TrisHCl  $\text{pH}$  7.6 ,  
[M] 50 millimolar (mM) NaCl
- 24.3 Activate with 1V of [M] 10 millimolar (mM) TrisHCl  $\text{pH}$  7.6 , [M] 1 Molarity (M) NaCl
- 24.4 Equilibrate with 3V of starting buffer
- 25 Load  80 mL of suspension and then washed with 100 ml [M] 50 millimolar (mM) NaCl [M] 10 millimolar (mM) TrisHCL,  300 mL of gradient elution (0-100%), 2 ml/min flow rate. Collected samples using fraction collector 2, every fraction 4 ml (use 10 ml glass tubes)
- 26 Place supernatant into channel A1 (was previously use for starting buffer, do not generate bubbles)
- 27 Place starting buffer in channel A2 (clean the tubing using the program mode)



- 28 Place elution buffer in channel B1 ([M] 10 millimolar (mM) TrisHCl  $\text{pH}$  7.6 , [M] 1 Molarity (M) NaCl)  
Collected samples using fraction collector 2, every fraction 4 ml (use 10 ml glass tubes);
- 29 Analyze the fractions eluted at 250-350 mM salt (20 RFU conductivity) though SDS-PAGE (stain with Coomassie).  
*Combine  $\text{A}$  10  $\mu\text{L}$  of each fraction with  $\text{A}$  10  $\mu\text{L}$  of 2X laemmli buffer and analyze fractions by SDS-PAGE with 4–20% gradient gels, followed by coomassie staining/destaining*
- 30 Measure A280/260 for the fractions containing single a-syn band, avoid collecting samples with A280/260 > 0.85
- 31 Combine the evaluated factions and measure total protein concentration using nanodrop.
- 32 Dialyze with  $\text{A}$  4 L of [M] 10 millimolar (mM) TrisHCl  $\text{pH}$  7.6 , [M] 50 millimolar (mM) NaCl  
(follow instruction for dialysis)



## Further purification














- 33 Repeat section 'Protein purification (anion-exchange chromatography)' for the further fractionation of the purified preparation [↺ go to step #23](#)



## Protein concentration



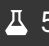



10m

- 34 Concentrate dialyzed protein sample to approximately [M] 30 mg/mL aliquot  
Prepare the ultra-concentration system
- 35 Use 50 mL ultra centrifugation units with 3K cutoff

- 36 Wash off the unit with miliQ water through centrifugation at  5000 x g at  4 °C for  00:05:00  5m, JA10 rotor
- 37 Load first  15 mL of the sample into ultracentrifugation unit (max load of the unit is approx.  15 mL )
- 38 Centrifuge at  5000 x g at  4 °C for  00:05:00, JA10 rotor  5m
- 39 Resuspend concentrated sample, add more of protein sample and concentrate until the total volume is ~  5 mL
- 40 Store at  -80 °C. Yield should be approximately  80 mg per 2 L culture



## Endotoxin removal

- 41 Follow instructions for  ToxinEraser™ Endotoxin Removal Kit Genscript Catalog #89233-330 with modifications  
For a more successful endotoxin removal, add  1 mL of  ToxinEraser™ Endotoxin Removal Resin Genscript Catalog #L00402 before the regeneration process
- 42 Collect the eluate into  5 mL endotoxin-free tube and save 2 aliquots (  10 µL and  50 µL ) for protein concentration and endotoxin measurements

## Endotoxin quantification

43

Follow instructions for Endotoxin detection kit LAL Genscript



Endotoxin detection kit LAL  
Genscript Catalog #95045-024