



Oct 13, 2022

DNA extraction for fermented plant based foods


Anique Ahmad¹, Arya Gautam¹, Tsedenia Deneke¹, Ahmed Shibl¹, Aashish Jha¹

¹New York University, Abu Dhabi

Anique Ahmad: aa5506@nyu.edu

Aashish Jha: jhaar@nyu.edu

1 Works for me

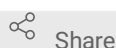

Share

dx.doi.org/10.17504/protocols.io.36wgqj9wkvk5/v1

nyuad.geneticheritagegroup

ABSTRACT

DNA Extraction for fermented plant based foods

DOI

dx.doi.org/10.17504/protocols.io.36wgqj9wkvk5/v1

PROTOCOL CITATION

Anique Ahmad, Arya Gautam, Tsedenia Deneke, Ahmed Shibl, Aashish Jha
2022. DNA extraction for fermented plant based foods. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.36wgqj9wkvk5/v1>

FUNDERS ACKNOWLEDGEMENT

The Bill and Melinda Gates Foundation

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 19, 2022

LAST MODIFIED

Oct 13, 2022

GUIDELINES

Disinfect mortar and pestle with 70% ethanol in between samples

If the goal is metagenomic sequencing, we recommend bypassing the mortar and pestle step. Simply suspend the samples in 1 ml of ZymoBIOMICS™ Lysis Solution, shake vigorously or vortex, and transfer 750 µL of the supernatant to a ZR BashingBead™ Lysis Tube (0.1 & 0.5 mm) and proceed to Step 2.

MATERIALS TEXT

ZymoBIOMICS DNA Extraction Kit

Fecal samples

Mortar and pestle

Liquid nitrogen

Ice box with dry ice

1.5 ml microcentrifuge tubes

Vortex

Omni Bead Ruptor Elite bead beater



Centrifuge



SAFETY WARNINGS

Handle liquid nitrogen with proper training and PPE

BEFORE STARTING

- Label the following:
 - ZR BashingBead™ Lysis Tubes (0.1 & 0.5 mm)
 - ZymoSpin™ III-F Filter in Collection Tube
 - Collection Tube
 - ZymoSpin™ II-CR Column in Collection Tube
 - ZymoSpin™ III-HRC Filter in Collection Tube
 - 2 Sets of 1.5 ml microcentrifuge tubes (not provided with kit)
- Spin the labelled ZR BashingBead™ Lysis Tubes (0.1 & 0.5 mm) for 10 seconds in a mini-centrifuge to ensure that the beads have settled at the bottom
- Include 1 control per batch and assign its position randomly.

- 1 Take  **250 mg** of food sample while on dry ice and transfer to mortar. Add liquid nitrogen until food sample is completely immersed. Using pestle vigorously homogenize sample until a powder forms. Transfer all of the sample to a *ZR BashingBead™ Lysis Tube (0.1 & 0.5 mm)*. Add  **750 µL** ZymoBIOMICS™ Lysis Solution to the tube. Cap tightly.

- 2 Secure in a Omni Bead Ruptor Elite bead beater fitted with a 2 ml tube holder assembly and ^{10m} process at max speed (30 m/s) for  **00:05:00** . Rest for  **00:05:00** .

- 3 Centrifuge the *ZR BashingBead™ Lysis Tubes (0.1 & 0.5 mm)* in a microcentrifuge at \geq **10000 x g, 00:01:00** . 1m
- 4 Transfer up to **600 µL** supernatant to the *Zymo-Spin™ III-F Filter* in the labelled *Collection Tube* and centrifuge at **8.000 x g, 00:01:00** . Discard the Zymo-Spin™ III-F Filter. 1m
- 5 Add **600 µL** of ZymoBIOMICS™ DNA Binding Buffer to the filtrate in the labelled *Collection Tube* from Step 4. Mix well.
- 5.1 Repeat Step 5 so that the final volume of ZymoBIOMICS™ DNA Binding Buffer added is **1200 µL** .
- 6 Transfer **800 µL** of the mixture from Step 5 to a *Zymo-Spin™ IICR Column* in a *Collection Tube* and centrifuge at **10.000 x g, 00:01:00** . 1m
- 6.1 Discard the flow through from the Collection Tube and repeat Step 6.
- 7 Add **400 µL** ZymoBIOMICS™ DNA Wash Buffer 1 to the *Zymo-Spin™ IICR Column* in a *new Collection Tube* and centrifuge at **10.000 x g, 00:01:00** . Discard the flow-through. 1m
- 8 Add **700 µL** ZymoBIOMICS™ DNA Wash Buffer 2 to the *Zymo-Spin™ IICR Column* in a *Collection Tube* and centrifuge at **10.000 x g, 00:01:00** . Discard the flow-through. 1m
- 9 Repeat with **200 µL** ZymoBIOMICS™ DNA Wash Buffer 2 to the *Zymo-Spin™ IICR Column* in a *Collection Tube* and centrifuge at **10.000 x g, 00:01:00** . 1m
- 10 Transfer the *Zymo-Spin™ IICR Column* to a clean 1.5 ml microcentrifuge tube and add **100 µL** (50 µl minimum) ZymoBIOMICS™ DNase/RNase Free Water directly to the column 6m

matrix and incubate for 🕒00:05:00 . Centrifuge at 🌀10.000 x g, 00:01:00 to elute the DNA

- 11 Place a *Zymo-Spin™ III-HRC Filter* in a new *Collection Tube* and add 📏600 µL 3m
ZymoBIOMICS™ HRC Prep Solution. Centrifuge at 🌀8.000 x g, 00:03:00 .
- 12 Transfer the eluted DNA (Step 10) to a prepared *Zymo-Spin™ III-HRC Filter* in a clean 1.5 ml ^{3m}
microcentrifuge tube and centrifuge at exactly 🌀16.000 x g, 00:03:00 .
- 13 The filtered DNA is now suitable for PCR and other downstream applications. Eluted DNA should be frozen (–30 to –15°C or –90 to –65°C)