




DEC 15, 2023

pH-Rodo pathogens to monitor phagocytosis in pMacs

Rebecca
Wallings¹

¹University of Florida

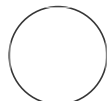

Rebecca Wallings

ABSTRACT

pH-Rodo pathogens to monitor phagocytosis in pMacs

OPEN  ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.q26g7pn18gwz/v1

Protocol Citation: Rebecca Wallings 2023. pH-Rodo pathogens to monitor phagocytosis in pMacs.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.q26g7pn18gwz/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Dec 14, 2023

Last Modified: Dec 15,
2023

PROTOCOL integer ID:
92331

Funders

Acknowledgement:

ASAP

Grant ID: ASAP-020527

- 1 Once pMacs have adhered to plate, resuspend pHrodo™ Green *E. coli* BioParticles™ in Invitrogen™ Live Cell Imaging Solution to a final concentration of 1mg/mL and sonicate in bath sonicator for 10 minutes.

Add 100 µL of *E. coli* BioParticles solution to each well (96-well plate).

Place cells in a live imaging chamber (37°C, 5% carbon dioxide, and 95% humidity) of an EVOS™ M7000 (Invitrogen) at 20 x magnification capture 4 images per well every 20 minutes for 5 hours.

Use brightfield images for automated-focus throughout the time course and use for normalization of GFP *E. coli* images.