



MAY 08, 2023

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.x54v9dj4zg3e/v1

Protocol Citation: Dan Dou, Alexander Boecker, Erika L.F. Holzbaur 2023. Culture and transfection of iPSC-derived neurons for live-imaging of axonal cargoes. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.x54v9dj4zg3e/v1>

MANUSCRIPT CITATION: Boecker, C.A., and Holzbaur, E.L.F. (2021). Hyperactive LRRK2 kinase impairs the trafficking of axonal autophagosomes. *Autophagy* 00, 1–3. Boecker, C.A., Olenick, M.A., Gallagher, E.R., Ward, M.E., and Holzbaur, E.L.F. (2020). ToolBox: Live Imaging of intracellular organelle transport in induced pluripotent stem cell-derived neurons. *Traffic* 21, 138–155. Fernandopulle, M.S., Prestil, R., Grunseich, C., Wang, C., Gan, L., and Ward, M.E. (2018). Transcription Factor-Mediated Differentiation of Human iPSCs into Neurons. *Curr. Protoc. Cell Biol.* 79, e51.


🌐 Culture and transfection of iPSC-derived neurons for live-imaging of axonal cargoes

Dan Dou^{1,2}, Alexander Boecker³, Erika L.F. Holzbaur^{1,2}

¹Department of Physiology, Perelman School of Medicine, University of Pennsylvania, Philadelphia, PA 19104, USA;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, USA;

³Department of Neurology, University Medical Center Goettingen, 37077 Goettingen, Germany

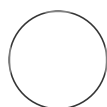

Dan Dou

ABSTRACT

Here, we plate, culture, and transfect human iPSC-derived excitatory glutamatergic neurons for the purpose of observing transport of axonal cargoes under spinning disk confocal microscopy. Protocol is largely as previously described (Boecker et al., 2020, 2021; Fernandopulle et al., 2018). For preceding differentiation of neurons, see “Protocol: Piggybac-mediated stable expression of NGN2 in iPSCs for differentiation into excitatory glutamatergic neurons” and “Protocol: iNeuron differentiation from human iPSCs.”

ATTACHMENTS

[550-1146.pdf](#)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Oct 13, 2022

Last Modified: May 08, 2023

PROTOCOL integer ID: 71294

Keywords: iPSC, iNeuron, live-imaging, axon

MATERIALS

Materials

Equipment	
35 mm Dish No. 1.5 Coverslip 20 mm Glass Diameter Uncoated	NAME
Coverslip	TYPE
Mattek	BRAND
p35g-1-5-20-c	SKU
https://www.mattek.com/store/p35g-1-5-20-c-case/	LINK

Reagents

- PLO (CATALOG)
- Borate buffer (CATALOG)
- BrainPhys media (CATALOG)
- NT-3 (CATALOG)
- BDNF (CATALOG)
- B-27 supplement (CATALOG)
- Mouse laminin (CATALOG)
- 5-Fluoro-2'-deoxyuridine
- Uridine

Culture and transfection of iPSCderived neurons for live-im...30m

1 In advance, prepare 10x PLO stock.



A	B
PLO	50 mg
0.1M borate buffer	50 mL

Note

Store 10x PLO stock at  -80 °C .

2



The day before plating, coat imaging dishes with 1x PLO solution (10x PLO stock diluted in ddH₂O).

Note

It is only necessary to fully coat the glass center of the imaging dish.

3



The day of plating, remove PLO solution from imaging dishes and wash twice with ddH₂O.

3.1


Add  2 mL of iNeuron culture media.



BrainPhys supplemented with

A	B
BDNF	10 ng/mL
NT-3	10 ng/mL
Laminin	1 µg/mL
B-27 supplement	1x


3.2

Place dishes in cell culture incubator for >  00:30:00 .

30m



4

Rapidly thaw cryopreserved iNeurons in  37 °C water bath.

Note

Retrieve vial to tissue culture hood when only a small amount of ice remains visible.

- 5 Centrifuge to remove freezing media and resuspend cell pellet in iNeuron culture media.



BrainPhys supplemented with

A	B
BDNF	10 ng/mL
NT-3	10 ng/mL
Laminin	1 µg/mL
B-27 supplement	1x

- 6 Count cells and plate 300k neurons per 35 mm imaging dish.



- 6.1 Add cells dropwise to the center area of the dish (so that they sink onto the glass, PLO-coated center).

- 7 For Piggybac-delivered NGN2 neurons, include 10 micromolar (µM) 5-Fluoro-2'-deoxyuridine and 10 micromolar (µM) uridine at the time of plating to prevent survival of mitotic cells.

Note

These drugs were removed 24 hours after plating.

- 8 Store neurons in cell culture incubator. Perform partial change of iNeuron media twice per week.





9

On DIV18 (~  72:00:00 prior to imaging), transfect iNeurons for imaging.

3d



Note

Transfection conditions may require optimization, but a typical transfection will use  4 μL Lipofectamine Stem and  1 μg of plasmid DNA. Plasmids with the PGK or EF1 α promoters express best in iNeurons.