



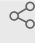
Jun 21, 2022

Inserting RFID into rats

a.beeson¹, Leah Solsberg-Woods², Gabriel J Barrero¹, Oksana Polesskaya¹, Abraham Palmer¹

¹University of California, San Diego; ²Wake Forest University

1 *Works for me*

 Share

dx.doi.org/10.17504/protocols.io.rm7vzyxzrlx1/v1

CGORD



Anita Bandrowski
SciCrunch Inc and UCSD

ABSTRACT

Insert RFID chips into rats

DOI

dx.doi.org/10.17504/protocols.io.rm7vzyxzrlx1/v1

PROTOCOL CITATION

a.beeson , Leah Solsberg-Woods, Gabriel J Barrero, Oksana Polesskaya, Abraham Palmer 2022. Inserting RFID into rats. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.rm7vzyxzrlx1/v1>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Mar 09, 2022

LAST MODIFIED

Jun 21, 2022

PROTOCOL INTEGER ID

59274

GUIDELINES

Save all old cage cards for verification of DOB. Weaned Pups will be housed in Coleman B113 for at least 2 days before transported to TSRB 155.

MATERIALS TEXT

- . Syringe with needle and non-sterile RFIDs (soaked in [M]75 % volume EtOH
- 🕒 **Overnight** or longer, then air dry before use)
- . Funnel bag

- 1 Disinfect needle with 70% EtOH and air dry
- 2 Carefully install one RFID into the clean needle
- 3 scruff rat skin on the back of neck with fingers
- 4 Inset the sharp tips of the needle under the skin, target the back of the neck, close to the head.
- 5 Release the RFID by pushing the plunger of the syringe
- 6 Scan and make sure the RFID is under the skin
- 7 upload the scanned RFID into a spreadsheet