


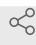

Jul 06, 2022

Isolation of ECs from Lymph node tissue for scRNAseq

Guido Krähenbühl¹¹MDC Berlin*In Development* Share

This protocol is published without a DOI.

AG Gerhardt



Guido Krähenbühl

ABSTRACT

Protocols describing the isolation of murine lymph node endothelial cells
Isolation of lymphatic endothelial cells for spectral flow cytometry and scRNA-seq
Basic analysis of spectral flow cytometry data of endothelial cells and scRNA-seq
Adapted from <https://doi.org/10.1016/j.xpro.2022.101267>

PROTOCOL CITATION

Guido Krähenbühl 2022. Isolation of ECs from Lymph node tissue for scRNAseq.
protocols.io
<https://protocols.io/view/isolation-of-ecs-from-lymph-node-tissue-for-scrnas-b98wr9xe>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

May 30, 2022

LAST MODIFIED

Jul 06, 2022

PROTOCOL INTEGER ID

63478

Preparation

- 1 Day before experiment check Antibodies to be used and prepare the following Reagents:

A	B	C
Dispase II	0.8 mg/mL	1 mL
Collagenase P	0.2 mg/mL	1 mL
DNase 1	0.1 mg/mL	1 ul
RPMI-1640		28 mL
Total		30 mL

Enzymatic Digestion solution for lymphnodes, store at 4C for 1 day (B is Final Concentration)

A	B	C
BSA	0.5%	2.5 g
EDTA, stock: 0.5M	2 mM	2 mL
PBS	-	498 mL
Total		500 mL

MACS buffer, store at 4C for up to 3 months (B is Final Concentration)


A	B	C
BSA	2%	10 mL
EDTA, stock: 0.5M	2 mM	2 mL
PBS	-	488 mL
Total		500 mL

FACS buffer, store at 4C for up to 3 months

- 2 Before starting prepare:
Dissection Tools
If Brain tissue is also taken, prepare perfusion kit
5mL RPMI-1640
Keep all reagents on ice

Tissue collection	1h
-------------------	----

- 3 Sacrifice mice using Isoflurine
(do not use cervical dislocation as this might affect perfusion of the brain)
- 4 Dissect inguinal, axillary, brachial and cervical lymphnodes

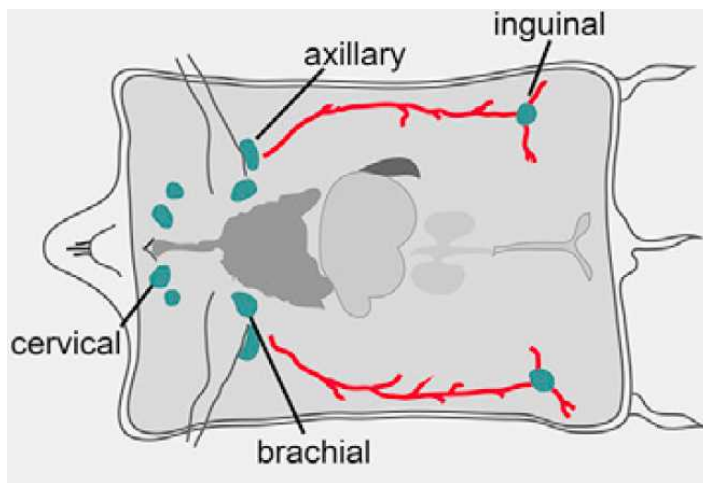

- 5 Perform transcardial perfusion via the left ventricle with **8 mL Ice-cold PBS** ice-cold at a perfusion rate of 2 mL/min for 5 min.
- 6 Dissect cervical lymph nodes and collect in **5 mL ice-cold RPMI-1640**
- 7 Collect remaining lymphnodes (axillary, brachial and inguinal) and add the to the cervical lymph nodes
- 8 **go to step #1 Repeat for a total of 5 mice to receive a sufficient number of cells (pool these samples)**

Tissue processing 35m

- 9 Remove the **RPMI-1640** and add **5 mL LN enzyme solution** at **37 °C** **00:15:00** ^{15m}
- 10 Aspirate and pipette the solution 10 times using a 5 mL serological pipette
- 11 Allow large fragments to settle down and transfer supernatant to 15m

20 mL ice-cold MACS buffer , add fresh 5 mL LN enzyme solution to remaining fragments and incubate at 37 °C 00:15:00
repeat 4 times or until fragments are fully digested
go to step #10

12 Centrifuge pooled supernatants at 300 x g, 4°C, 00:05:00

5m

13 Resuspend in 5 mL MACS buffer

14 Centrifuge at 300 x g, 4°C, 00:05:00

15 Resuspend in 500 µL FACS buffer

erythrocyte lysis

16 *

Lyse erythrocytes with lysis buffer
(can result in loss of sensitive endothelial cells)



FACS Sort 30m

17 Incubation with blocking antibody for 10 min at 4°C (CD16/CD32, 1:100 in FACS buffer)

18 Add Antibodies, mix and incubate at 4 °C 00:30:00 in the dark

30m

CD31 Antibody	<input type="text"/>
CD45 Antibody	<input type="text"/>
Antibody 3	<input type="text"/>
Antibody 4	<input type="text"/>

- 19 Wash cells twice with FACS buffer by centrifugation at  **300 rcf, 4°C, 00:05:00** and 5m
resuspend in  **300 µL FACS Buffer**
- 20 add viability dye immediately before analysis (7-AAD at 1:100, DAPI at 1:100)
- 21 Gating and sorting