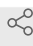




Sep 06, 2022

Amplifying the target genomic region by PCR

 In 1 collectionHanqin Li¹, Yogendra Verma¹, Dirk Hockemeyer¹, Frank Soldner²¹University of California, Berkeley; ²Albert Einstein College of Medicine1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.b4nvqve6

Devin E Snyder

ABSTRACT

This protocol describes a procedure for amplifying targeted genomic regions using a PCR reaction

DOI

dx.doi.org/10.17504/protocols.io.b4nvqve6

PROTOCOL CITATION

Hanqin Li, Yogendra Verma, Dirk Hockemeyer, Frank Soldner 2022. Amplifying the target genomic region by PCR. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b4nvqve6>



FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's
Grant ID: ASAP-000486

COLLECTIONS

**Genotyping by next generation sequencing**

KEYWORDS

ASAPCRN

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 03, 2022

LAST MODIFIED

Sep 06, 2022

PROTOCOL INTEGER ID

57781

PARENT PROTOCOLS

Part of collection

[Genotyping by next generation sequencing](#)

MATERIALS TEXT

Item	Vendor	Catalog #
dNTP	NEB	N0447L
10x HF buffer	NEB	B0518S
Phusion DNA polymerase	NEB	M0530S
Microseal® 'B' Adhesive Seals	Biorad	MSB1001
96-well PCR plate	GeneMate	T-3152-1

- 1 For a 150bp paired-end or single-read sequencing experiment, design primers in a way that the region of interest, like CRISPR cutting site or specific mutations, locates within 130bp from at least one of the primers
- 2 Add NGS adapters to 5' end of primers following the instructions of your local NGS facility
- 3 20 µl PCR reaction for each sample

3.1 PCR Reaction - setup

A	B
Ultrapure H ₂ O	10.6 µl
5x HF buffer	4 µl
2.5 mM dNTP	1.6 µl
10 µM primer Forward	0.5 µl
10 µM primer Reverse	0.5 µl
DMSO	0.6 µl
Titan DNA polymerase or Phusion	0.2 µl
Crude cell lysate	2 µl

- 4 Calculate the amount of each component needed for all samples, also count a negative control and positive control
- 5 In a pre-PCR area, mix all components except the crude cell lysate to a microcentrifuge tube or a 15 ml conical tube to make master mix
- 6 Aliquot 18 µl to each 200 µl microcentrifuge tube or each well in a 96-well PCR plate. Use reservoir and multi-channel pipet if there are many samples.
- 7 Add 2 µl crude cell lysate to the reaction. In one reaction, use 2 µl H₂O instead as negative control. In another reaction, use 2 µl previously validated crude cell lysate as positive control.
- 8 Cap the tubes or seal the plates properly with an adhesive seal.
- 9 Shake the tubes or plates vigorously to mix
- 10 Briefly spin the samples
- 11 Use the following touch-down PCR protocol in a thermocycler

11.1 Touch-down PCR protocol

A	B
98°C	3min
98°C	30s
70°C (touch down, 1°C/cycle)	30s
72°C	30s
Go to 2	12 cycles in total
98°C	30s
58°C	30s
72°C	30s
Go to 6	23 cycles in total
72°C	7min
4°C or 12°C	forever