



DEC 21, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.e6nvwjdjdlmk/v1

Protocol Citation: Robert Edwards, Shweta Jain 2023. Live Imaging of Primary Mouse Neuron Cultures. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.e6nvwjdjdlmk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


Protocol status: Working
 We use this protocol and it's working

Created: Aug 09, 2023

Live Imaging of Primary Mouse Neuron Cultures

Robert Edwards¹, Shweta Jain¹

¹University of California San Francisco

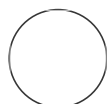








kelsey.barcomb

ABSTRACT

This protocol describes live imaging of primary neuron cultures. Included are methods for preparing hippocampal or dopamine neuronal cultures from neonatal mouse brain tissue. The imaging described involves labeling of dopamine neurons with a fluorescent DAT ligand and using virally encoded pHlourin sensors to measure vesicular release of neurotransmitter as the neurons are electrically stimulated. This technique was used in Jain et al., 2023 to compare vesicular release in neurons between various transgenic knockout mouse lines.

PROTOCOL MATERIALS

-  Basal Medium Eagle (BME) Thermo Fisher Catalog #21010046 Step 20
-  GDNF Recombinant Human Protein Thermo Fisher Catalog #PHC7045 Step 20
-  HyClone Fetal Bovine Serum Fisher Scientific Catalog #SH3039603 In 2 steps
-  D-()-Glucose solution Merck MilliporeSigma (Sigma-Aldrich) Catalog #G8769
- Step 8
-  B27 supplement without retinoic acid (50x) Gibco, ThermoFisher Catalog #17504044
- In 3 steps
-  Neurobasal™ Medium Gibco - Thermo Fisher Catalog #21103049 In 2 steps
-  Glutamax (100x) Gibco - Thermo Fischer Catalog #35050-061 In 2 steps

Last Modified: Dec 21,
2023


PROTOCOL integer ID:
86286

Keywords: ASAPCRN, Cell
Culture, Neuron, Mouse, Live
Imaging, Microscopy, Vesicle
Release

Funders
Acknowledgement:
ASAP-CRN
Grant ID: 020529

Preparation of Hippocampal Cultures

- 1 Coat coverslips with Poly-L-lysine
- 2 Collect mouse pups at postnatal day 0
 - 2.1 For Jain et al., 2023, cultures were prepared from WT, β 3A KO, β 3B KO, and *mocha* mice
- 3 Decapitate and remove brain, dissect out hippocampi from each hemisphere into Hank's balanced salt solution (HBSS) containing 10 mM HEPES and 20 mM glucose
- 4 Digest with papain (20 units/mL) in papain digestion buffer (HBSS with 20 mM CaCl_2 , 5 mM EDTA, 0.2 mg/mL L-Cysteine, 10 mM HEPES at pH 7.4) for 15 minutes
- 5 Wash hippocampi three times with HBSS containing 10 mM HEPES and 20 mM glucose

- 6 Triturate hippocampi to get a single cell suspension
- 7 Using a hemocytometer, count the density of cells in the suspension
- 8 Plate cells onto the poly-L-lysine-coated coverslip at a density of 350 cells/mm² in Minimal Essential Medium containing:
- 1X B27 supplement
⊗ B27 supplement without retinoic acid (50x) Gibco, ThermoFisher Catalog #17504044
 - 5% Fetal Bovine Serum (FBS)
⊗ HyClone Fetal Bovine Serum Fisher Scientific Catalog #SH3039603
 - 21 mM glucose ⊗ D-(+)-Glucose solution Merck MilliporeSigma (Sigma-Aldrich) Catalog #G8769
- 9 Store cultures in incubator at  37 °C
- 10 After day 1 *in vitro* (DIV1), replace 3/4 of the medium with Neurobasal medium
⊗ Neurobasal™ Medium Gibco - Thermo Fisher Catalog #21103049 containing:
- 1X B27 supplement
⊗ B27 supplement without retinoic acid (50x) Gibco, ThermoFisher Catalog #17504044
 - GlutaMAX ⊗ Glutamax (100x) Gibco - Thermo Fischer Catalog #35050-061
- 11 Infect cells at DIV4-5 using lentiviral vectors

11.1 For Jain et al., 2023, viruses encoding VGLUT2-pH, VMAT2-pH, or VGAT-pH were used

12 At DIV 6-7, inhibit glial proliferation by treating cultures with 4 μ M cytosine arabinoside (Ara-C)

Preparation of Dopamine Neuron Cultures

13 Coat coverslips with poly-L-lysine and laminin; plate astrocyte monolayer onto coverslips







14 Collect mouse pups at postnatal day 0

15 Decapitate and remove brain, dissect out midbrain (including ventral tegmental area and substantia nigra) from each hemisphere into Hank's balanced salt solution (HBSS) containing 10 mM HEPES and 20 mM glucose

16 Digest with papain (20 units/mL) in papain digestion buffer (HBSS with 20 mM CaCl_2 , 5 mM EDTA, 0.2 mg/mL L-Cysteine, 10 mM HEPES at pH 7.4) for 15 minutes


17 Wash digested tissue three times with HBSS containing 10 mM HEPES and 20 mM glucose

18 Triturate tissue to get a single cell suspension

- 19 Using a hemocytometer, count the density of cells in the suspension
- 20 Plate cells onto the pre-prepared coverslips at a density of 1000 cells/mm² in medium consisting of 60% Neurobasal Medium  Neurobasal™ Medium Gibco - Thermo Fisher Catalog #21103049, 30 % Basal Eagle Medium  Basal Medium Eagle (BME) Thermo Fisher Catalog #21010046, 10% FBS
-  HyClone Fetal Bovine Serum Fisher Scientific Catalog #SH3039603 containing:
- 1X B27 supplement
 -  B27 supplement without retinoic acid (50x) Gibco, ThermoFisher Catalog #17504044
 - 2 mM GlutaMAX  Glutamax (100x) Gibco - Thermo Fischer Catalog #35050-061
 - 10 ng/mL glial cell line-derived neurotrophic factor (GDNF)
 -  GDNF Recombinant Human Protein Thermo Fisher Catalog #PHC7045
 - 1X penicillin/streptomycin
- 21 At DIV 2-4, infect with lentivirus
- 21.1 For Jain et al., 2023, cells were infected with lentivirus encoding VMAT2-pH or VGLUT2-pH

Live Imaging

- 22 Maintain cultures until ready for imaging. Live imaging of cultures is performed at DIV14-16 for hippocampal cultures and DIV13-15 for dopamine neuron cultures

- 23** Prepare Tyrode's buffer:
- 119 mM NaCl
 - 25 mM HEPES
 - 2 mM CaCl_2
 - 2 mM MgCl_2
 - 2.5 mM KCl
 - 30 mM glucose
- Adjust pH to 7.4
- 24** For midbrain cultures, label dopamine neurons, incubate dopamine neuron cultures in Tyrode's buffer (119 mM NaCl, 25 mM HEPES, 2 mM CaCl_2 , 2 mM MgCl_2 , 2.5 mM KCl, and 30 mM
- 24.1** Make Tyrode's +JHC1-64 by adding the fluorescent DAT ligand JHC1-64 at 30 nM to Tyrode's buffer
- 24.2** Incubate dopamine neuron culture in Tyrode's +JHC1-64 for 5 minutes at  Room temperature
- 24.3** Rinse neurons in Tyrode's buffer (no JHC 1-64)
- 25** Add fresh Tyrode's buffer to culture dishes, mount coverslips in a perfusion and stimulation chamber of a TE300 inverted epifluorescence microscope
- 26** Monitor fluorescence:
- For pHlourin, use 470/40 nm excitation and 525/70 nm emission bandpass filters
 - For JHC1-64, use 545/25 nm excitation and 605/70 nm emission bandpass filters
- Use the red channel to identify dopamine neurons then switch to the pHlourin channel for experiments

- 27 Place platinum-iridium stimulating electrodes near the dopamine neuron being observed
- 28 Acquire images, typically at 1 Hz except for Fluo-5F image where images are acquired at 10 Hz
- 29 While acquiring images, evoke action potentials using 1 ms bipolar current pulses through the stimulating electrode at 5-10 V/cm