



AUG 03, 2023

OPEN  ACCESS**DOI:**

dx.doi.org/10.17504/protocols.io.6qpvr45wpgmk/v1

Protocol Citation: Oyedele Adeyi, Andrew Nelson, Beth Thompson, Laura Niedernhofer 2023. Human FFPE Liver Pathology Assessment -- University of Minnesota Human TMC.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.6qpvr45wpgmk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


Protocol status: Working
We use this protocol and it's working

Human FFPE Liver Pathology Assessment -- University of Minnesota Human TMC

Oyedele Adeyi¹, Andrew Nelson¹, Beth Thompson¹,
Laura Niedernhofer¹

¹University of Minnesota Medical School, Minneapolis, MN USA

Cellular Senescence Network (SenNet) Method Development
Community

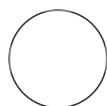

Mickayla DuFresne-To

University of Minnesota - Twin Cities

ABSTRACT

Purpose: Describe the parameters evaluated per FFPE human liver sample received/obtained for the University of Minnesota SenNet TMC. These parameters are used to evaluate whether the tissue can be characterized as "normal" vs "diseased" prior to further protocols.

Created: Apr 25, 2023

Last Modified: Aug 03, 2023

PROTOCOL integer ID:
81034

Keywords: Pathology, Liver, Human, Sample evaluation, UMN, UMinnesota, University of Minnesota, Niedernhofer, Tissue, FFPE

- 1 Prepare two slides with serial sections of the FFPE liver sample – 1 section per slide.
- 2 Section A -- Stain with H&E.
Section B -- Stain with trichrome.
- 3 Mount coverslip.
- 4 Schedule evaluation with pathologist (ideally with a strong liver background).
- 5 Request assessment and record of the following scores:
 1. Steatosis, macro %
 2. Steatosis, score (0, 1, 2, 3)
 3. Ballooning (0, 1, 2)
 4. Inflammation (0, 1, 2, 3)
 5. Fibrosis (0, 1a, 1b, 1c, 2, 3, 4)
 6. Any comments on preservation/quality of the sample
 7. Other/Additional identifications (e.g. the presence of NASH or other disease states)
- 6 Record these values and consider them prior to moving the sample forward for experimentation processing.

- 7 Inventory and archive H&E and trichrome stains at room temperature – Can be scanned if future images are required.

Inventory and archive FFPE blocks at room temperature.