




JUL 25, 2023

🌐 Immunostaining of iPSC-derived neurons for quantification of synaptic proteins

Dan
Dou¹, Erika Holzbaur¹

¹University of Pennsylvania

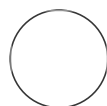

Dan Dou

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.8epv5x91ng1b/v1

Protocol Citation: Dan Dou, Erika Holzbaur 2023. Immunostaining of iPSC-derived neurons for quantification of synaptic proteins. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.8epv5x91ng1b/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jul 20, 2023

Last Modified: Jul 25, 2023

PROTOCOL integer ID:
85305

Keywords: ASAPCRN

ABSTRACT

Here, we fix, permeabilize, and stain human iPSC-derived neurons for the purpose of observing and quantifying somal proteins of interest. For preceding culture of neurons, see "Protocol: Culture and transfection of iPSC-derived neurons for live-imaging of axonal cargoes."

MATERIALS

35 mm imaging dish, 20 mm glass diameter (Mattek, p35g-1-5-20-c)
18mm circular coverglass (Electron Microscopy Sciences, Cat# 72222-01)
Prolong Gold Antifade Mountant (Thermo Fisher, P36930)

1 At DIV14, fix human iNeurons in 4% paraformaldehyde supplemented with 4% sucrose for 15

minutes at 37 degrees C

- 2 Wash four times with PBS
- 3 Permeabilize for 15 minutes in 0.2% Triton-X in PBS
- 4 Block for 1 hour with 5% goat serum and 1% BSA in PBS
- 5 Incubate in primary antibody diluted in blocking solution at room temperature for 1 hour
- 6 Wash three times with PBS
- 7 Incubate in secondary antibody diluted in blocking solution for 1 hour at room temperature
- 8 Wash three times with PBS
- 9 Remove PBS and add 40 μ L Prolong Gold (Thermo Fisher). Using forceps, apply coverglass.

