




JUL 25, 2023

Co-immunoprecipitation using GFP-trap

Dan
Dou¹, Erika Holzbaur¹

¹University of Pennsylvania

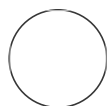

Dan Dou

ABSTRACT

Here, we describe performing co-immunoprecipitation experiments in HEK293 cells using GFP-trap beads

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.6qpvr36o2vmk/v1

Protocol Citation: Dan Dou, Erika Holzbaur 2023. Co-immunoprecipitation using GFP-trap. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.6qpvr36o2vmk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jul 20, 2023

Last Modified: Jul 25, 2023

PROTOCOL integer ID:
85310

Keywords: ASAPCRN

- 1 See "Protocol: HEK293 cell culture for co-immunoprecipitation experiments" for preceding culture and transfection. For the present study, HEK cells were transfected with EGFP-tagged

bait proteins. Prey proteins were either endogenous or transfected with HA- or SNAP-tags.

2 Prepare wash buffer:

A	B
Tris-HCl pH 7.5	10 mM
NaCl	150 mM
EDTA	0.5 mM
Triton-X (optional)	0.4%

Note

Triton-X should be included for prey proteins where appreciable non-specific binding is observed in the EGFP vector condition.

3 Prepare lysis buffer. Lysis buffer composition differed for experiments needing lambda phosphatase treatment

3.1 Lysis buffer for non-lambda phosphatase experiment:

A	B
Tris-HCl pH 7.5	10 mM
NaCl	150 mM
EDTA	0.5 mM
NP-40	0.5%
PMSF	1 mM
TAME	0.01 mg/mL
Leupeptin	0.01 mg/mL
Pepstatin A	0.001 mg/mL

3.2 Lysis buffer compatible with lambda phosphatase experiment:

A	B
1x NEBuffer for Protein MettaloPhosphatases (New England BioLabs)	1x
NP-40	0.5%
Leupeptin	0.01 mg/mL
ddH ₂ O	To desired volume

4 24 hours after transfection, wash HEK cells twice in ice-cold PBS and lyse in appropriate lysis buffer. Use 600 µL lysis buffer per experimental condition (pooled across the three 10cm dishes).

5 Clarify lysates at 10 x g at 4 degrees C for 10 minutes

6 Wash 25 µL GFP-trap beads per experimental condition in 500 µL wash buffer, in low protein binding Eppendorf tubes under rotating agitation.

Note

For large protein complexes, GFP-Trap Magnetic Particles M-270 should be used instead of GFP-Trap Magnetic Agarose beads

7 Equilibrate beads in lysis buffer for 5 min at 4 degrees C under rotating agitation

8 For lambda phosphatase experiments: add 60 µL of 10 mM MnCl₂ and 24 µL lambda phosphatase (2,400 units; New England BioLabs) for a final reaction volume of 600 µL per experimental condition.

Note

For conditions that are lambda phosphatase-negative, add 24 μL ddH₂O instead

- 8.1** Incubate at 30 degrees C for 30 minutes
- 9** Incubate beads with clarified lysate (or lambda phosphatase-treated lysate) for 1 hour at 4 degrees C under rotating agitation
- 10** Wash beads three times for 5 min in wash buffer at 4 degrees C under rotating agitation
- 11** Resuspend beads in 60 μL denaturing buffer, and boil for 10 minutes to release bound protein