




JUN 21, 2023

🌐 LB Agar

sgoodwin¹

¹Sheffield iGEM

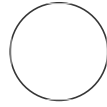

sgoodwin

ABSTRACT

Protocol for producing LB agar growth medium

OPEN ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.dm6gp38xjvzp/v1

Protocol Citation: sgoodwin 2023. LB Agar. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.dm6gp38xjvzp/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: In development

We are still developing and optimizing this protocol

Created: Jun 21, 2023

Last Modified: Jun 21, 2023

PROTOCOL integer ID: 83813

LB Agar

3h 11m

1

Add  1.5 g LB powder,  1.5 g agarose powder to conical flask

5m

2

Dilute with  200 mL distilled water. Ensure fully dissolved

5m

3

Seal conical flask with sponge bung and foil

1m

4

Sterilise growth medium by autoclaving

3h