



MAY 14, 2023

ATAC-seq from nuclei from frozen tissue

Eric RA Pederson¹

¹Uppsala University

methods



Eric RA Pederson
Uppsala University

ABSTRACT

Modified ATAC-seq method for frozen tissue, in this case brain tissue.

GUIDELINES

The Tn5 enzyme used in this experiment was from a locally produced batch from the Protein Science Facility at the Karolinska Institutet in Stockholm, which is first discussed in Picelli and others (2014). Thus, if the Tn5 is purchased through a company it may react differently.

MATERIALS

- NEBNext High-Fidelity 2X PCR Master Mix - 250 rxns New England Biolabs Catalog #M0541L
- Digitonin, 40ul Promega Catalog #G9441
- Tn5 transposase with Nextera adapters loaded Contributed by users Catalog #UC-Macro-Tn5-Nextera adapter
- Zymo DNA Clean & Concentrator - 5 Zymo Research Catalog #D4014
- MACS SmartStrainers 30um Miltenyi Biotec Catalog #130-098-458
- Halt™ Protease Inhibitor Cocktail, EDTA-Free (100X) Thermo Fisher Catalog #78437
- SYBR™ Green I Nucleic Acid Gel Stain - 10,000X concentrate in DMSO Thermo Fisher Catalog #S7563
- Protease Inhibitor Tablets cOmplete Mini EDTA free Roche Catalog #11836170001

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.bp2l6b4nkgqe/v1

Protocol Citation: Eric RA Pederson 2023. ATAC-seq from nuclei from frozen tissue. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bp2l6b4nkgqe/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Sep 14, 2020

Last Modified: May 14, 2023

PROTOCOL integer ID:
42096

Keywords: ATAC-seq, frozen tissue



KIMBLE 2mL Glass Dounce Tissue Grinder Set Merck MilliporeSigma (Sigma-Aldrich) Catalog #D8938



2 ml LoBind Tubes Eppendorf Catalog
#0030108078



1.5 mL LoBind tubes Eppendorf Catalog
#022431021



Falcon® 100 mm TC-treated Cell Culture Dish Corning Catalog
#353003

Equipment

DNA LoBind Tubes

NAME

Microcentrifuge tubes

TYPE

Eppendorf

BRAND

0030108051

SKU

<https://online-shop.eppendorf.com/OC-en/Laboratory-Consumables-44512/Tubes-44515/DNA-LoBind-Tubes-PF-56252.html>

LINK



Qubit 2.0 Fluorometer Contributed by users



Qubit® dsDNA HS Assay Kit Thermo Fisher Scientific Catalog #Q32854

Equipment

Countess 3 FL Automated Cell Counter

NAME

Automated Cell Counter

TYPE

Thermofisher scientific

BRAND

AMQAF2000

SKU

<https://www.thermofisher.com/th/en/home/life-science/cell-analysis/cell-analysis-instruments/automated-cell-counters/models/countess-3-fl.html>

LINK

Equipment	
Thermomixer C	NAME
Eppendorf	BRAND
2231000667	SKU
https://www.pipette.com/2231000667-Promotion-Eppendorf-ThermoMixer-C-with-24x1-5-mL-SmartBlock-and-ThermoTop	LINK

Equipment	
Dounce, 2ml	NAME
Homogenizer	TYPE
KIMBLE	BRAND
432-0250	SKU
https://se.vwr.com/store/product/561196/homogenisator-dounce-kimble	LINK
2 ml	SPECIFICATIONS

Equipment	
Bio RS-24 Mini-rotator	NAME
mini-rotator	TYPE
BioSan	BRAND
RS-24	SKU
https://biosan.lv/products/-bio-rs-24-mini-rotator-for-test-tubes-with-timer/	LINK

Equipment	
4200 TapeStation System	NAME
Electrophoresis tool for DNA and RNA sample quality control.	TYPE
TapeStation Instruments	BRAND
G2991AA	SKU
https://www.agilent.com/en/product/automated-electrophoresis/tapestation-systems/tapestation-instruments/4200-tapestation-system-228263	LINK

- ⊗ 1 M Calcium Chloride (CaCl₂) Fisher Scientific Catalog #BP510
- ⊗ Mg(Ac)₂ Contributed by users
- ⊗ Tris-HCl pH 7.5 Contributed by users
- ⊗ Sodium Chloride Fisher Scientific Catalog #S271
- ⊗ Magnesium Chloride Fisher Scientific Catalog #AC223210010
- ⊗ Molecular grade water nuclease-free Contributed by users
- ⊗ NN-Dimethylformamide (DMF) solution Merck MilliporeSigma (Sigma-Aldrich) Catalog #D4551
- ⊗ Triton X-100 Merck MilliporeSigma (Sigma-Aldrich) Catalog #T8787-50ML
- ⊗ DTT Merck MilliporeSigma (Sigma-Aldrich) Catalog #D0632
- ⊗ Glycerol, 1000ml Promega Catalog #H5433
- ⊗ EDTA (0.5 M), pH 8.0 Life Technologies Catalog #AM9260G
- ⊗ 2-Mercaptoethanol Merck MilliporeSigma (Sigma-Aldrich)
- ⊗ High Sensitivity D1000 Reagents Agilent Technologies Catalog #5067-5585
- ⊗ High Sensitivity D1000 ScreenTape Agilent Technologies Catalog #5067-5584

—

Amount	Reagents
3 ml	1 M CaCl ₂
0.6 ml	3 M Mg(Ac) ₂
6 ml	1 M Tris pH 7.9
89.2 ml	molecular grade water

6x Homogenization Buffer Stable Solution.

Buffer recipes;

A	B
Amount	Reagents
2.27084 ml	6x Homogenization Buffer stable
3.78 ml	100mM PMSF
0.28 ul	14.3 M B-mercaptoethanol
1/2 tablet	Protease inhibitor (cOmplete Mini)

6x Homogenization Buffer Unstable Solution.

Amount	Reagents
500ul	1M Tris-Hcl ph 7.5
100ul	5M NaCl
150ul	1M MgCl ₂
49.25ml	H ₂ O

1x Homogenization Buffer Unstable Solution

Amount	Reagents
500ul	1M Tris-Hcl ph 7.5
100ul	5M NaCl
150ul	1M MgCl ₂
49.25ml	H ₂ O

ATAC-RSB

Amount	Reagents
8 ml	ATAC-RSB
8 ul	10% Tween

ATAC-RSB + 10% Tween

2X TD buffer:

Amount	Reagents
2 ml	1 M Tris-Hcl pH 7.5
1 ml	1 M MgCl ₂
20 ml	100% Dimethyl Formamide

2X TD buffer

(before the addition of dimethyl formamide, adjust the pH to 7.6 with 100% acetic acid)

DF buffer

Amount	Reagents
100 mM	HEPES (pH 7.2)
200 mM	NaCl
0.2 mM	EDTA
2 mM	DTT
2.00%	Triton X-100
20.00%	Glycerol

2X DF buffer

Primer	Sequence	Amount
Tn5-A primer	TCGTCGGCAGCGTCAGATGTGTATAAGAGACAG	16 ul
Tn5-rev primer	CTGTCTCTTATACACATCT	16 ul

Tube A

Primer	Sequence	Amount
Tn5-B primer	GTCTCGTGGGCTCGGAGATGTGTATAAGAGACAG	16 ul
Tn5-rev primer	CTGTCTCTTATACACATCT	16 ul

Tube B

SAFETY WARNINGS

- Regular lab safety rules apply. Ensure the use of goggles during the usage of dry ice.

BEFORE START INSTRUCTIONS

The pre-experimentation steps are important
Also it is important that all the primers have been ordered and reconstituted beforehand.

Pre-experimentation

- 1 All steps should be performed on ice or at 4 °C .

1m

2 Pre-chill all Dounces and pestles to 4 °C in a fridge or on ice

10m

3 Pre-chill all tubes.

10m

For each sample you are processing, you will need:

(i) One 2 mL 2 ml LoBind Tubes Eppendorf Catalog #0030108078 per sample

(ii) Three 1.5 mL LoBind tubes Eppendorf Catalog #022431021 per sample

(iii) one PCR tube per sample

(iv) One 50 ml Falcon tube Contributed by users for filtration step per sample

4 Prepare buffers.

1h 30m

- i) 6x Homogenization Buffer Stable Solution.
- ii) 6x Homogenization Buffer Unstable Solution.
- iii) 1x Homogenization Buffer Unstable Solution.
- iv) ATAC-RSB
- v) ATAC-RSB + 10% Tween
- vi) 2X TD buffer
- vii) 2X DF buffer

4.1 i) 6x Homogenization Buffer stable Solution.

20m

A	B
Amount	Reagents
3 ml	1 M CaCl ₂
0.6 ml	3 M Mg(Ac) ₂
6 ml	1 M Tris pH 7.9
89.2 ml	molecular grade water


6x Homogenization Buffer Stable

4.2 ii) 6x Homogenization Buffer Unstable Solution.

20m

Amount	Reagents
2.27084 ml	6x Homogenization Buffer stable
3.78 ml	100mM PMSF
0.28 ul	14.3 M B-mercaptoethanol
1/2 tablet	Protease inhibitor (cOmplete Mini)

6x Homogenization Buffer Unstable Solution

 cOmplete™, Mini, EDTA-free (Protease Inhibitor) Roche Catalog
##11836170001)

4.3 ii) 1x Homogenization Buffer Unstable Solution.

20m

Amount	Reagents
2.166645 ml	6x Homogenization Buffer Unstable Solution
4.16 ml	1 M sucrose
2.6 ul	500 mM EDTA
130 ul	10.00% NP10
6.540755 ml	H2O

1x Homogenization Buffer Unstable Solution

4.4 iii) ATAC-RSB:

20m

A	B
Amount	Reagents
500ul	1M Tris-Hcl ph 7.5
100ul	5M NaCl
150ul	1M MgCl2
49.25ml	H2O

ATAC-RSB

4.5 iv) ATAC-RSB + 10% Tween

20m

Amount	Reagents
--------	----------

8 ml	ATAC-RSB
8 ul	10% Tween

ATAC-RSB + 10% Tween

4.6 v) 2X TD buffer

20m

A	B
Amount	Reagents
2 ml	1 M Tris-Hcl pH 7.5
1 ml	1 M MgCl ₂
20 ml	100% Dimethyl Formamide

2X TD buffer

(before the addition of dimethyl formamide, adjust the pH to 7.6 with 100% acetic acid)

4.7 vi) 2X DF buffer

20m

A	B
Amount	Reagents
100 mM	HEPES (pH 7.2)
200 mM	NaCl
0.2 mM	EDTA
2 mM	DTT
2.00%	Triton X-100
20.00%	Glycerol

2X DF buffer

5 Tn5 assembly reaction

2h

5.1

1h

A	B	C
Primer	Sequence	Amount
Tn5-A primer	TCGTCGGCAGCGTCAGATGTGTATAAGAGACAG	16 ul
Tn5-rev primer	CTGTCTCTTATACACATCT	16 ul

Tube A

A	B	C
Primer	Sequence	Amount
Tn5-B primer	GTCTCGTGGGCTCGGAGATGTGTATAAGAGACAG	16 ul
Tn5-rev primer	CTGTCTCTTATACACATCT	16 ul

Tube B

🔥 95 °C to 🔥 25 °C
cooling -0.1°C/second

~1 hour on PCR machine
Olgio Annealing program

5.2

TN5 Assembly

1h

A	B
Amount	Reagent
25 ul	Tn5
15.5 ul	Tube A
15.5 ul	Tube B

A	B
33 ul	2X DF buffer

Tn5 Assembly Mix



Room temperature



01:00:00

Tn5 from the Protein Science Facility at the Karolinska Institutet in Stockholm (Picelli *et al.*, 2014)

or



Tn5 transposase with Nextera adapters loaded Contributed by users Catalog #UC-Macro-Tn5-Nextera adapter

Nuclei extraction and filtration

- 6 Materials for the cutting of tissue on dry ice. 1m

- 6.1 1m
 - Gloves
 - White warm gloves
 - Dry ice
 - Blades and handle
 - Forceps
 - Ethanol spray
 - Cell culture dish


- 7 In the most sterile way possible, cut a small piece of tissue, half a pea size or so, and leave it in the petri dish with a marking on the lid, in the dry ice. Weigh it and cut again if needed. 20m

Make sure to use the ethanol to clean everything and be careful not to cut yourself.

- 8 Add 2 mL 1X HB buffer into the dounce, which is sitting in the ice. 2m

Add 0.2 µL 1 M DTT and and 1 µL 100X Halt protease inhibitor

- 9 Place 20 mg frozen tissue into a pre-chilled 2 ml Dounce containing 1 ml cold 1x HB and let thaw 5m

for  00:05:00 .

10 Dounce with “A” loose pestle until resistance goes away (~10 strokes).

2m

Put the A pestle into the beaker of water


11 Dounce with “B” tight pestle for 20 strokes.

2m

Put the B pestle into the beaker of water


12 Pour everything from the dounce into a 30 um MACS smartstrainer which is sitting on top of a labelled 50 ml falcon tube sitting in ice.

2m


 MACS SmartStrainers 30um Miltenyi Biotec Catalog #130-098-458

13 Let it drip through for 10-15 minutes.

15m

14 Transfer to a labelled  2 mL Lobind Eppendorf tube, already cold from sitting in ice.

5m

15 To pellet the nuclei, centrifuge  900 rpm, 4°C, 00:10:00

10m






15.1 Transfer the supernatant to a new tube without disturbing the pellet

10m

Repeat the centrifugation

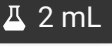
16 Discard supernatant

5m


- 17** Resuspend the nuclei in  200 µL of ATAC-RSB + 10% Tween 5m
- 18** Count the amount of cells using the cell counter in the cell lab. 20m
- 18.1** Put  10 µL Trypton blue onto a piece of parafilm and mix with  10 µL of the sample
- 18.2** Take  10 µL of the mix and pipette into a cell counting cell
- 18.3** Measure on the cell counting machine, making sure to adjust for the smaller size of the nuclei.
- The machine will think they are dead cells.
- 19** Calculate the amount of nuclei to use for the next step;
(500000 nuclei) 5m
- 19.1** ATAC-seq requires 50,000 cells for the experiment. So for example if the number of nuclei is;
- 4.76×10^7
- then the calculation will look like this;
- $50000 / 47600 = 1.10$
- 20** Turn on the thermomixer to  37 °C 1m

ATAC-seq

27m


- 21** Add the calculated amount of cells into a  2 mL lobind eppendorf tube with 1 ml ATAC-RSB + 10% Tween.

5m

So taking the example above one would take 1.10 ul of nuclei in the  200 μ L ATAC-RSB + 10% Tween to the new tube.

- 22** Centrifuge nuclei for  00:10:00  900 rcf, 4°C, 00:10:00

10m

- 22.1** Here one can also save the nuclei for later. Simply spin the tube down with the new tube of diluted nuclei, discard the supernatant, add  300 μ L of the and put in the -20° freezer.

10m

 00:10:00  900 rcf, 4°C, 00:10:00

- 23** Discard most of the supernatant. Leave a bit in the bottom to ensure that there is enough for the reaction and that the pellet stays intact

2m



- 24** Add the reaction mix to the pellet and remaining water and resuspend the pellet in the mix

5m



A	B
Amount	Reagents
25 ul	2X TD buffer
16.5 ul	1X PBS (cold)
0.5 ul	1% digitonin
0.5 ul	10% Tween -20
2.5 ul	TN5 assembly

Reaction Mix

- 25** Put the tube into the thermomixer at  1000 rpm, 37°C, 00:30:00 30m
- 26** Take the tubes out and immediately proceed to the column clean up 2m
- 27** Use the Zymo DNA clean & concentrator to clean the nuclei 30m
 Zymo DNA Clean & Concentrator - 5 Zymo Research Catalog #D4014
- 27.1** Add 2-7 volumes of the DNA binding buffer to each volume of DNA sample
So in this case the volume is approximately  50 µL , so add  300 µL of DNA binding buffer
- 27.2** Transfer mixture to a provided Zymo-spin column in a collection tube
- 27.3** Centrifuge  10000 x g, 00:00:30
Discard supernatant
- 27.4** Add  200 µL DNA wash buffer to the column
Centrifuge  10000 x g, 00:00:30
- 27.5** repeat the wash step

27.6 Transfer the column to a 1.5 low-bind eppendorf tube.

Add  21 µL 2 DNA Elution buffer to the column.

Centrifuge  10000 x g,
00:00:30

28 Can stop here and start the next day if necessary.
Leave samples in 4°

29 Set up PCR;

20m

A	B
Amount	Reagent
2.5 ul	Primer AD1
2.5 ul	Primer AD2.#
25 ul	NEBNext Master Mix
20 ul	sample

PCR



NEBNext High-Fidelity 2X PCR Master Mix - 250 rxnsNew England Biolabs Catalog
#M0541L

29.1 PCR program - ATAC-seq pre-amplification

40m

A	B	C
Temperature	Time	Cycle
72°	5 minutes	1
98°	30 seconds	-

A	B	C
98°	10 seconds	5
63°	30 seconds	-
72°	1 minute	-
72	5 minutes	1
4	infinite	-

ATAC-seq pre-amplification

- 30** Remove PCR tubes from the machine and put immediately onto ice.

1m

- 31** Proceed immediately to the qPCR amplification to determine additional cycles step

1m


- 32** qPCR amplification to determine additional cycles


1h 30m

prepare a mix or master mix

A	B	C
Reagent	1X	6.5X
Molecular grade water	3.76 ul	24.5 ul
Primer AD1	0.5 ul	6.5 ul
Primer AD2.#	0.5 ul	0.5 ul *
25x SYBR green (in DMSO)	0.24 ul	1.56 ul
2x NEBNext Master Mix	5 ul	32.5 ul

qPCR mix

 SYBR™ Green I Nucleic Acid Gel Stain - 10,000X concentrate in DMSO Thermo Fisher Catalog #S7563

 NEBNext High-Fidelity 2X PCR Master Mix - 250 rxns New England Biolabs Catalog #M0541L

32.1 qPCR program

A	B	C
Temperature	Time	Cycle
98°	30 seconds	1
98°	10 seconds	20
63°	30 seconds	-
72°	1 minute	-

qPCR program

32.2 qPCR extra setup options

1m

- 6 machine not 6
- standard
- SYBR green
- no melt curve
- turn off ROX
- Fill in sample list

33 Determine the required amount of cycles that each sample needs in addition

5m

33.1 Looking at the final amplification curve, determine the max fluorescence where the graph plateaus

33.2 Determine 1/3 of that number.

For example if the max is 1400000 then 1/3 would be 433333

33.3 Check the graph for how many cycles line up with this number from the curve.

In the case above it would be 6 because of where the line was

34 Continue the PCR with the additional amount of cycles skipping the beginning parts of the program

34.1 # of additional cycles

40m

A	B	C
Tempurature	Time	Cycle
98°	10 seconds	Based on calculation
63°	30 seconds	-
72°	1 minute	-
4°	Infinite	1

Additional cycles


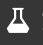




35 Take the tubes out and immediately proceed to the column clean up

1m

36 Use the Zymo DNA Clean & Concentrator to clean the nuclei

30m



- 36.1** Add 2-7 volumes of the DNA binding buffer to each volume of DNA sample 1m
- So in this case the volume is approximately 50 ul, so add 300 ul of DNA binding buffer
- 36.2** Transfer mixture to a provided Zymo-spin column in a collection tube 1m
- 36.3** Centrifuge  10000 x g, 00:00:30 1m
- Discard supernatant
- 36.4** Add  200 μ L DNA wash buffer to the column 2m
- Centrifuge  10000 x g, 00:00:30
- 36.5** Repeat the wash step 2m
- 36.6** Transfer the column to a 1.5 low-bind Eppendorf tube. 1m
-  Add 21 ul **Sterile water** buffer to the column.
-  Centrifuge  10000 x g, 00:00:30
- 37** Quality control 1h
- TapeStation (DS1000) and Qubit quantification