



DEC 14, 2022

WORKS FOR ME

1



nCoV-2019 Illumina Miniseq sequencing protocol (2,000bp amplicon)

Forked from [nCoV-2019 sequencing protocol \(RAPID barcoding, 1200bp amplicon\)](#)

DOI

dx.doi.org/10.17504/protocols.io.6qpvrpdjbgmk/v1Bruno Gomez-Gil¹, juli.enciso²¹CIAD AC Mazatlan; ²CIAD

juli.enciso: Julissa Enciso Ibarra;



Bruno Gomez-Gil

CIAD

COMMENTS 1

DISCLAIMER

It work for us, any modification is up to you.

ABSTRACT

This is a fork of the protocol <https://dx.doi.org/10.17504/protocols.io.bh7hj9j6> but modified for tiled 2000bp amplicons, tagmentation with Nextera XT, indexing, and sequencing with the Illumina Miniseq platform.

It has already produced very good sequences.

Much of this protocol is base on this paper: <https://doi.org/10.1093/biomethods/bpaa014>

DOI

dx.doi.org/10.17504/protocols.io.6qpvrpdjbgmk/v1

PROTOCOL CITATION

Bruno Gomez-Gil, juli.enciso 2022. nCoV-2019 Illumina Miniseq sequencing protocol (2,000bp amplicon). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.6qpvrpdjbgmk/v1>



FUNDERS ACKNOWLEDGEMENT

CONACYT

Grant ID: 321122

FORK NOTE

FORK FROM

[Forked from nCoV-2019 sequencing protocol \(RAPID barcoding, 1200bp amplicon\), Nikki Freed](#)

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Feb 28, 2021

LAST MODIFIED

Dec 14, 2022

PROTOCOL INTEGER ID

47758

GUIDELINES

Tested with high viral copy numbers (<30 Ct).

Citation: Bruno Gomez-Gil, juli.enciso nCoV-2019 Illumina Miniseq sequencing protocol (2,000bp amplicon)

<https://dx.doi.org/10.17504/protocols.io.6qpvrpdjbgmk/v1>

This is an open access protocol distributed under the terms of the **Creative Commons Attribution License** (<https://creativecommons.org/licenses/by/4.0/>), which permits unrestricted use, distribution, and reproduction in any medium,

MATERIALS TEXT

cDNA

- [GoScript™ Reverse Transcriptase Kit](#) by Promega.

Multiplex DNA

-

Library preparation

- [Nextera XT](#) DNA Library Preparation Kit
-

SAFETY WARNINGS

Please follow standard health and safety guidelines when working with COVID-19 patient samples.

DISCLAIMER

It work for us, any modification is up to you.

Sample preparation

10m

- 1 Dilute the sample depending on the Ct values, this will reduce the likelihood of PCR-inhibition.

Ct range	Sample	Water
----------	--------	-------

12-15	2 µL	198 µL
16-18	2 µL	18 µL
>18	no dilution	

cDNA preparation

2h

- 2 We use the GoScript™ Reverse Transcriptase Kit **Promega Catalog #A5001**
Mix (pipetting) the following components in Eppendorf tube.

5m

Component	Volume
-----------	--------

Nuclease-free water	4 µL
GoScript™ Reaction Buffer, Random Primer	2 µL
GoScript™ Enzyme Mix	0.4 µL
Final volume	10 µL

Note

Note

- 3 Prepare the **mastermix** on ice, mix by pipetting.

Component	Volume
-----------	--------

GoScript™ Reverse Transcription Mix 10 µL
 RNA 3 µL
Final volume 13 µL

1h 21m

4 Incubate the reaction as follows:

Anneal primers 25 °C for 00:05:00
 Extension 42 °C for 01:00:00
 Inactivation 70 °C for 00:15:00

Snap cool in a prechilled metal rack or on ice 00:01:00

Note

Note

Primer pool preparation

5 **PRIMERS** for this protocol are described in protocol [nCoV-2019 sequencing protocol \(RAPID barcoding, 1200bp amplicon\)](#) and are the 2,000 bp option. We selected for 2,000 bp amplicons so they could be more easily tagged in the library preparation. Two pools are prepared, Pool1 with 18 primers and Pool2 with 16 primers.

POOL1. Average Tm 61.01 °C

A	B	C	D	E	F	G	H	I	J	K
Name	Sequence	Direction	Start	End	Length	Product Siz	Tm	%GC	Hairpin Tm	Pa
SARSCoV_2000_01_LEFT	ACCAACCAACTTTTCGATCTCTTGT	forward	31	54	24	2049	60.7	41.7	None	Nc
SARSCoV_2000_01_RIGHT	ACACCACCTGTAATGTAGGCCA	reverse	2058	2079	22	2049	61.4	50	39	Nc
SARSCoV_2000_03_LEFT	TCGCACAAATGTCTACTTAGCTGT	forward	3772	3795	24	1814	60.6	41.7	None	0.7
SARSCoV_2000_03_RIGHT	GTGTGCCCATGTACATAACAGCT	reverse	5563	5585	23	1814	61.2	47.8	42.7	0.7
SARSCoV_2000_05_LEFT	CAATCATGCAATTGTTTTTCAGCTATTT	forward	7299	7328	30	1825	60.4	30	33	Nc
SARSCoV_2000_05_RIGHT	CGTGTGTCAGGGCGTAACTTT	reverse	9102	9123	22	1825	61.6	50	None	Nc
SARSCoV_2000_07_LEFT	GGACGTACCATATTGGGTAGTGC	forward	10886	10908	23	1834	60.8	52.2	43.8	Nc
SARSCoV_2000_07_RIGHT	TCTGTCGTAGTGCAACAGGACT	reverse	12698	12719	22	1834	61.3	50	41.4	Nc
SARSCoV_2000_09_LEFT	ACCACTTCAGAGAGCTAGGTGT	forward	14477	14498	22	1925	60.5	50	None	Nc
SARSCoV_2000_09_RIGHT	ACAACCTGGAGCATTGCAACA	reverse	16380	16401	22	1925	61.5	45.5	None	Nc
SARSCoV_2000_11_LEFT	TGGCATACCTAAGGACATGACCT	forward	18168	18190	23	1814	60.9	47.8	38.3	Nc
SARSCoV_2000_11_RIGHT	CAGTGAGTGGTGACAAAATCGT	reverse	19960	19981	22	1814	61.6	50	38.7	Nc
SARSCoV_2000_13_LEFT	TCCTCAGTTTTACATTCAACTCAGGA	forward	21695	21720	26	1937	60.2	38.5	45.2	Nc
SARSCoV_2000_13_RIGHT	TGACTAGCTACACTACGTGCCC	reverse	23610	23631	22	1937	61.5	54.5	37	Nc
SARSCoV_2000_15_LEFT	AGGAGTCAAATTACATTACACATAAAC	forward	25360	25389	30	1805	60.1	30	None	Nc
SARSCoV_2000_15_RIGHT	ACTGCTACTGGAATGGTCTGTGT	reverse	27142	27164	23	1805	61.6	47.8	None	Nc
SARSCoV_2000_17_LEFT	ACTTGTACGCCTAAACGAACA	forward	27873	27894	22	1918	60.7	45.5	36.7	Nc
SARSCoV_2000_17_RIGHT	TAGGCAGCTCTCCCTAGCATTG	reverse	29769	29790	22	1918	61.6	54.5	45.3	Nc




Pool1

POOL2. Average Tm  61.05 °C

A	B	C	D	E	F	G	H	I	J	K
Name	Sequence	Direction	Start	End	Length	Product Siz	Tm	%GC	Hairpin Tm	Pa
SARSCoV_2000_02_LEFT	AGGCCGCTATAACAATACTAGATGGA	forward	1956	1981	26	1923	61.3	42.3	None	Nc
SARSCoV_2000_02_RIGHT	CAGCGATCTTTTGTCAACTTGCT	reverse	3855	3878	24	1923	60.8	41.7	None	Nc
SARSCoV_2000_04_LEFT	TCAACATGCCAATTTAGATTCTTGCA	forward	5473	5498	26	1929	60.3	34.6	None	8.1
SARSCoV_2000_04_RIGHT	GCTGAAATCGGGGCCATTGTGA	reverse	7380	7401	22	1929	61.5	50	None	8.1
SARSCoV_2000_06_LEFT	GCTGCTGAATGTACAATTTTAAAGATG	forward	9011	9039	29	2003	61.1	34.5	None	Nc
SARSCoV_2000_06_RIGHT	AACCAGTGGTGTGTACCCTTGA	reverse	10992	11013	22	2003	61.5	50	47	Nc
SARSCoV_2000_08_LEFT	TCACCTAATTTAGCATGGCCTCTT	forward	12620	12643	24	1956	60.1	41.7	None	Nc
SARSCoV_2000_08_RIGHT	CAGGGTCAGCAGCATACACAAG	reverse	14554	14575	22	1956	61.5	54.5	None	Nc
SARSCoV_2000_10_LEFT	TGCATACGTAGACCATTCTTATGTTGT	forward	16291	16317	27	1985	60.8	37	32.8	Nc
SARSCoV_2000_10_RIGHT	GCTTCTTCGCGGGTGATAAACA	reverse	18254	18275	22	1985	61.5	50	None	Nc
SARSCoV_2000_12_LEFT	GGACTACAAAAGAGATGCTCCAGC	forward	19878	19901	24	1920	61.5	50	42.5	11
SARSCoV_2000_12_RIGHT	ACCTCTTAGTACCATTGGTCCCA	reverse	21775	21797	23	1920	60.5	47.8	37.1	11
SARSCoV_2000_14_LEFT	GCTGAACATGTCAACAACATCATATGA	forward	23519	23544	26	1973	60.1	38.5	35.2	Nc
SARSCoV_2000_14_RIGHT	TGCAGTAGCGCGAACAATACT	reverse	25470	25491	22	1973	61.4	45.5	46.2	Nc
SARSCoV_2000_16_LEFT	TCTTATTACAAATTGGGAGCTTCGCA	forward	27051	27076	26	1995	61.3	38.5	37.6	Nc
SARSCoV_2000_16_RIGHT	GCTTCTTAGAAGCCTCAGCAGC	reverse	29024	29045	22	1995	61.6	54.5	43.6	Nc

Pool2.

6 PRIMER STOCKS (100 micromolar (μM))

If you have ordered each primer independently and need to generate primer pool stocks: add  5 μL of each primer from Pool 1 to a  1.5 mL Eppendorf labeled "Pool 1 (100μM)" and each primer from Pool 2 to a  1.5 mL Eppendorf labelled "Pool 2 (100μM)". These are your 100 micromolar (μM) stocks of each primer pool.


Note



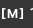

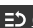

7 WORKING PRIMERS (10 micromolar (μM))

Dilute the primer stocks 1:10 in molecular grade water, to generate 10μM primer working stocks. It is recommend that multiple aliquots of each primer pool are made to in case of degradation or contamination.


Note





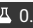

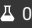
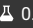











8 POOLING OF PRIMERS

Label two  1.5 mL Eppendorf tubes, one as **POOL1** and the another as **POOL2**.


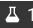
- 8.1 Add  10 μL of each of the primers of **set 1**  [go to step #7](#) to the Eppendorf tube labelled **POOL1**, final concentration will be  10 micromolar (μM) each
- 8.2 Add  10 μL of each of the primers of **set 2**  [go to step #7](#) to the Eppendorf tube labelled **POOL2**, final concentration will be  10 micromolar (μM) each

Multiplex PCR

- 9 In the PCR hood set up the multiplex mastermix reaction as follows in 2 0.2mL PCR tubes. We use the  KAPA Taq HotStart **Sigma Aldrich Catalog #KK1510**

Component	Pool 1	Pool 2	Final Concentration
5X Kapa HotStart Buffer	 2.5 μL	 2.5 μL	 1 X
25 mM MgCl ₂	 0.75 μL	 0.75 μL	 1.5 millimolar (mM)
10 mM dNTPs	 0.25 μL	 0.25 μL	 0.2 millimolar (mM) each
Primer Pool 1 or 2 (10 μM ea)	 0.6 μL	 0.6 μL	 0.5 micromolar (μM)
Kapa Taq HotStart Polymerase	 0.1 μL	 0.1 μL	 0.5 U
Nuclease-free water	 5.8 μL	 5.8 μL	
Final mastermix volume	 10.0 μL	 10.0 μL	

Note




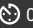



- 10 In the **extraction and sample addition cabinet** add  2.5 μL cDNA to each tube and mix well by pipetting.
The final volume will be  12.5 μL

Note

- 11 Pulse centrifuge the tubes to collect the contents at the bottom of the tube.

- 12 Set-up the following program on the thermal cycler:

5m 45s

Step	Temperature	Time	Cycles
Heat Activation	 98 °C	 00:00:30	1
Denaturation	 95 °C	 00:00:15	25-35
Annealing and Extension	 60 °C	 00:05:00	25-35
Hold	 4 °C	Indefinite	1

Note

Note

Expected result

Pooling and PCR quantification

- 13 Amplicon quantification to make an equimolar mixture.

Note

- 13.1 Put $1\ \mu\text{L}$ of each pool in a Nanodrop or similar spectrometer and quantify the DNA concentration.

- 13.2 Label a $1.5\ \text{mL}$ Eppendorf tube for each sample and make a **equimolar mix** with the two pools.
Calculate to achieve a final concentration of $50\ \text{Mass Percent}$

- 14 Quantify DNA using a **Qubit** or other method.
Quantification using Nanodrop is not recommended for a good estimation of the final pool.

Note

- 14.1 Protocol

Protocol



NAME

DNA quantification using the Qubit fluorometer


CREATED BY


Nikki Freed

PREVIEW

- 14.1.1 Prepare a mastermix of Qubit™ working solution for the required number of samples and standards. The Qubit dsDNA kit requires 2 standards for calibration (see note below).

Per sample:

Qubit® dsDNA HS Reagent  1 µL

Qubit® dsDNA HS Buffer  199 µL

Note

- 14.1.2 Label the tube lids. Do not label the side of the tube as this could interfere with the sample reading.

Note

- 14.1.3 Aliquot Qubit™ working solution to each tube:
- standard tubes requires 190µL of Qubit™ working solution
 - sample tubes require anywhere from 180–199µL (depending how much sample you wish to add).

The final volume in each tube must be 200µL once sample/standard has been added.

- 14.1.4 Add 10µL of standard to the appropriate tube.


- 14.1.5 Add 1–20µL of each user sample to the appropriate tube.



Note

- 14.1.6 Mix each tube vigorously by vortexing for 3–5 seconds.
- 14.1.7 Allow all tubes to incubate at room temperature for 2 minutes, then proceed to “Read standards and samples”.
- 14.1.8 On the Home screen of the Qubit™ 3 Fluorometer, press DNA, then select 1X dsDNA HS as the assay type. The Read standards screen is displayed. Press Read Standards to proceed.
- Note**
- 14.1.9 Insert the tube containing Standard #1 into the sample chamber, close the lid, then press Read standard. When the reading is complete (~3 seconds), remove Standard #1.
- 14.1.10 Insert the tube containing Standard #2 into the sample chamber, close the lid, then press Read standard. When the reading is complete, remove Standard #2.
- 14.1.11 The instrument displays the results on the Read standard screen. For information on interpreting the calibration results, refer to the Qubit™ Fluorometer User Guide, available for download at thermofisher.com/qubit.
- 14.1.12 Press Run samples.
- 14.1.13 On the assay screen, select the sample volume and units:
- Press the + or – buttons on the wheel, or anywhere on the wheel itself, to select the sample volume added to the assay tube (from 1–20µL).
 - From the unit dropdown menu, select the units for the output sample concentration (in this case choose ng/µL).
- 14.1.14 Insert a sample tube into the sample chamber, close the lid, then press Read tube. When the reading is complete (~3 seconds), remove the sample tube.
- 14.1.15 **The top value (in large font) is the concentration of the original sample and the bottom value is the dilution concentration.** For information on interpreting the sample results, refer to the Qubit™ Fluorometer User Guide.
- 14.1.16 Repeat step 14 until all samples have been read.
- 14.1.17 Carefully **record all results** and store run file from the Qubit on a memory stick.

- 14.1.18 All negative controls should ideally be 'too low' to read on the Qubit machine, but MUST be < 1ng per ul. If your negative controls >1ng per ul, considerable contamination has occurred and you must redo previous steps.

Normalisation


- 15 Label a  0.2 mL PCR tube for each sample.

- 15.1 Adjust the amount of DNA in the tube to be  100 ng total per sample in  7.5 µL molecular grade water.




Note

Tagmentation

- 16  Nextera XT DNA Library Preparation Kit **Illumina, Inc. Catalog #FC-131-1024** We have reduced the amount of reagents used per reaction.

- 16.1 Label a  1.5 mL Eppendorf tube for each sample and add in this order:




Component	Volume
-----------	--------


- | | |
|---|---|
| 1. Tagment DNA buffer (TD) |  3.75 µL |
| 2. cDNA  2 ng/µl |  3 µL |

Pipette 10 times to mix

- 16.2 Add **Amplicon Tagment Mix**  0.5 µL and pipette 10 times to mix

- 16.3 Centrifuge at  280 x g, 20°C, 00:01:00

- 16.4 Incubate in thermal cycler
 55 °C for  00:10:00
 10 °C indefinitely

- 16.5 Add  2 µL of Neutralize Tagment Buffer (NT).
Pipette 10 times to mix

- 16.6 Centrifuge at  280 x g, 20°C, 00:01:00

16.7 Incubate at Room temperature for 00:05:00 0 °C

5m

16.8 Preserve the samples at 4-8 °C until use.

Indexing

30m

17 La siguiente reacción requiere de la enzima polimerasa 2x Ampigene HS Taq Mix Catalog # ENZ-NUC101-0200 y de los sets de Nextera XT. Cada set contiene 96 combinaciones de TAGs con esto podemos alcanzar 384 muestras usando los Nextera XT Index Kit v2 (Sets A,B,C y D) Catalog # 20027213;20027214;20027215;20027216.

17.1 A cada tubo de reacción agregar:

Component	Volume
Nuclease-free water	10 µL
2x Ampigene HS Taq Mix.	10 µL
Nextera XT index i5	1 µL
Nextera XT index i7	1 µL
ADN Tagmentado	3 µL
Final volume	25 µL

17.2 Centrifuge at 280 x g, 20°C, 00:01:00

1m

17.3 Set-up the following program on the thermal cycler:

Cover 100 °C

Step	Temperature	Time	Cycles
Heat Activation	72 °C	00:03:00	1
Denaturation	95 °C	00:00:30	1
Annealing and Extension	95 °C	00:00:10	14
	55 °C	00:00:30	
Hold	72 °C	00:00:30	1
	4 °C	Indefinite	

4m 40s

18

Purificación final y pooling

19 Al volumen que se tiene en el tubo agregar 0.8X de perlas magnéticas Ampure XP.

19.1 Mezclar y spin-down.

- 19.2 Llevar al magneto hasta formar el pellet y desechar el sobrenadante.
- 19.3 Agregar 150 uL de etanol 80% en posición contraria al pellet. Esperar 30 segundos.
Desechar el etanol.
- 19.4 Retirar el exceso de etanol y secar las perlas por 5 min.
- 19.5 Resuspender las perlas en 26 uL de Agua libre de nucleasas.
Mix y spin down. Incubar por 5 min a TA.
- 19.6 Llevar al magneto hasta formar el pellet y transferir 25 uL del sobrenadante a un nuevo tubo previamente etiquetado.
- 19.7 Cuantificar 2 uL por Qubit HS y analizar los tamaños mediante electroforesis en gel de agarosa al 1.0%.

Secuenciación de bibliotecas en plataforma Illumina Miniseq

- 20 A partir de la concentración en ng/uL determinada por Qubit y obtenido el tamaño aproximado del fragmento, llevar cada una de las librerías a una concentración de 4 nM.
Transferir 5 uL de cada librerías a 4 nM a un tubo previamente etiquetado para obtener el pool final.
Seguir el protocolo Library Denaturing and miniseq Sample Loading del kit