

VERSION 2

NOV 08, 2022

 SHAREWORKS FOR ME  Assembly: Chronic recoverable Neuropixels in mice V.2**This protocol is published without a DOI.**Emily A Aery Jones¹¹Stanford University EEmily A Aery Jones
Stanford University

COMMENTS 0

ABSTRACT

This protocol collection explains how to build a low-cost, lightweight system to implant Neuropixels 1.0 probes into mice, record during freely moving behavior, then recover the probe for future use. This protocol explains how to 3D print components, sharpen, solder, and test the probe, and prepare components for surgery.

Chronic recoverable implant based off of Juavinett et al, 2019 (<https://elifesciences.org/articles/47188>), then modified by Janna Aarse, then modified by Emily Aery Jones.

See also other chronic recoverable designs:

* for rats: Luo & Bondy et al, 2020 (<https://elifesciences.org/articles/59716>)

* Neuropixels 2.0: Steinmetz, Aydin, Lebedeva, Okun, & Pachitariu et al, 2021 (<https://www.science.org/doi/10.1126/science.abf4588>)

* overview: Neuropixels 2021 course lecture 3.7 by Yoh Isogai (https://www.youtube.com/watch?v=7cZqYIiGvBQ&list=PLfhWmWntvjl64ti_a-MzHlwqwEU0ZIALb&index=18)

PROTOCOL CITATION

Emily A Aery Jones 2022. Assembly: Chronic recoverable Neuropixels in mice. **protocols.io**
<https://protocols.io/view/assembly-chronic-recoverable-neuropixels-in-mice-cizcuf2w>
Version created by Emily A Aery Jones



KEYWORDS

electrophysiology, silicon probe, entorhinal cortex, probe sharpening, ground screw

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Nov 08, 2022

LAST MODIFIED

Nov 08, 2022

PROTOCOL INTEGER ID

72452

GUIDELINES

- These steps can be done in any order and simultaneously, except that printing and sharpening must happen before gluing, and soldering must happen before testing.
- The probe shank is fragile and hard to see. Secure the probe in a closed container whenever not in use, or secure in the Sensapex holder to a clamp or other immobilizing device with the shank pointed well clear of any items the experimenter might reach for to prevent accidental breakage.

MATERIALS TEXT

Consumables:

3D printing:

- Formlabs black resin (RS-F2-GPBK-04)
- Isopropyl alcohol

Probe box:

- non-hardening modeling clay
- plastic case, >4cm tall and >10cm wide (e.g. Amazon B01M9JU210)
- Neuropixels probe with metal dovetail

Soldering:

- Silver wire (WPI AGW1030)
- Solder
- Soldering flux
- Precision applicator brushes (Parkell S379)
- Gold pins (Digikey ED1058-ND)
- Ground screws (00-90 1/8" stainless steel screws, Amazon B002SG89S4)

Testing:

- Saline & DH2O

Assembling:

- Assembly screws (00-90 1/8" brass round head screws, McMaster-Carr 92453A854)

- Hex Nuts (00-90 brass, McMaster-Carr 92736A112)
- Loctite glue
- 1" Headbar
- Sterile 200uL pipette tips

Equipment:

3D printing:

- Formlabs Form 2 printer
- Formlabs finish kit

Sharpening:

- Narishige EG-40 grindstone
- Sensapex holder with 0.89mm hex key
- Thin metal rod for Sensapex holder (see "Before start")

Soldering:

- Soldering iron
- Forceps
- Soldering clamps

Testing:

- Recording system
- Beaker
- Clamp or micromanipulator

Assembling:

- Tiny flathead screwdriver
- Sterilized scissors
- Plastic jewelry divider box for holding printed parts

BEFORE STARTING

- If your probes do not have a metal cap, you can 3D print the file provided at https://github.com/emilyasterjones/chronic_NPX_mouse/tree/main/probe_holders, or machine according to this file: https://e2f49b05-3010-45e2-a7bf-b8f67811425a.filesusr.com/ugd/832f20_1617705787d84f7795ec85aa08630ea1.pdf. To attach the cap, secure the probe in the modeling clay, apply superglue to dovetail cap, position, press together for  00:00:30, and allow to cure  Overnight.
- The grindstone only accepts thin metal rods, not those compatible with stereotaxes. Machine a compatible rod using the file provided at https://github.com/emilyasterjones/chronic_NPX_mouse/tree/main/probe_holders.

4h 40m

3D print components

10m

- 1 Build a print file for the following pieces per mouse: 1 each of body piece, back and front flex cable holders, and dome, plus 2 wings. Print 1 headstage holder per recording rig. To re-use explanted probes, print everything

except for the body piece, which is permanently affixed to the probe. Files located at https://github.com/emilyasterjones/chronic_NPX_mouse.

Note

These components have been tested at 50um resolution on a Form2 printer with black resin. Different printers and resolutions (e.g. ESD resin on a Form3 printer) may also work.

Note

You only need to print the headstage holder if you will record freely moving. See *Freely moving recording* protocol for how to wire and assemble the headstage holder.

Note

For freely moving recordings, the headstage holder must pull upwards on the implant, but is attached only to the flex cable holder, which is during surgery affixed to the skull only through a bit of dental cement connecting it to the dome piece. This attachment is secure enough, but for a more secure attachment, you can opt to print the "with_screw" versions of the flex holder and body piece, then affix the flex holder to the body piece with a screw prior to surgery. Note that this design prevents the body piece from being mounted on the posterior of the skull as you will not be able to visualize the probe through the body piece. Read the *Implant surgery* protocol for more detail.

- 1.1 Place all build files into a single print file. Orient each component so it is well-supported, with supports attaching to non-interface points. These are: rounded hooks of wings, top of body piece (where hex nut slot is), flat backs of flex cable holders and headstage holder, and any side of dome. See *single_mouse_print.form* file for example using Formlabs system.

4h

- 2 Print the file.

- 3 Remove prints from the build platform. Remove liquid resin from prints and cure according to manufacturer instructions for your printer and resin.

10m

Note

For Formlabs black resin: use Finish Kit to wash in 2 isopropyl alcohol baths, manually agitating for a few


minutes in each bath, to completely remove uncured resin. Allow to dry. No UV cure required. See full instructions at <https://support.formlabs.com/s/article/Form-2-Basic-Finishing-Steps>.

- 4 Remove supports with wire cutters for fine surfaces & twisting print & raft apart for larger surfaces.

20m

Build probe box

- 5 Along the base inside the plastic case, place a thick (>1cm) piece of modeling clay in a ~3cm strip about one-quarter of the way from the top of the box. The clay should be thick enough so that you can press a probe into it without worrying about the shank hitting the box and wide enough to securely hold the PCB board of the probe.
- 6 From an empty probe box, remove the foam interior. Cut the probe holding foam strip out of this. Tape this strip to the top of the box. This will grip any metal rod attached to the Sensapex holder to secure the probe during soldering, gluing, or waiting to be mounted to the stereotax or clamp.

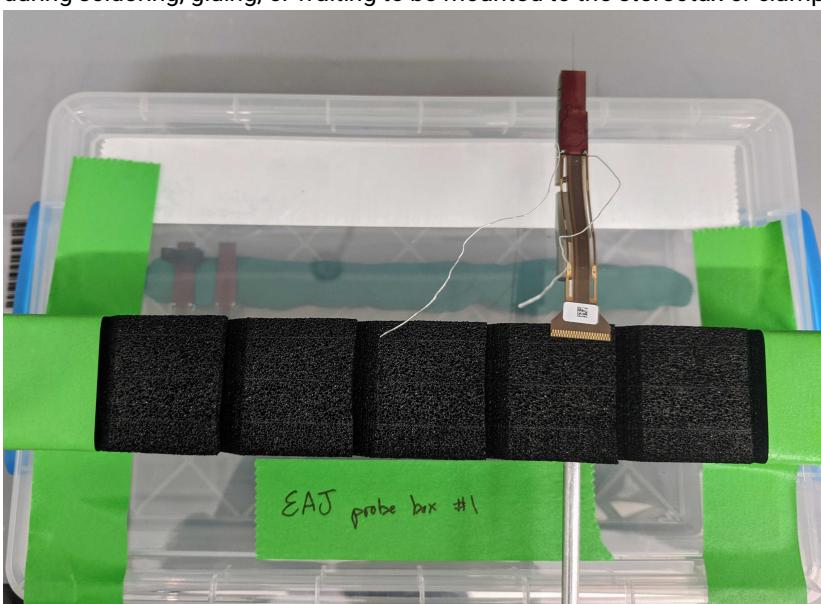

A probe box with modeling clay along the interior and foam gripper along the top.

Sharpen the probe


- 7 Load the probe into the Sensapex holder. Screw the holder into the thin metal rod.

20m

- 8 Sharpening probes can improve cell yield and allow you to puncture dura without removing it. Sharpen probe according to this protocol: <https://github.com/cortex-lab/neuropixels/wiki/Sharpening>



Probe mounted to grindstone at 20 degrees.

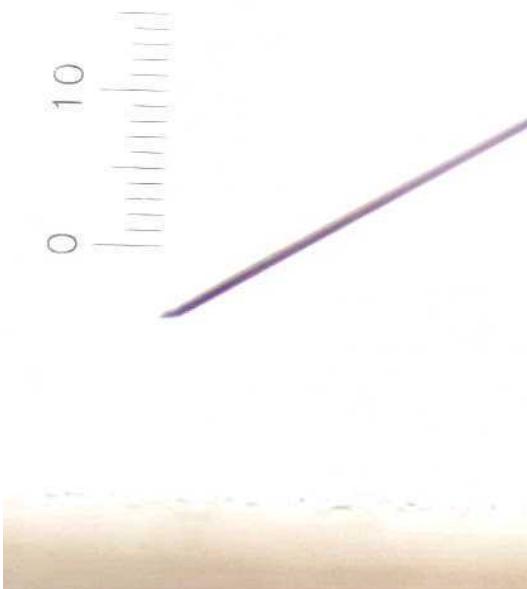

Sharpened probe

Solder probe and ground screw

- 9 Thread silver wire between the ground and reference pads on the back of the PCB board, then down to a few

mm below the PCB. Solder. Keep the iron cool (~315C/600F) and don't heat for longer than 4s.

Note

Total length covered by wires & pins, between top of screw & pads on PCB, should be at least 3cm, but not much longer. This is about the distance from the PCB pads to the tip of shank.

Note

With the Sensapex holder facing away from you and the probe pointed down, ground pad is on the right and reference pad is on the left (ref=left). Shorting ground to reference is not necessary for internally referenced recordings (just convenient, in case you swap which pad is which or want to externally reference later on). Shorting ground to reference works for most experimenters using external referencing.

Safety information

Lead solder creates toxic fumes. Have a ventilation fan or downdraft table pulling melted solder fumes away from the user.

Safety information

Soldering flux is toxic. Dispose of all components which touched the flux (e.g. weighboats, precision applicators) appropriately.

10 Solder gold male pin to the end of wire.

11 Apply flux just under the head of the ground screw. Wrap a loop of wire around this & twist to tighten. Solder closed.

Note

You will need one screw per mouse, or two if keeping reference and ground separate.

- 12 Solder gold female pin to the end of this wire.

3m


Test signal

- 13 Mount the probe in a clamp or micromanipulator and submerge shank into saline. Clip the ground wire into the saline on the side of the beaker.
- 14 Plug in the headstage. Run BIST tests. Observe the signal and noise level on SpikeGLX to confirm your soldering is good and the probe is functional.

10m

Assemble

- 15 Superglue nuts into slots on body piece. Wait for glue to dry  00:00:30. Attach wings and screw to affix.
- 16 Superglue the back of the probe to the 3D printed body piece. The wings extend 6mm beyond the body piece and you'll need some room for dental cement, so set the piece at [your target depth]-3mm away from the base of the shank. Press together for  00:00:30. Allow to cure  Overnight.

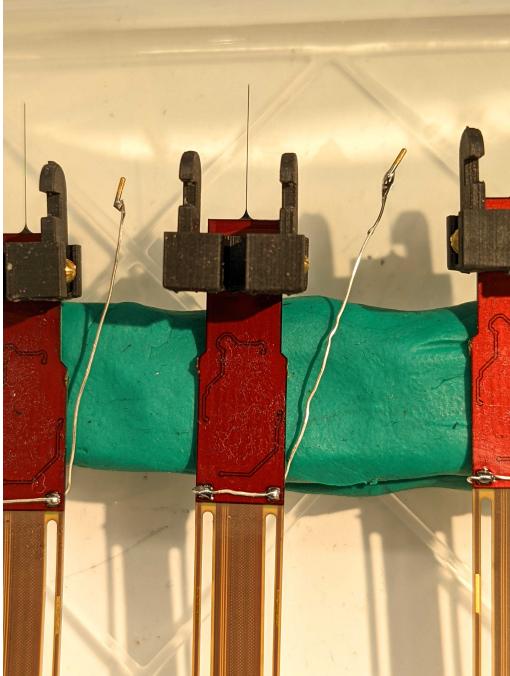

Soldered and glued probes mounted on modeling clay inside plastic storage case.

Note

Body piece can be attached centered to reach central structures or allow targeting to either hemisphere, or can be attached to one side to allow implant to be more centered on the skull when targeting lateral structures. Center-attached body pieces work well to reach MEC (3.3mm from midline).

17 Mark the center of each headplate (1" long) with a lab marker.

18 Slice the pipette tips into 1mm diameter, 0.5mm tall circles to create wells for the craniotomy.



Implant components not mounted to probe. Top row: well, headbar, ground screw. Bottom row: flex cable holder with tab slot, flex cable holder without tab slot, dome.

1h

Build the headstage holder

- 19 Insert the LEDs through the ends of the headstage holder arms. Wire: LED short lead => resistor => ground, LED long lead => switched power. Insert a coin cell battery and flip the switch to check the connection.

Note

Green LED is much dimmer than red LED and so does not need a resistor, to try to balance the illumination levels.

- 20 Epoxy the battery breakout to the back of the headstage holder. Mount headstage into slot and optionally secure with tape. Plug into Omnetics connector and secure this connection with tape (this connection easily comes loose).



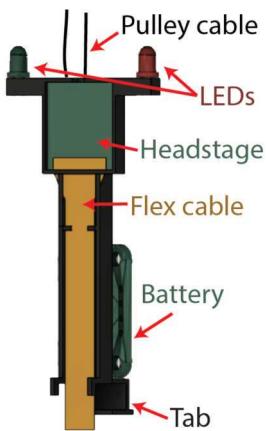


Diagram of headstage holder, version with LEDs mounted on top. Pulley cable loops into slots on top of holder to connect to counterweight, 12mm coin battery breakout board glued to back near base, and tab inserts into tab slot to attach to flex cable holder on mouse.

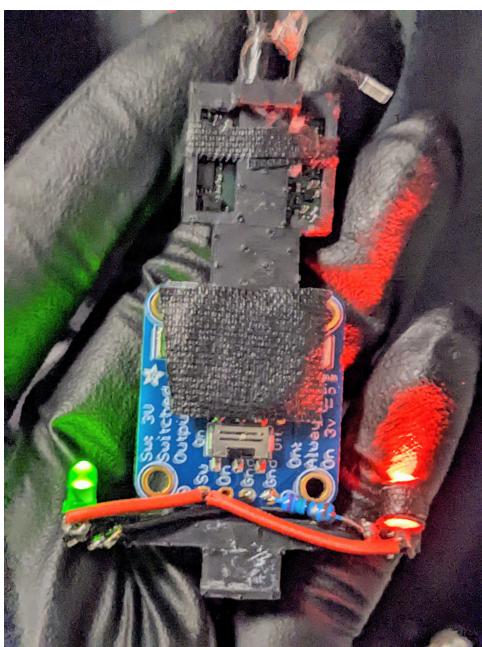

Headstage holder from the back. LEDs are wired to switched power (red wires) and ground (black wire/resistor). In this design, red LED is surrounded by tape to reduce diffusion; battery is covered with tape to reduce reflection from overhead lights.

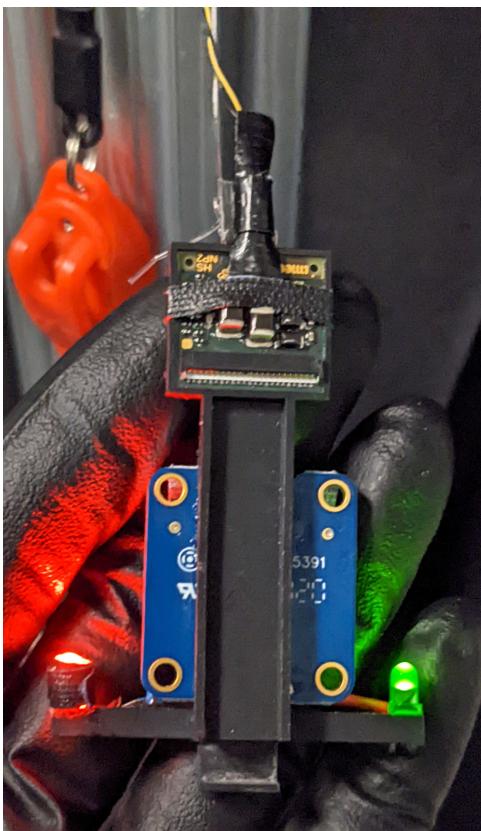

Heastage holder from the front. Pulley cable (clear) emerges from behind the headstage.