

Aug 29, 2024

Immunofluorescence and confocal microscopy

DOI

dx.doi.org/10.17504/protocols.io.6qpvr8p1olmk/v1

Elias Adriaenssens¹

¹Sascha Martens lab, University of Vienna, Max Perutz Labs - Vienna



Elias Adriaenssens

Sascha Martens lab, University of Vienna, Max Perutz Labs - ...

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.6qpvr8p1olmk/v1

Protocol Citation: Elias Adriaenssens 2024. Immunofluorescence and confocal microscopy. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.6qpvr8p1olmk/v1>

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: July 02, 2024

Last Modified: August 29, 2024

Protocol Integer ID: 103082

Keywords: ASAPCRN

Funders Acknowledgement:

**Aligning Science Across
Parkinson's (ASAP)**

Grant ID: ASAP-000350

**Marie Skłodowska-Curie
MSCA Postdoctoral
fellowship**

Grant ID: 101062916



Abstract


This protocol details the process of immunofluorescence and confocal microscopy.

Materials

28906, Thermo Fisher Scientific

 Pierce™ 16% Formaldehyde (w/v) Methanol-free **Thermo Fisher Scientific Catalog #28906**

0.1% (v/v) Triton X-100 (9002-93-1, Sigma-Aldrich)













5% (v/v) BSA (9048-46-8, Sigma-Aldrich)  Bovine Serum Albumin [BSA] **Fisher Scientific Catalog #9048-46-8**

DAPI Fluoromount-G mounting medium (0100-20, Southern Biotech)

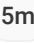

Steps

2h 15m

- 1 Seed the cells on glass coverslips (12 mm #1.5) at a concentration of 100.000 cells/well, and after treatment with Rapalog for the indicated time, fix in 4% paraformaldehyde (28906, Thermo Fisher Scientific) for  00:10:00 at  Room temperature .  10m
- 2 After washing with PBS, permeabilize the cells with 0.1% (v/v) Triton X-100 (9002-93-1, Sigma-Aldrich) in PBS for  00:05:00 .  5m
- 3 Perform the blocking with blocking buffer (5% (v/v) BSA (9048-46-8, Sigma-Aldrich) and 0.05% (v/v) Triton X-100 diluted in PBS) for  01:00:00 at  Room temperature .  1h
- 4 Dilute the primary and secondary antibodies in blocking buffer and incubate for  01:00:00 at  Room temperature with three PBS washing steps in between.  1h
- 5 Mount the cells on microscopy slides in DAPI Fluoromount-G mounting medium (0100-20, Southern Biotech), which stains the nuclei, and store at  4 °C until use.
- 6 Perform the confocal microscopy with a Zeiss LSM700 laser scanning confocal microscopy with Plan-Apochromat 40×/1.30 Oil DIC, WD 0.21 mm objective.