


Sep 16, 2024

# Micro Volume Purification of His-Tagged Proteins with IMAC Ni-Charged Resin

DOI

[dx.doi.org/10.17504/protocols.io.3byl49yejgo5/v1](https://dx.doi.org/10.17504/protocols.io.3byl49yejgo5/v1)

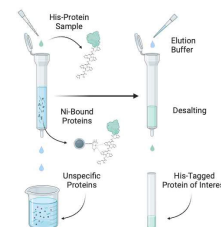

Syon Schlecht<sup>1,2</sup>, Emily Gunderson<sup>2</sup>, Ruthie Fowler<sup>1,2</sup>, Takara Aguilar<sup>3,2</sup>

<sup>1</sup>Georgia Institute of Technology; <sup>2</sup>Empower College and Career Center; <sup>3</sup>University of Georgia

GA Tech School of Chemi...



Syon Schlecht

Georgia Insitute of Technology

OPEN  ACCESS



DOI: [dx.doi.org/10.17504/protocols.io.3byl49yejgo5/v1](https://dx.doi.org/10.17504/protocols.io.3byl49yejgo5/v1)

**Protocol Citation:** Syon Schlecht, Emily Gunderson, Ruthie Fowler, Takara Aguilar 2024. Micro Volume Purification of His-Tagged Proteins with IMAC Ni-Charged Resin. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.3byl49yejgo5/v1>

**License:** This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**We use this protocol and it's working**

**Created:** September 09, 2024

**Last Modified:** September 16, 2024

**Protocol Integer ID:** 107141

## Disclaimer

### DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to **protocols.io** is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with **protocols.io**, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.







## Abstract

This protocol describes a micro volume purification method for His-tagged proteins using IMAC with Ni-charged resin in a mini spin chromatography column. His-tagged proteins are selective bound to a Ni-charged resin due to the interaction between the polyhistidine tag and  $\text{Ni}^{2+}$  ions and eluted using imidazole-containing buffers. Desalting is the performed using SEC where salts and imidazole are fractionated, while larger molecules like the target protein are excluded. This results in the isolation of purified His-tagged proteins with a high yield and minimal contaminants, suitable for downstream applications.

## Image Attribution

Created in BioRender. Schlecht, S. (2024) BioRender.com/d98d031










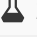


## Protocol materials

-  Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** In [11 steps](#)
-  Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories** In [6 steps](#)
-  Profinity™ Ni-charged IMAC Resin **Bio-Rad Laboratories Catalog #1560131** Step 1
-  Equilibration Buffer Step 5
-  Wash Buffer Step 9
-  Elution Buffer Step 10





## Preparing and Equilibrating Column

6m

- 1 Resuspend  
 Profinity™ Ni-charged IMAC Resin **Bio-Rad Laboratories Catalog #1560131** and  
pipette  200 µL into a  
 Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories**
- 2 Centrifuge the  Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** at  
 1000 x g, 00:02:00 and discard any collected buffer
- 3 Wash the  Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** by adding  
 200 µL of distilled water
- 4 Centrifuge the  Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** at  
 1000 x g, 00:02:00 and discard any water collected
- 5 Equilibrate the column by adding  200 µL of  Equilibration Buffer to the  
 Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories**




### Note

The equilibration buffer is 20 mM sodium phosphate, 300 mM NaCl, and 5 mM imidazole.

- 6 Centrifuge the  Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** at  
 1000 x g, 00:02:00 and discard any buffer collected and cap the column

## Sample Binding and Elution of His-Proteins

20m

- 7 From the soluble fraction of lysed cells containing the  His-Sample add  600 µL to  
the  Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** and gently mix for

20m



00:20:00

8 Centrifuge the Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** at 1000 x g, 00:02:00 and discard any flowthrough 2m

9 Add 600 µL of Wash Buffer to the Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** and centrifuge at 1000 x g, 00:02:00 and discard any collected buffer 2m

#### Note

The wash buffer is 20 mM sodium phosphate, 300 mM NaCl and 10 mM imidazole.

10 Add 100 µL Elution Buffer to the Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** and centrifuge at 1000 x g, 00:02:00 , and **COLLECT** the flowthrough, **containing** the His-Sample 2m

#### Note

The elution buffer is 20 mM sodium phosphate, 300 mM NaCl, 250 mM imidazole, pH 8.0

#### Note

The Mini Bio-Spin® Chromatography Column **Bio-Rad Laboratories** can now be disposed of

## Desalting of His-Sample 2m

11 Snap the bottom and remove the top of a Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories** and allow the 5m



packing buffer to drain for 00:05:00

#### Note

The Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories** has a fractionation range of 1000 to 6000 Da, so any remaining salts or small proteins will be fractionated.

- 12 Centrifuge the Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories** at 1000 x g, 00:02:00 and discard any collected buffer

2m

- 13 Add the approx. 100 µL of the collected His-Sample to the Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories**

#### Note

The maximum load volume of the column is 100 µL

- 14 Centrifuge the Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories** at 1000 x g, 00:05:00 and any excluded liquid is the purified His-Sample

5m

- 15 Repeat up to 3 times with the same Micro Bio-Spin® Column with Bio-Gel® P-6 **Bio-Rad Laboratories** to obtain the desired His-Sample volume



## Protocol references

Protein Expression and Purification Series.

[https://www.bio-rad.com/webroot/web/pdf/lse/literature/pepsi\\_hr\\_1665067.pdf](https://www.bio-rad.com/webroot/web/pdf/lse/literature/pepsi_hr_1665067.pdf)

Schlecht, S., Gunderson, E., Fowler, R., & Aguilar, T. (2024). Synthesis and Characterization of Naproxen-Salicylate Derivatives as Potential Dual-Targeted Inhibitors of Dihydrofolate Reductase. *Advances in Biological Chemistry*, 14(04), 87–102. <https://doi.org/10.4236/abc.2024.144008>

His-tagged Proteins – Production and Purification. ThermoFisher Scientific.

<https://www.thermofisher.com/us/en/home/life-science/protein-biology/protein-biology-learning-center/protein-biology-resource-library/pierce-protein-methods/his-tagged-proteins-production-purification.html>