




MAR 22, 2023

Water Extractable Organic Matter (WEOM)

In 2 collections

maggie.bowman¹

¹Environmental Molecular Sciences Laboratory

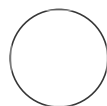

nicholas.sconzo

ABSTRACT

The method covers the extraction of water extractable organic matter for soil organic matter characterization using FT-ICR-MS, LCMS, and TOC/TN.

MATERIALS

- Genese Scientific Centrifuge Tubes 50 mL (9 per sample)
- Genese Scientific Centrifuge Tubes 15 mL (3 per sample)
- Air dried soil
- Ultrapure DI water
- 30 mL Pipet and auto-pipettor
- 1000 µL and 2 µL pipets
- Shaker (800 rpm)
- Centrifuge
- Sample Labels for each analysis
- 0.45 um PES Whatan Syringe Filters
- Syringe

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.ewov1o4oylr2/v1

Protocol Citation: maggie.bowman 2023. Water Extractable Organic Matter (WEOM). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.ewov1o4oylr2/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working




Created: Oct 06, 2022


Last Modified: Mar 22, 2023


PROTOCOL integer ID:
70969

Extraction

- 1 Weigh out 6 g of air-dried soil in triplicate into a labeled 50 mL falcon tube

2 Add  30 mL of Ultrapure DI (maintaining a  2 g to  10 mL ratio) using large pipette

3 Place samples onto orbital shaker for  2 rpm, 02:00:00

4 Once completed centrifuge samples at  6000 rpm, 00:08:00

8m

4.1 If samples are still cloudy, centrifuge a second time


5 Transfer the supernatant to two centrifuge tubes and do not discard the soil pellet

5.1  5 mL in a  15 mL centrifuge tube for FT-ICR-MS

5.2 ~  20 mL in a  50 mL centrifuge tube for TOC/TN and LC-MS

6 Acidify the FT-ICR-MS subsample using  2 μ L of concentrated phosphoric acid (37%)

7 Freeze the TOC/TN and LC-MS sample on the side at  -20 °C until filtration

8 Freeze FT-ICR-MS samples on their side in the  -80 °C until SPE


8.1 SPE Protocol: [Gilson Automated SPE Protocol](#)


Sample Filtering (LC-MS and TOC/TN)n




9 Thaw ~  20 mL sample prior to filtering

10 Using a new clean syringe for each sample, pull up  50 mL up Ultrapure DI water

11 Attach a new 0.45 um filter to the syringe and press the water through the filter to rinse the filter.

12 Remove the filter without touching the tip and pull up the ~  20 mL of thawed sample

13 Re attach the rinsed filter and discard the first  1 - 2 mL of sample filtered

- 14 Filter the sample into two falcon tubes, two  15 mL tube labeled for TOC and LC-MS respectively.
- 15 Filter  10 mL of solution into the LC-MS tube and all remaining sample into the TOC/TN tube.
- 16 If the filter becomes difficult to press sample through, rinse a new filter using a clean (DI ONLY syringe using steps 3 and 4) and attach the rinsed filter to the syringe with sample. Discard the first  1 - 2 mL and continue filtering into the respective tubes.
- 17 Store filtered samples in the freezer prior to analyses.