



JUN 13, 2023

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.eq2ly7b1mlx9/v1

Protocol Citation: michela.deleidi, Bianca Marchetti, Federico Bertoli, Carmela Giachino, [carmela.giachino](mailto:carmela.giachino@unict.it) 2023. Rotarod test.

protocols.io
<https://dx.doi.org/10.17504/protocols.io.eq2ly7b1mlx9/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Mar 31, 2023

Last Modified: Jun 13, 2023

PROTOCOL integer ID:
79823

Keywords: ROTAROD

Rotarod test

Federico

michela.deleidi^{1,2}, Bianca Marchetti^{3,4}, Bertoli^{1,2},
Carmela
Giachino³, [carmela.giachino](mailto:carmela.giachino@unict.it)⁵

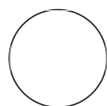
¹Mitochondria and Inflammation in Neurodegenerative Diseases, DZNE, Tübingen-Germany;

²Hertie Institute for Clinical Brain Research, University of Tübingen;

³Neuropharmacology Laboratory, Oasi Research Institute-IRCCS, Troina, Italy;

⁴Biomedical and Biotechnological Sciences, Pharmacology Section, University of Catania-Italy;

⁵Institut Imagine, INSERM UMR1163 Paris Cite' University, 24 boulevard du Montparnasse, 75015, Paris, France.

[carmela.giachino](mailto:carmela.giachino@unict.it)

ABSTRACT

Rotarod test was used to evaluate the effect of LPS treatment on motor performance during aging process between genotypes.

ATTACHMENTS

[k4dvb48if.docx](#)

- 1 To test motor performance, we used the accelerated Rotarod (five-lane accelerating rotarod; Ugo Basile, Comerio, Italy), a motor behavioral test widely used to assess motor deficit in neurodegenerative disease models in rodents. Mice have to keep their balance on a horizontal rotating rod (diameter, 3 cm) and rotation speed is progressively increased every 30 sec by 4 rpm.
- 2 Before the test, WT and G2019S within the different treatment groups were housed five per cage and acclimated to a 12h shift in light/dark cycle so that the exercise occurred during the animal normal wake period.
- 3 To start the trials, the mice (five mice are tested at the same time) are placed on a rotating rod; when the mice fall down or when 5 min are completed, a switch is activated that automatically stops a timer. On the day of testing, the mice perform 5 trials, separated by an interval of 30 min between each trial.
- 4 The Rotarod performance was assessed at day -7, 0 days, and during NaCL/LPS injection paradigm, at monthly intervals, for the determination of motor deficit (*Supplementary Figure 1*).