




Aug 24, 2022

Preparation and imaging of enriched Golgi from GolgiTAG-IP using Transmission Electron Microscopy

Rotimi Fasimoye¹, Alan Prescott², Dario R Alessi¹¹Medical Research Council Protein Phosphorylation and Ubiquitylation Unit, School of Life Sciences, University of Dundee, Dow Street, Dundee DD1 5EH, UK;²Dundee Imaging Facility, School of Life Sciences, University of Dundee, Dow Street, Dundee DD1 5EH, UK

1 Works for me

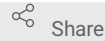

Share

dx.doi.org/10.17504/protocols.io.x54v9y9nqg3e/v1

asap

Dario R Alessi

ABSTRACT

Transmission electron microscopy (TEM) is a tool used to image, in good resolution, the structure of organelles within the cell. Available protocols are designed to image structures within fixed intact cells. Here, we described a protocol where Golgi, isolated from cells by GolgiTAG immunoprecipitation (IP) (as described in [ddx.doi.org/10.17504/protocols.io.6qpvrjrogm/v1](https://dx.doi.org/10.17504/protocols.io.6qpvrjrogm/v1)), can be prepared and imaged using TEM. This protocol can also be used to image any organelles isolated using various organelle-IP protocols that are available.

ATTACHMENTS

[ivgnbgrdf.docx](#)

DOI

dx.doi.org/10.17504/protocols.io.x54v9y9nqg3e/v1

PROTOCOL CITATION

Rotimi Fasimoye, Alan Prescott, Dario R Alessi 2022. Preparation and imaging of enriched Golgi from GolgiTAG-IP using Transmission Electron Microscopy. **protocols.io**
<https://protocols.io/view/preparation-and-imaging-of-enriched-golgi-from-gol-cfm6tk9e>

FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's (ASAP) initiative

Grant ID: ASAP-000463

UK Medical Research Council

Grant ID: MC_UU_00018/1

KEYWORDS

Enriched Golgi imaging, GolgiTAG-IP, Transmission Electron Microscopy

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 22, 2022

LAST MODIFIED

Aug 24, 2022

OWNERSHIP HISTORY

Aug 22, 2022 | Urmilas

Aug 23, 2022 | Dario R Alessi

PROTOCOL INTEGER ID






69022

MATERIALS TEXT




Materials:

- [☒ Paraformaldehyde Sigma](#)
- [Aldrich Catalog #158127](#)
- [☒ Glutaraldehyde EM Grade 25% Sigma](#)
- [Aldrich Catalog #G5882-50ML](#)
- [☒ Sodium cacodylate trihydrate Merck Millipore](#)
- [Sigma Catalog #C0250](#)
- Osmium tetroxide (Agar. Catalog # R1023)
- [☒ Sodium ferrocyanide decahydrate Sigma](#)
- [Aldrich Catalog #13425](#)
- Tannic acid (BDH. Catalog # 30337)
- Uranyl acetate (Agar. Catalog # R1260A)
- [☒ Durcupan™ ACM Sigma](#)
- [Aldrich Catalog #44611](#)
- [☒ Lead\(II\) citrate tribasic trihydrate Sigma](#)
- [Aldrich Catalog #15326](#)
- [☒ Sodium citrate monobasic Sigma](#)
- [Aldrich Catalog #71497](#)
- [☒ Sodium hydroxide \(NaOH\) Sigma](#)
- [Aldrich Catalog #S5881](#)
- Potassium chloride (KCl) (VWR. Catalog# 267264)
- [☒ Potassium phosphate monobasic Sigma –](#)
- [Aldrich Catalog #P5379](#)
- Calcium chloride (Calbiochem. Catalog# 208291)
- [☒ \(R\)-\(-\)-Propylene oxide Sigma](#)
- [Aldrich Catalog #540048](#)
- [☒ Ethyl alcohol, 200 proof, anhydrous, ≥99.5% Sigma](#)
- [Aldrich Catalog #459836](#)

Buffer:

- KPBS ([M]**136 millimolar (mM)** KCl, [M]**10 millimolar (mM)** KH₂PO₄. Adjust to pH 7.25 with KOH).
- [M]**0.4 Molarity (M)** sodium cacodylate buffer ( **17.12 g** sodium cacodylate dissolved in water. Adjust to pH7.4 with HCl).
- Reynold lead citrate ( **1.33 g** lead citrate,  **1.76 g** sodium citrate;  **8 mL** 1N NaOH in  **50 mL** water).


Equipment:

-  [Pierce™ Anti-HA Magnetic Beads Thermo](#)
- **Fisher Catalog #88837**
-  [DynaMag™- Spin Magnet Thermo](#)
- **Fisher Catalog #12320D**
- Microcentrifuge with thermostat (VWR Micro Star 17R. S/N 42209232)
- Scalpel
- P1000 pipette tips
- Pasteur pipette
- Oven (TAAB. I064)
- Leica UCT ultramicrotome (Leica)
-  [Clear glass with solid top black phenolic 14B rubber lined cap VWR international](#)
- **Ltd Catalog #215-3571**
- JEOL 1200EX TEM using SIS III camera




Method

2d 5h 4m

1 Perform GolgiTAG-IP as previously described in dx.doi.org/10.17504/protocols.io.6qpvrjdjrogmk/v1

2 

1h

After last wash of beads with KPBS, fix the beads in  **1 mL** solution containing 4% (w/v) paraformaldehyde and 2.5% (v/v) glutaraldehyde in [M]**0.1 Molarity (M)** sodium cacodylate buffer dilute in water for  **01:00:00** at  **Room temperature** .

3 Pellet beads with magnet DynaMag™ Spin Magnet (or equivalent) and remove the supernatant.

4 

Wash beads pellet three times in  **1 mL** [M]**0.1 Molarity (M)** sodium cacodylate buffer, pelleting the beads

with magnet after each wash.




1m

After third wash, spin down beads in tabletop centrifuge at **1000 x g** for **00:01:00**. Pellet should be tightly packed.



Use a needle to gently dislodge the pellet and use a Pasteur pipette to remove the pellet to a clean glass vial. (Note: All steps are carried out in a glass vial).

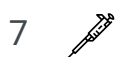

1h

To ensure rapid dehydration and embedding, cut pellet with scalpel into approximately 1mm³ pieces, then post-fix in **1 mL** 1% (w/v) Osmium tetroxide with 1.5% (w/v) sodium ferrocyanide in **0.1 Molarity (M)** sodium cacodylate buffer for **01:00:00** at **Room temperature**.



Wash pellets three times in **0.1 Molarity (M)** sodium cacodylate buffer, by adding the buffer, allowing the pellets to settle and gently taking out supernatant with P1000 tip.



1m

Wash pellets three times in **0.1 Molarity (M)** sodium cacodylate buffer for **00:01:00** (1/3).




1m

Wash pellets three times in **0.1 Molarity (M)** sodium cacodylate buffer for **00:01:00** (2/3).



1m

Wash pellets three times in **0.1 Molarity (M)** sodium cacodylate buffer for **00:01:00** (3/3).

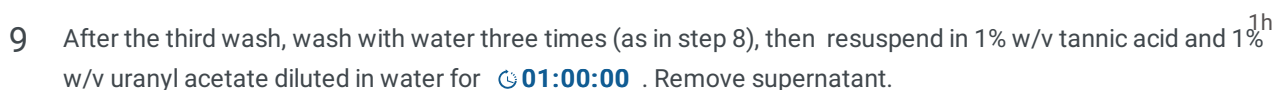

1h

10



1h 10m

Dehydrate pellet through alcohol series:

50% ethanol for ⌚ 00:10:00

70% ethanol for ⌚ 00:10:00

80% ethanol for ⌚ 00:10:00

90% ethanol for ⌚ 00:10:00

95% ethanol for ⌚ 00:10:00

then twice in 100% ethanol for ⌚ 00:10:00 (1/2)

100% ethanol for ⌚ 00:10:00 (2/2)

Use a P1000 pipette to gently remove ethanol after each series.

11 Further dehydrate in 100% propylene oxide twice for 10 min.

11.1 Further dehydrate in 100% propylene oxide twice for ⌚ 00:10:00 (1/2).

10m

11.2 Further dehydrate in 100% propylene oxide twice for ⌚ 00:10:00 (2/2).

10m

12



10m

Resuspend pellet into 50% v/v propylene oxide and 50% v/v Durcupan resin ⌚ Overnight .

13 Open glass vial to allow propylene oxide to evaporate (This will take approximately ⌚ 01:00:00).

1h

14 To embed pellet in resin, resuspend in 100% Durcupan resin.

15 Embed polymerise resin in an oven at 🔥 60 °C for ⌚ 48:00:00 .

2d

Cut sections of polymerised embedded resin with Leica UCT ultramicrotome.

16

Sections should be 70-100nm thick

17

Stain sections with 3% aqueous uranyl acetate for 00:20:00 .

20m

18

Stain with Reynold lead citrate for 00:20:00 .

20m

19



Image stained sections on JEOL 1200EX TEM using SIS III camera.

Note: Due to the cutting and many centrifugation steps, it is expected that the magnetic particles dissociate from the plastic bead core.)

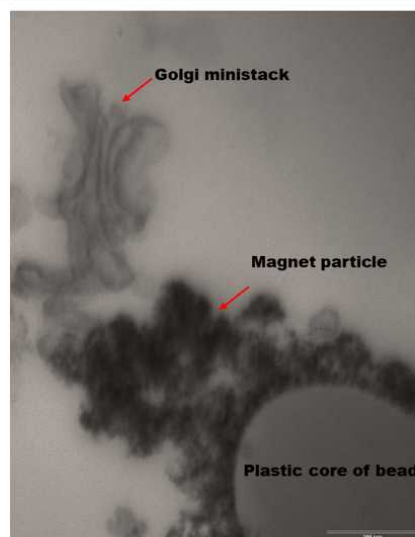


Figure 1: Transmission electron micrograph of enriched Golgi. Intact Golgi ministack is shown to be captured by GolgiTAG-IP. Scale bar 200nm