



VERSION 2

APR 28, 2023

Flaviviruses (West Nile, Zika, Dengue) NS2B/NS3 Fluorescence Dose Response V.2

Haim

Barr^{1,2}, Noa Lahav^{1,2}

¹Weizmann Institute of Science; ²ASAP Drug Discovery Consortium



ASAP Discovery

DISCLAIMER

OPEN ACCESS

Protocol Citation: Haim Barr, Noa Lahav 2023. Flaviviruses (West Nile, Zika, Dengue) NS2B/NS3 Fluorescence Dose Response. **protocols.io** <https://protocols.io/view/flaviviruses-west-nile-zika-dengue-ns2b-ns3-fluore-cthtwj6n>Version created by ASAP Discovery

MANUSCRIPT CITATION:

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Apr 28, 2023

Last Modified: Apr 28, 2023

PROTOCOL integer ID:
81171

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

ABSTRACT

This is a **functional, biochemical assay** used to identify treatments for viral infectious diseases related to viral **Flaviviridae infection**, (specifically **West Nile, Zika, and Dengue**) and targets the conserved **NS2B/NS3 protein**.

Utilizing a direct enzyme activity measurement method, the experiment was performed in a 384-well plate reading the fluorescence intensity. This assay tested the mode of action of inhibition.

It was developed at the Weizmann Institute of Science, as a part of the ASAP Drug Discovery Consortium.

Experiment Concentrations (From Stock to Assay)

A	B	C	D	E

A	B	C	D	E
Reagent	Stock	Concentration Loaded into GNF	Final Concentration in Assay Plate	Units
Substrate	10000	10	5	μM
DENV NS2B/NS3	217000	200	100	nM
ZIKV NS2B/NS3	225000	200	100	nM
WNV NS2B/NS3	222000	200	100	nM

Assay Buffer

A	B	C	D	E
Reagent	Stock	Concentration Loaded into GNF	Final Concentration in Assay Plate	Units
HEPES (pH 7.3)	20	10	10	mM
NaCl	100	50	50	mM
Glycerol	50	5	5	%
Igepal	10	0.05	0.05	%
TCEP	1000	1	1	mM

Compound Plate Design for Dose Response:

Total assay volume: 20 μL

Compounds top assay concentration: 60 μM

Dilution factor: 2

Dose response points: 12

Number of replicates: 2

Backfill with DMSO: yes

Compounds Plate Design for 2-Point Assay:

Total Assay Volume: 16 μL

Compounds Assay Concentration: 100 μM and 50 μM

Dilution Factor: 2

Dose Response Points: 2






Number of Replicates: 2

Backfill with DMSO: Yes

MATERIALS

Assay Buffer Reagents (Concentration listed are Stock Solution Concentrations)


1. [M] 20 millimolar (mM)

1.  HEPES 1M Solution pH 7.3 Fisher Scientific Catalog #AAJ16924K2
2. [M] 100 millimolar (mM)
 -  Sodium Chloride Merck MilliporeSigma (Sigma-Aldrich) Catalog #S9888
3. [M] 50 % volume
 -  Glycerol - for molecular biology, ≥99% Merck MilliporeSigma (Sigma-Aldrich) Catalog #G5516
4. [M] 10 Mass Percent
 -  IGEPAL-CA630 Merck MilliporeSigma (Sigma-Aldrich) Catalog #I3021 SIGMA-ALDRICH
5. [M] 1000 millimolar (mM)
 -  Tris(2-carboxyethyl)phosphine hydrochloride Merck MilliporeSigma (Sigma-Aldrich) Catalog #75259

*Note: all components are added fresh to the assay buffer before each experiment

Additional Reagents:

West Nile Virus (WNV) Reagents:

- [M] 222000 nanomolar (nM) WNV NS2B/NS3 Enzyme
 - WNV NS2B/NS3 was originally [M] 222000 nanomolar (nM) and was diluted to [M] 200 nanomolar (nM) with freshly made **Assay Buffer** before each experiment
- [M] 10000 nanomolar (nM) WNV Enzyme Substrate
 - Enzyme Substrate was
 -  Boc-Gly-Arg-Arg-AMC acetate salt Biosynth Catalog #FB110553
 - Substrate stock was created by dissolving the substrate in **DMSO** to create [M] 10 millimolar (mM) Substrate Stock
 - Before each experiment, the Substrate Stock was diluted again to be [M] 10 micromolar (μM) Substrate before every experiment with freshly made **Assay Buffer**

Zika (ZIKV) Virus Reagents:

- [M] 225000 nanomolar (nM) ZIKV NS2B/NS3 Enzyme
 - WNV NS2B/NS3 was originally [M] 225000 nanomolar (nM) and was diluted to [M] 200 nanomolar (nM) with freshly made **Assay Buffer** before each experiment
- [M] 10000 nanomolar (nM) ZIKV Enzyme Substrate
 - Enzyme Substrate was



Boc-Gly-Arg-Arg-AMC acetate salt Biosynth Catalog
#FB110553

- Substrate stock was created by dissolving the substrate in **DMSO** to create
[M] 10 millimolar (mM) Substrate Stock
Before each experiment, the **Substrate**
Stock was diluted again to be [M] 10 micromolar (μ M) Substrate before every
experiment with freshly made **Assay Buffer**

Dengue (DENV) Reagents:

[M] 217000 nanomolar (nM) DENV NS2B/NS3
Enzyme

- WNV NS2B/NS3 was originally [M] 217000 nanomolar (nM) and was diluted to
[M] 200 nanomolar (nM) with freshly made **Assay Buffer** before each
experiment

[M] 10000 nanomolar (nM) WNV Enzyme
Substrate

- Enzyme Substrate** was



Bz-Nle-KRR-AMC (hydrochloride) Cayman Chemical Company Catalog
#27710

- Substrate Stock** was created by dissolving the substrate in **DMSO** to create
[M] 10 millimolar (mM) Substrate Stock
Before each experiment, the **Substrate**
Stock was diluted again to be [M] 10 micromolar (μ M) Substrate before every
experiment with freshly made **Assay Buffer**

SAFETY WARNINGS



Please be sure to wear proper Personal Protective Equipment (PPE) while performing this experiment.

BEFORE START INSTRUCTIONS

Note: Inhibitor compounds stock concentration is 20 mM. Compounds are pre-dispensed into 384 plates and stored at -200°C until use.

Determine which Flavivirus is needed and prepare solutions

- Determine which Flavivirus is needed and prepare solutions based on the materials section.



A	B	C	D	E

A	B	C	D	E
Reagent	Stock	Loaded into GNF	Final in assay plate	units
Choice of NS2B/NS3 Enzyme Protein				
DENV NS2B/NS3	217000	200	100	nM
ZIKV NS2B/NS3	225000	200	100	nM
WNV NS2B/NS3	222000	200	100	nM
Choice of Viral Substrate				
WNV/ZIKV Substrate	10000	10	5	μM
DENV Substrate	10000	10	5	μM
Assay buffer				
HEPES pH=7.3	20	10	10	mM
NaCl	100	50	50	mM
Glycerol	50	5	5	%
Igepal	10	0.05	0.05	%
TCEP	1000	1	1	mM

Prepare 384-well Plate for experiment



2h 31m



2 **OPEN** the EQUIcon Software and **SELECT** the "Flavivirus dispense 7,8 C" Program

3 **PRIME** the GNF Washer/Dispenser II (GNF) with  3 mL Ehtanol and  3 mL Dionized Water




4 **CONFIRM** that the GNF had accurately dispensed Ethanol and Water


4.1 **WEIGH** the plate and **RECORD**

4.2 **DISPENSE**  3 mL Ehtanol and  3 mL Dionized Water into a plate


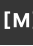
4.3 **WEIGH** the plate and **RECORD**. Determine if the GNF Washer/Dispenser II had accurately dispensed  3 g Dionized Water and  2.367 g Ethanol



5 **CONNECT Assay Buffer** to 7C and your **Flavivirus NS2B/NS3** to position 8C of the GNF Washer/Dispenser II.


5.1 **PRIME** the GNF with  300 μ L Assay Buffer and with  300 μ L  200 nanomolar (nM) Flavivirus NS2B/NS3 respectively.

6 **DISPENSE**  10 μ L Assay Buffer to columns **1 and 23** using the 7C position of the GNF



- **Note:** These columns will be the inhibitor control columns (Containing: substrate + assay buffer + DMSO, no compounds)

7 **DISPENSE**  10 μ L  200 nanomolar (nM) Flavivirus NS2B/NS3 to columns 2 through 22 and column 24 using the 8C position of the GNF.

- **Note:**  200 nanomolar (nM) Flavivirus NS2B/NS3 is two times the assay concentration. The final concentration of the Flavivirus NS2B/NS3 is  100 nanomolar (nM) during the assay.
- Columns 2 and 24 are **neutral control columns** (Contain: Enzyme + substrate + DMSO, no compounds)

8 **CENTRIFUGE**  1500 rpm, Room temperature, 00:01:00 plate to remove bubbles

1m

9 **INCUBATE** plate  02:00:00 at  Room temperature

2h



⚠ Make sure the plate is protected from light!

During Incubation: PREPARE the GNF to dispense the Flavivirus Substrate


9.1 **EMPTY** 7C of the GNF.



9.2 **WASH** 7C tubing in **Assay Buffer**. Discard used Assay buffer

9.3 **PRIME** 7C of the GNF with  300 μ L  10 micromolar (μ M) Flavivirus Substrate

10 **DISPENSE**  10 μ L  10 nanomolar (nM) Flavivirus Substrate to Columns 1 through 23 (the full plate)

- **Note:**  10 nanomolar (nM) Flavivirus Substrate is two times the assay concentration. The final concentration of the Flavivirus Substrate is  5 nanomolar (nM) during the assay.

11 **CENTRIFUGE**  1500 rpm, Room temperature, 00:01:00 plate to remove bubbles

12 **INCUBATE** plate for  00:30:00 at  Room temperature

30m

⚠ Make sure the plate is protected from light!

Recommended: Clean GNF during incubation

Read Plate Flourescence

13 **READ** and **RECORD** the plate Relative fluorescence units (RFU) via the "**Flavivirus protocol**" on the **PERAstar FS Control Software**.

Expected result

gain 300 should yield ~20,000 RFU in full reaction; 7000 RFU in Buffer control

14