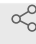





Sep 12, 2022

Intracellular neuromelanin quantification

[miquel.vila](#)¹¹Vall d'Hebron Research Institute1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.yxmvm27e6g3p/v1 [joan.compte](#)

ABSTRACT

Protocol for quantifying intracellular neuromelanin

DOI

dx.doi.org/10.17504/protocols.io.yxmvm27e6g3p/v1

PROTOCOL CITATION

miquel.vila 2022. Intracellular neuromelanin quantification. **protocols.io**
<https://protocols.io/view/intracellular-neuromelanin-quantification-cgg4ttyw>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 12, 2022

LAST MODIFIED

Sep 12, 2022

PROTOCOL INTEGER ID

69884

- 1 Scan midbrain section using the Panoramic Midi II FL, HQ Scientific 60x
- 2 Acquire SNpc images at 63x with CaseViewer

- 3 Upload individual images at Image J software
- 4 Click on Adjust canvas size and adjust it at 1596x1198
- 5 Click on Invert image
- 6 With the free hand selections tool, draft a neuromelanin-pigmented neuron (excluding the nucleus) and measure the optical density (pixel brightness) and the area
- 7 With the free hand selections tool, draft the neuromelanin pigment of the neuron and measure the optical density (pixel brightness) and the area
- 8 With the free hand selections tool, draft 15-25 non-pigmented neurons (excluding the nucleus), measure the optical density and calculate the mean value
- 9 Normalize (i.e., subtract) the values of the neuromelanin-pigmented neuron's optical density with the mean value of the optical density of the non-pigmented neurons
- 10 Extra measurement: calculate the percentage of the neuronal area occupied by the pigment dividing the neuromelanin pigment area by the neuron's area