




Jun 28, 2021

ELISA for assessing the burden of Neutrophil Extracellular Traps (NETs) in clinical serum samples.

kathryn.hally^{1,2,3}¹Department of Surgery and Anaesthesia, The University of Otago, Wellington, New Zealand;²School of Biological Sciences, Victoria University of Wellington, Wellington, New Zealand;³Wellington Cardiovascular Research Group, Wellington, New Zealand.

1 Works for me

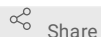

Share

dx.doi.org/10.17504/protocols.io.bmp7k5rn

The University of Otago, Wellington, New Zealand

kathryn.hally

ABSTRACT

This protocol steps through the ELISA procedure for detecting three NET-specific biomarkers: 1) myeloperoxidase-DNA (MPO-DNA) complexes, 2) neutrophil elastase-DNA (NE-DNA) complexes, and 3) citrullinated histone H3 (H3(cit)).

DOI

dx.doi.org/10.17504/protocols.io.bmp7k5rn

PROTOCOL CITATION

kathryn.hally 2021. ELISA for assessing the burden of Neutrophil Extracellular Traps (NETs) in clinical serum samples.. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bmp7k5rn>

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Hally KE, Parker OM, Brunton-O'Sullivan MM, Harding SA, Larsen PD. Linking Neutrophil Extracellular Traps and Platelet Activation: A Composite Biomarker Score for Predicting Outcomes after Acute Myocardial Infarction. *Thromb Haemost.* 2021 May 13. doi: 10.1055/s-0041-1728763. Epub ahead of print. PMID: 33984869.

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 24, 2020

LAST MODIFIED

Jun 28, 2021

PROTOCOL INTEGER ID

42463

Preparing ELISA standards.

6h

6h

1

Supernatant from PMA-stimulated neutrophils will be used as the standard for the MPO-DNA and H3(cit) ELISAs.

Carry out the following to obtain NET-rich supernatant.

Isolating granulocytes:

- From a healthy volunteer, draw **10 mL blood** into EDTA blood collection tubes.
- To two 15 mL centrifuge tubes, add **5 mL Polymorphprep** .
[Polymorphprep Axis-Shield Catalog #1114683](#)
- Carefully layer **5 mL undiluted blood** on top of each Polymorphprep layer.
- Centrifuge these tubes **650 x g, 25°C, 00:30:00 , with soft start and brakeless deceleration** .
- Following centrifugation, two cell layers will have formed which are distinct from the packed red blood cell fraction. Aspirate off and discard the first cell layer (PBMCs).
- Carefully aspirate off the remaining granulocyte layer and pool in a 50 mL centrifuge tube and add 1x PBS to **25 mL** .
- Centrifuge at **100 x g, 4°C, 00:20:00** .
- Resuspend the granulocyte pellet in **1 mL ice-cold ddH2O** **On ice** . This will lyse any contaminating red blood cells.
- After **00:00:30** , add **24 mL 1x PBS** .
- Centrifuge at **100 x g, 4°C, 00:20:00** .
- Resuspend in **1 mL cell culture media** , count and adjust cell count to 2×10^6 cells/mL. Cell culture media used here: 10% FCS, 2 mM L-glutamate, 100 U/mL penicillin, 100 µg/mL streptomycin, 0.01M HEPES buffer, 0.1% β-mercaptoethanol, 0.01 nM non-essential amino acids.

Diluting phorbol 12-myristate 13-acetate (PMA):

- PMA is provided as a lyophilized powder that is reconstituted in DMSO. Dilute PMA to a final concentration of **200 Nanomolar (nM)** using cell culture media. Ensure that the concentration of DMSO does not exceed 0.1% in this solution.

[Phorbol 12-myristate 13-acetate \(PMA\) Sigma](#)

Aldrich Catalog #P1585

PMA stimulation:

- Add isolated granulocytes with diluted PMA in a 1:1 ratio (50 µL: 50 µL) into wells within a 96-well non-tissue culture treated plate.
- Incubate for **04:00:00** at **37 °C** and **5% CO2** .
- After incubation, centrifuge at **400 x g, 25°C, 00:20:00** .
- Aspirate off the NET-rich supernatant, aliquot and store at **-80 °C** until required. Be careful not to disturb the neutrophil monolayer, only collect cell-free supernatant.

Repeat:

- Repeat this protocol for a total of 10 times. Each time, draw blood from a new healthy donor.

2 3d
For the NE-DNA ELISA, a number of unknown clinical serum samples should be screened to determine a number of samples (preferably > 5) that show the highest level of endogenous NE-DNA complexes.

The NE-DNA ELISA procedure described below can be used to determine these samples.

3 30m
Preparation and storage of NET standards:

- For MPO-DNA and H3(cit) ELISAs.

For each of 10 healthy volunteers, NET-rich supernatant from PMA-stimulated neutrophils has been aliquoted and stored. Before proceeding with NET ELISAs, thaw an aliquot from each healthy volunteer and pool.

Aliquot out pooled NET-rich supernatant in **150 µl aliquots** and store at **-80 °C** until required.

- For NE-DNA ELISA.

Having identified a number of clinical serum samples with high endogenous levels of NE-DNA complexes, thaw an aliquot from each sample and pool.

Aliquot out pooled serum in **150 µl aliquots** and store at **-80 °C** until required.

Solutions to be prepared for ELISAs. 1h 10m

4 10m
For Day 1:

- Coating buffer: **100 mL 1x PBS**

Phosphate buffered saline Sigma

Aldrich Catalog #P4417-50TAB

5 30m
For Day 2:

- Wash buffer (0.1% Tween 20 in 1X PBS): **1.5 mL Tween 20** + **1498.5 mL 1X PBS**, adjust to **pH 7.4**

Tween 20 Sigma

Aldrich Catalog #P1379

- Assay diluent (1% BSA in 1X PBS): **0.8 g BSA powder** + **80 mL 1x PBS**, filter before use

- Blocking buffer (5% BSA in 1X PBS): **4.0 g BSA powder** + **80 mL 1x PBS**, filter before use

Bovine serum albumin MP

Biomedicals Catalog #0219989680

For Day 3:

- Stop solution (1.8M H₂SO₄): **5.5 mL sulfuric acid (H₂SO₄)** + **44.5 mL 1x PBS**

Sulfuric acid (H₂SO₄) Sigma

Aldrich Catalog #258105

- Substrate buffer (0.11M acetate buffer): **1.5 g sodium acetate** + **100 mL ddH₂O**, adjust to **pH 5.5** with 100% acetic acid. Can be stored for up to two weeks refrigerated.

Sodium acetate trihydrate Sigma

Aldrich Catalog #S8625

- 6 mg/mL TMB solution: **9 mg TMB** + **1.5 mL DMSO**. Protect against light exposure. Can be stored for up to one month at room temperature.

3355-Tetramethylbenzidine (TMB) powder Sigma

Aldrich Catalog #860336

- 3% hydrogen peroxide: **10 µL 30% hydrogen peroxide** + **90 µL ddH₂O**. Can be stored for up to one week refrigerated.

- TMB substrate solution: **42 mL substrate buffer** + **700 µL TMB solution** + **42 µL 3% hydrogen peroxide**. Protect against light exposure.

ELISA day 1.

25m

- Label three 96-well ELISA plates as follows: MPO-DNA, NE-DNA, H3(cit).

Nunc MaxiSorp™ flat-bottom Thermo

Fisher Catalog #44-2404-21

- Add **200 µL coating buffer** to each well and incubate for **00:01:00**.
- Following this, flick off coating buffer and bang each plate face down on paper towels.

- Dilute each capture antibody using coating buffer, following the table below:

Type of ELISA	Capture antibody	Dilution	Antibody amount	Coating buffer amount	Final volume
MPO-DNA	MPO	1:300	35 µL	10,465 µL	10,500 µL
NE-DNA	Elastase	1:300	35 µL	10,465 µL	10,500 µL
H3(cit)	Histone	1:100	105 µL	10,395 µL	10,500 µL

[Anti-histone antibody (clone H11-4) is supplied as a part of the Cell Death Detection ELISA]

Mouse anti-human myeloperoxidase (MPO) clone 4A4 Bio-rad




Laboratories Catalog #0400-0002

☒ Mouse anti-human neutrophil elastase clone ELA10-101.5 Life

Technologies Catalog #MA1-10608

☒ Cell Death Detection






ELISA Roche Catalog #11544675001

- Add  **100 µl appropriate diluted capture antibody** to each well in the ELISA plate.
- Seal each plate with a cover slip and incubate  **Overnight** at  **4 °C**.

ELISA day 2. 4h 20m

3h 20m

9

- Thoroughly remove the capture antibody by flicking off the solution and banging each plate face down on paper towels.
- Add  **200 µl wash buffer** to each well and incubate for  **00:01:00**. Flick and bang again.
- Repeat another four times (five washes in total).
- Add  **200 µl blocking buffer** to each well, cover and incubate for  **03:00:00** at  **Room temperature**.

30m

10

- During this incubation, dilute all unknown clinical serum samples with assay diluent as follows:

A	B	C	D
ELISA	Dilution (% of final volume)	Amount of serum	Amount of assay diluent
MPO-DNA	5%	6 µL	114 µL
NE-DNA	40%	48 µL	72 µL
H3(cit)	20%	24 µL	96 µL

- Also during this incubation, dilute NET standards with assay diluent as follows:

A	B	C	D	E	F	G
ELISA	Type of standard	Dilution of Standard 1 (% of final volume)	Amount of standard	Amount of assay diluent	Serial dilution	Concentration of standards
<i>MPO-DNA</i>	Supernatant from PMA-stimulated neutrophils	5%	12 µL	228 µL	1:1	5.00, 2.50, 1.25, 0.63, 0.31, 0.16
<i>NE-DNA</i>	Pooled clinical serum samples with high endogenous levels of target antigen	50%	120 µL	120 µL	1:1	50.00, 25.00, 12.50, 6.25, 3.13, 1.56, 0.78
<i>H3(cit)</i>	Supernatant from PMA-stimulated neutrophils	40%	96 µL	144 µL	1:1	40.00, 20.00, 10.00, 5.00, 2.50, 1.25

All sample dilutions should be empirically established according to spike-and-recovery and parallelism experiments. Likewise, the concentrations for the standard should be empirically established to ensure saturation and acceptable recovery. We provide these dilutions as a starting point.

11

30m

- After three hours of blocking, wash each plate five times with wash buffer as described above.
- Add **50 µl standard** or **50 µl diluted serum** to each well, and ensure that each standard/serum sample is run in duplicate.
- Cover and incubate **Overnight** at **4 °C**.

ELISA Day 3. 6h 40m

12

2h 30m

- Dilute each detection antibody using assay diluent, following the table below:

Type of ELISA	Detection antibody	Dilution	Antibody amount	Coating buffer amount	Final volume
<i>MPO-DNA and NE-DNA</i>	DNA-Peroxidase	1:200	105 µL	20,895 µL	21,000 µL
<i>H3(cit)</i>	H3(cit)	1:4000	2.5 µL	9,997.5 µL	10,000 µL

[Anti-DNA-peroxidase (clone MCA-33) is supplied as a part of the Cell Death Detection ELISA]




All antibody dilutions should be empirically established to ensure maximal signal-to-noise ratio, but we provide these dilutions as a starting point.

[Cell Death Detection](#)

ELISA Roche Catalog #11544675001

⊗ Anti-histone H3 (citrulline R2 R8 R17)

antibody Abcam Catalog #ab5103

- Wash each plate five times with wash buffer and add  **100 µl diluted detection antibody** to each well in the ELISA plate.
- Cover and incubate for  **02:00:00** at  **Room temperature** in the dark.

13

30m

- During this incubation, make up the second detection antibody for the H3(cit) ELISA using assay diluent, as follows:

A	B	C	D	E	F
Type of ELISA	Detection antibody	Dilution	Antibody amount	Coating buffer amount	Final volume
H3(cit)	Goat anti-rabbit IgG-HRP	1:2500	4.2 µL	10,495.8 µL	10,500 µL

All antibody dilutions should be empirically established to ensure maximal signal-to-noise ratio, but we provide these dilutions as a starting point.

⊗ Goat anti-rabbit IgG HRP Bio-rad





Laboratories Catalog #1706515

- Wash each plate five times with wash buffer.

14

1h 10m








FOR MPO-DNA and NE-DNA ELISA plates:


- Add  **100 µl TMB substrate solution** to each well, cover and incubate for  **00:45:00** at  **Room temperature** in the dark.
- Stop the colour change reaction by adding  **100 µl stop solution** .
- Read the optical density at 450 nm.

15

2h 30m

For the H3(cit) ELISA plate:

- Add  **100 µl second detection antibody** to each well, cover and incubate for  **01:00:00** at  **Room temperature** in the dark.
- Wash the plate five times with wash buffer.
- Add  **100 µl TMB substrate solution** to each well, cover and incubate for  **00:45:00** at  **Room temperature** in the dark.
- Stop the colour change reaction by adding  **100 µl stop solution** .

- 
- Read the optical density at 450 nm.