



May 23, 2022

Aggregated aSyn Dot Blot assay

rabdelha¹¹California Institute of Technology

2

dx.doi.org/10.17504/protocols.io.261gen2xdg47/v1 rabdelha

This protocol details aggregated aSyn Dot Blot assay.

[432-920.docx](#)

DOI

dx.doi.org/10.17504/protocols.io.261gen2xdg47/v1

rabdelha 2022. Aggregated aSyn Dot Blot assay. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.261gen2xdg47/v1>



Aggregated aSyn Dot Blot assay, Proteinase K resistant aSyn histology

_____ protocol ,



May 12, 2022

May 23, 2022


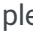
May 12,
2022 maria.sMay 17,
2022 rabdelha

62478

Aggregated aSyn Dot Blot assay

- 1 Normalize protein samples to equal concentrations between  100 ng -  1000 ng/μl in PBS.

2 

Carefully spot  **1.5 µL** each sample onto dry nitrocellulose ( **0.45 µm** pores) membrane.

3 Allow moisture to dry for a minute or two.

4 Block membrane in 5% skim milk in TBST.

5 

Wash membrane 1X in TBST.

6 Apply anti-aggregated aSyn antibody in 5% skim milk in TBST (1:1000 for Abcam fibril.aggregate specific aSyn antibody: MJFR 14-6-4-2).

7  

1h

Incubate at  **Room temperature** for 1-2 hours, or  **Overnight** at  **4 °C** .

8 

Wash 3x with TBST.

8.1 Wash with TBST (1/3).

8.2 Wash with TBST (2/3).

8.3 Wash with TBST (3/3).

9



Incubate with anti-rabbit IgG (or appropriate secondary) at 1:1000 for 1-2 hours in TBST.

10



Wash 3x with TBST.

10.1 Wash with TBST (1/3).

10.2 Wash with TBST (2/3).

10.3 Wash with TBST (3/3).

11

Develop with Biorad Chemiluminescent substrate, or similar.

Proteinase K resistant aSyn histology

1h 10m

12



Certainly can be modified for either DAB or Fluorescent staining, I have DAB protocol too

Wash 50um brain sections 1x in PBS.

13

Treat with Proteinase K (NEB, 1:4000, stock at **20 mg/mL**) for **00:10:00** at **Room temperature**.

10m

14



Wash 2x with IHC buffer (3% NGS, 0.3% TX100, in PBS).

14.1 Wash with IHC buffer (3% NGS, 0.3% TX100, in PBS) (1/2).

14.2 Wash with IHC buffer (3% NGS, 0.3% TX100, in PBS) (2/2).

15 Block with Mouse-on-Mouse kit (Vector Lab, PK-2200, just first step in instructions) for 1h
🕒 01:00:00 .

16 

Wash 2x with IHC buffer.

16.1 Wash with IHC buffer (1/2).

16.2 Wash with IHC buffer (2/2).

17  

Incubate 🕒 **Overnight** at 🌡 **4 °C** with primary antibody 1:250 – 1:1000 (anti-synuclein BD Biosci 610787).

18 

Wash 3x with PBS.

18.1 Wash with PBS (1/3).

18.2 Wash with PBS (2/3).

18.3 Wash with PBS (3/3).

19



Incubate with desired secondary (anti-Mouse IgG AF488 is my choice at 1:1000).

20



Wash 3x with PBS.

20.1 Wash with PBS (1/3).

20.2 Wash with PBS (2/3).

20.3 Wash with PBS (3/3).

21 Stain with Nissl and Mount with Profade with DAPI.