



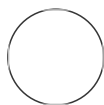
VERSION 1

NOV 14, 2022



SHARE

IN DEVELOPMENT

 Goga Lab RT-qPCR protocol: QuantStudio6 V.1**This protocol is published without a DOI.**Jeremy.williams¹¹UCSF

Jeremy.williams

COMMENTS 0

ABSTRACT

Guidelines for preparing RT-qPCR samples for QuantStudio 6 located in HSW7 lab space.

ATTACHMENTS

[GogaLab-RTqPCR-Preparation.xlsx](#)

PROTOCOL CITATION

Jeremy.williams 2022. Goga Lab RT-qPCR protocol: QuantStudio6 . **protocols.io**
<https://protocols.io/view/goga-lab-rt-qpcr-protocol-quantstudio6-ci9puh5n>



LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Nov 14, 2022

LAST MODIFIED

Nov 14, 2022

PROTOCOL INTEGER ID

72719

GUIDELINES

- *start with full tip boxes and use tips in coordination with your plate map, so you never get lost
- *watch 1uL volumes in the pipet tip like a hawk - major source of variability
- *cap and gently vortex your *mixture per primer set* every 12 replicates or so, sometimes things separate

MATERIALS TEXT


PowerUp SYBR Green Master Mix
Your generated cDNA samples
100uM single primer stocks

BEFORE STARTING

*Thaw on ice (leave time!) and keep all reagents on ice through preparation. Prepare the plate on ice.

ATTACHMENTS

[GogaLab-RTqPCR-Preparation.xlsx](#)

- 1  GogaLab-RTqPCR-Preparation.xlsx Dilute stock IDT primers to 100uM
- 2 Mix your forward and reverse primer pairs together, to a final dilution of 10uM forward and 10uM reverse. For example, add 10uL each of forward and reverse primers to 80uL PCR-quality DI for 100uL final volume.
- 3 Mix *total reagent volumes required* (see excel spreadsheet, green) for DI and PowerUP

- 4 Mix *mixtures per primer set* (blue) volumes together.
- 5 Pipet 1uL diluted cDNA into respective qPCR well, aiming for the sidewall of each respective well.
- 6 Add 19uL *mixture per primer set* into each respective well.
- 7 Seal plate, and spin down 1000rpm for 1 minute. Use bacterial, not tissue culture, centrifuge.
- 8 Load plate into centrifuge and proceed using QuantStudio software suite.