

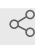


Aug 01, 2022

IB BIO SL IA M23 - A study on the conditions that affect yeast fermentation

Marcy Roca Cugat¹¹Universitat de Girona

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.261geny6jg47/v1

Marcy Roca Cugat

Universitat de Girona, Universitat Autònoma de Barcelona, In...

ABSTRACT

This protocol is designed to aid me in the research for my International Baccalaureate Biology Internal Assignment

DOI

dx.doi.org/10.17504/protocols.io.261geny6jg47/v1

PROTOCOL CITATION

Marcy Roca Cugat 2022. IB BIO SL IA M23 - A study on the conditions that affect yeast fermentation . **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.261geny6jg47/v1>



KEYWORDS

biology, fermentation, yeast, ib, ia

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 01, 2022

LAST MODIFIED

Aug 01, 2022

PROTOCOL INTEGER ID

68006

GUIDELINES

None

MATERIALS TEXT

- Flask (3 mouths)
- Thermometer
- Holed plugs x2
- Glass tube
- Rubber tube
- Gas recovery flask
- Basin

- Water
- Sugar
- Yeast

SAFETY WARNINGS

None

Preparation

5m 30s

1 Check that all reagents and materials are ready.

30s

2 Assemble the instrumentation

5m

2.1 One end of the flask must be plugged with a holed stopper and a thermometer.
The other end should be plugged with a holed stopper and a glass tube.
The middle opening should be covered with plastic and secured with a rubber band.

2.2 A long and shallow recipient must be filled with water as much as possible.
The gas recovery flask should also be filled with water as much as possible.
Then, the latter should be inverted and suspended inside the former.

2.3 The tube should be connected from the glass tube to the gas recovery flask.

3 The block of compressed yeast must be cut in 4 parts

- 4 Several factors must be checked and accounted for: pressure, temperature and initial weight of the full gas recovery flask.
- 5 One quarter of the yeast block must be dissolved in sugary water (concentration undetermined, volume undetermined, temperature undetermined). This should then be placed in the flask.
- 6 Wait for fermentation: 🕒 **00:30:00** 30m
- 7 Weigh the gas recovery flask. The difference from the first measurement is the gas generated. Convert it to moles of $\text{CO}_2 + \text{H}_2\text{O}$ for a universal measurement