



DEC 12, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.eq2lyn25pvx9/v1](https://dx.doi.org/10.17504/protocols.io.eq2lyn25pvx9/v1)

**Protocol Citation:** Jinho Lee, Jose Ceja Navarro, Koei Chin, Heidi S Feiler, Christopher Corless 2023. NanoString GeoMx DSP TMA-TNP protein assay.  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.eq2lyn25pvx9/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jun 02, 2022


## 🌐 NanoString GeoMx DSP TMA-TNP protein assay

Heidi S

Jinho Lee<sup>1</sup>, Jose Ceja Navarro<sup>2</sup>, Koei Chin<sup>1</sup>, Feiler<sup>1</sup>, Christopher Corless<sup>2</sup>

<sup>1</sup>Knight Cancer Institute, OHSU; <sup>2</sup>Knight Diagnostics Laboratory, OHSU

NCIHTAN

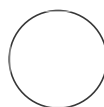

Jinho Lee

### ABSTRACT

This protocol outlines the NanoString GeoMx DSP protein assay that was applied in the Human Tumor Atlas Network (HTAN) Tissue MicroArray (TMA) -TransNetwork Project (TNP).

The TMA-TNP evaluates various characterization and analytics methodologies on a large array of breast tumor samples representing a broad spectrum of disease state and subtype. A commercially available anonymized breast tumor TMA was purchased and serial sections were distributed. Participating HTAN Centers characterized the FFPE specimens using various imaging platforms and generated a spatially resolved cell type/state census using each center's method of choice. Data was recorded in a common repository to enable joint analysis.

The protocol that immediately precedes this one for TMA-TNP Phase 3 (when this NanoString GeoMx DSP protein assay protocol was first tested) can be found at: [dx.doi.org/10.17504/protocols.io.6qpvr6992vmk/v1](https://dx.doi.org/10.17504/protocols.io.6qpvr6992vmk/v1). It describes FFPE block serial sectioning, slide processing and TMA sample distribution.

In this protocol, the DSP protein assay was performed with the following human marker panels: Immune cell profiling core, Immune activation status, Immune-Oncology (IO) Drug target, Cell death, Pan-tumor, MAPK and PI3K/AKT panels including OHSU custom panel (Cell cycle and DNA damage). A total of 87 protein targets were evaluated in the TMA samples. Three compartments (Immune, Tumor and Stroma) in each TMA core were analyzed to determine cell-to-cell interactions in the tissues.

**Last Modified:** Dec 12, 2023

**PROTOCOL integer ID:** 63771

**Keywords:** NanoString, GeoMx, Digital Spatial Profiler, HTAN, TMA, DSP, cell profiling, cell type

## MATERIALS

### 1. Materials

Slide Baking Oven  
Antigen retrieval pressure cooker  
Humidity chamber  
Thermocycler  
Centrifuge (up to 2,000g)

### 2. Reagents

Xylene  
Ethanol  
TBS  
TBS-T  
Citrate buffer pH6.0  
10% Neutral Buffered Formalin (NBF)  
DEPC-treated water

### 3. Nanostring reagents

GeoMx Instrument Buffer kit / Collection plate (Item no. 100474)  
FFPE slide prep kit (Buffer W / Buffer S) (Item no. 121300312)

Human Immune cell profiling core panel (Item no. 121300101)  
Human Immune activation status panel (Item no. 121300103)  
Human IO Drug target panel (Item no. 121300102)  
Human Cell death panel (Item no. 121300112)  
Human MAPK panel (Item no. 121300111)  
Human PI3K/AKT panel (Item no. 121300113)  
Human Pan-tumor panel (Item no. 121300105)  
OHSU SQ-35430 custom protein panel

Solid tumor Morphology kit (Anti-PanCK / Anti-CD45 / SYTO-13) (Item no. 121300301)

Hybridization CodeSets (A-H) (Item no. 121300401)

Master kit / Hybridization buffer (Item no. 100052)

## FFPE slide sample preparation

1h 59m

### 1 FFPE slide preparation




- 1.1 Bake FFPE slides at  $60\text{ }^{\circ}\text{C}$  for 01:00:00 . 1h
- 1.2 Deparaffinize by sequential incubation with Xylene (00:03:00 twice), 100% EtOH (00:01:00), 95% EtOH (00:01:00) and 70% EtOH (00:01:00). 9m
- 1.3 Briefly rinse the slides with  $\text{diH}_2\text{O}$  (distilled water) and incubate with 1x Citrate antigen retrieval buffer (Epitope retrieval buffer, pH 6.0) for 00:15:00 at high pressure in a pressure cooker. 15m
- 1.4 Take out the slides from the pressure cooker and let them cool down to room temperature around 00:25:00 . 25m
- 1.5 Wash the slides 5 times with TBS-T for 00:02:00 each. 10m

## Antibody panel incubation

1h 59m

### 2 DSP protein panel and visualization marker incubation






- 2.1 Mark off the entire tissue section with a hydrophobic pen to create a reagent boundary.

- 2.2 Block the slides for  01:00:00 at room temperature using the Nanostring blocking buffer (Buffer W) 1h
- 2.3 Freshly prepare each commercial protein panel in 1:25 dilution and each visualization marker (PanCK and CD45) in 1:40 dilution with Buffer W to a volume of 220ul per slide.
- For the custom antibody panel, the final concentration at 250ng/ml per antibody is calculated and adjusted to the final volume of 220ul.
- 2.4 Incubate slides with all antibody panels and visualization markers for  Overnight at  4 °C 16h

## Post fixation and DNA staining

1h 18m

### 3 Post fixation

- 3.1 After antibody incubation, wash the slides 3 times with TBS-T for  00:03:00 each. 9m
- 3.2 Fix the slides with 10% NBF (Neutral buffered Formalin) for  00:30:00 at RT. 30m
- 3.3 Wash 3 times with TBS-T for  00:03:00 each, then stain the samples with SYTO-13, prepared 1:10 dilution in TBS for  00:15:00 at RT. 18m
- 3.4 Slides are briefly washed 2 times with TBS-T for  00:03:00 each. 6m

## Loading samples to DSP

### 4 DSP running preparation

- 4.1 Scrape the hydrophobic barrier off with a scalpel or a straight-edged razor.
- 4.2 Place slides on the GeoMx slide tray and clean the back of the slide with 70% EtOH.
- 4.3 Once the gasket is sealed, place 6 mL of Buffer S on slides.
- 4.4 Load the GeoMx slide tray onto DSP platform.

## GeoMx operation

6h 30m

### 5 Slide scan, ROI (Region of Interest) selection and AOI (Area of Interest) segmentation.

- 5.1 Log onto GeoMx software and start with **"New / Continue Run"**.
- 5.2 After loading slides with the collection plate information, DSP is ready to scan slides.

30m

- 5.3** A slide scan name is created and panel/visualization marker information is selected as below: **10m**
- In the **Probe Reagent Kit** field, in any order select **Human Immune Cell Profiling Protein Core, Human Immune Activation Status Protein, Human IO Drug Target Protein, Human Cell Death Protein, Human MAPK Signaling Protein, Human Pan-Tumor Protein, Human PI3K/AKT Signaling Protein, and SQ-35430 OHSU Custom**
1. Select the FITC/525 nm, Cy5/568 nm and Texas Red/615 nm channels.
  2. For **FITC/525 nm**, select **SYTO 13** as fluorophore, DNA as biological target and enter 50 as exposure time.
  3. For **Cy5/568 nm**, select **Alexa 532** as fluorophore, **PanCK** as biological target and enter 300 as exposure time.
  4. For **Texas Red/615 nm**, select **Alexa 594** as fluorophore, **CD45** as biological target and enter 300 as exposure time.
  5. Select FITC/525 nm as focus channel
- 5.4** When the scan area for each slide has been adjusted with sensitivity setting, select **Scan**. **30m**
- 6** ROI selection and AOI segmentation
- 6.1** After scanning is done, each color channel intensity is adjusted to show visualization markers also with tissue or cell line property. **10m**
- 6.2** Each ROI is determined and selected by pathologist's guide, and drawn with circle (maximum 660 radius), rectangle (maximum 660x785um) or polygonal shape (maximum 660x785um). **30m**
- TNP TMA slide contains total 88 cores.
- Due to the limitation of scan area in the slide loading slot, only 44 cores per slide can be scanned and collected.
- Two TMA slides (88 cores were embedded and slightly shifted to either left or right side of slide to cover half of 88 cores in each slide) are required to collect all cores.
- 6.3** In segment menu, 3 segmentation classes (Immune, Tumor and Stroma) are added and parameters are set in the following order: **10m**

For Immune segmentation, Alexa 594 (CD45) is set to positive ("+") and others set to ignore ("0") for the first collection

For Tumor segmentation, Alexa 532 (PanCK) is set to positive ("+") and others set to ignore ("0") for the second collection

For Stroma segmentation, FITC 525 (SYTO 13) is set to positive ("+") and others are set to ignore ("0") for the third collection

Then click **Generate Segments**.

- 6.4** Once all segments are automatically generated, each channel parameter needs to be manually adjusted with pathologist's input to confirm if the segmentation is correctly done. **1h 30m**

**Caution:** less than 20 cells in each segment is removed from collection due to threshold for low signal.

- 6.5** Once all AOI segmentation is complete, Exit Scan Workspace button icon is clicked to approve ROI selection and samples are collected in 96-well plate. **3h**

## Hybridization

**18h 10m**

### 7 Hybridization

- 7.1** After completion of sample collection, the 96-well plate is finalized, removed from DSP and transferred to the PCR thermocycler. **1h**  
DNA oligo samples in 96-well plate are completely dehydrated at **60 °C** for **01:00:00** in PCR thermocycler.

**Caution:** the plate lid of PCR thermocycler should be opened completely to avoid the contamination of DNA oligo from evaporation during this dry step.

- 7.2** Add **10 µL** of DEPC-treated water to each well and resuspend the DNA oligo samples for **00:30:00** at RT. **30m**

- 7.3** (1) Add and mix well **12 µL** of each Probe A per panel to a 1.5ml Eppendorf tube (labeled Probe A tube) **10m**

Since 7 commercial panels and 1 custom panel (consists of 2 Probe A ) are used, the total volume is  $\text{108 } \mu\text{L}$  ( $= 12\text{ul} \times 9 \text{ Probe A}$ ) in the Probe A tube.

(2) Add and mix well  $\text{12 } \mu\text{L}$  of Probe B (Universal probe) with  $\text{87 } \mu\text{L}$  of DEPC water in a fresh 1.5ml Eppendorf tube (labeled Probe B tube, total volume is  $\text{99 } \mu\text{L}$  ).

(3) Calculate and add  $\text{80 } \mu\text{L}$  of hybridization buffer per row into a fresh 1.5ml Eppendorf tube (labeled Master mix tube). For 44 core samples, we collected a total 8 rows in the collection plate so total ( $= 80 \times 8 \text{ rows}$ )  $\text{640 } \mu\text{L}$  of hybridization buffer was added.

The Probe/Master Mix is made by adding  $\text{64 } \mu\text{L}$  from the Probe A tube and  $\text{64 } \mu\text{L}$  from the Probe B into the Master mix tube ( $64\text{ul} + 64\text{ul} + 640\text{ul}$ ).

**7.4** The Hybridization CodeSet for each row (A-H) is taken out from  $\text{-80 } ^\circ\text{C}$  and thawed at RT 10m  
 $\text{84 } \mu\text{L}$  of Probe / Master mix is added to each hybridization Codeset and mixed well.

**7.5** Prepare the hybridization plate separately (96-well plate) 20m  
 $\text{8 } \mu\text{L}$  of hybridization CodeSet mix was added to each well per each row in the hybridization plate.  
 $\text{3 } \mu\text{L}$  of DNA oligo samples from the rehydrated plate (out of 10ul) was mixed into each hybridization Codeset ( $8\text{ul} + 3\text{ul} = \text{total } 11\text{ul} / \text{well}$ ).

**7.6** Seal the hybridization plate with aluminum foil using the microplate heat sealer and briefly spin it down in a centrifuge at 2000g.

**7.7** Samples are hybridized at  $\text{67 } ^\circ\text{C}$  for  $\text{16:00:00}$  in PCR Thermocycler and then kept at  $\text{4 } ^\circ\text{C}$  until loading in **nCounter MAX system**. 16h

## nCounter reading

5h 45m

### 8 nCounter preparation and reading



- 8.1** In the DSP server, the collection plate is finalized and the library preparation file for nCounter load is downloaded. 15m  
The sample loading volume is determined (it will vary) according to the area size calculation of total ROI/AOI in each column.  
The Samples from each row (A-H) are collected and transferred into a 12-well strip PCR tube using a 12-channel multi-pipette.
- 8.2** Briefly, the samples are spun down in the 12-well PCR strip and then transferred to Nanostring's MA Prep-station to load samples onto the cartridge with the standard sensitivity setting. 3h
- 8.3** The cartridge is transferred to nCounter to read the counts with defined **CDF (CodeSet Design)** setting downloaded from the DSP plate information. 2h 30m
- 8.4** The **RCC (Reporter Code Count) file** from nCounter reading is downloaded to your PC and imported into the DSP server.

## DSP data analysis

### 9 QC DSP data and analysis

- 9.1** Select and queue the slides to analyze using **"New Analysis"** in the DSP server.
- 9.2** Determine the New Analysis file name and save it in the designated folder.
- 9.3** Open an analysis file and perform the **QC** with preset parameters.

**9.4** QC passed samples are processed and the QC file (CSV) is downloaded to a PC for further analysis.

**Comment:** All documents related with GeoMx DSP run can be found at  
<https://university.nanostring.com/page/document-library>