



VERSION 2

NOV 08, 2023

OPEN ACCESS



Protocol Citation: Elisabeth Rebboah 2023. Protocol to isolate and fix nuclei from flash frozen mouse liver for IGVF. [protocols.io](https://protocols.io/view/protocol-to-isolate-and-fix-nuclei-from-flash-froz-c4puyvnnw) <https://protocols.io/view/protocol-to-isolate-and-fix-nuclei-from-flash-froz-c4puyvnnw> Version created by Elisabeth Rebboah

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working


Created: Nov 07, 2023

Last Modified: Nov 08, 2023

Protocol to isolate and fix nuclei from flash frozen mouse liver for IGVF V.2

Elisabeth Rebboah¹

¹University of California, Irvine

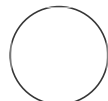

Elisabeth Rebboah
University of California, Irvine

ABSTRACT

This protocol describes isolation of nuclei from 10 week old mouse liver (tissue ID: 05) from 8 founder strains (B6J, AJ, 129S1J, NZOJ, WSBJ, NODJ, PWKJ, and CASTJ), preparation of a single nucleus suspension, and fixation for 1. single nucleus RNA-seq using the Parse Biosciences protocol (Split-seq) and 2. single nucleus RNA-seq + ATAC-seq using the SHARE-seq protocol. We process 1 rep from each strain per day; e.g. female rep 1 across all 8 strains. For 8 samples, this protocol takes about 4 hours from start to finish.

The results are 2 aliquots of fixed single-nucleus suspensions for Parse per each of the 8 samples at $\geq 2,500$ nuclei/ul, and 1 fixed nuclei pellet pooled across all 8 strains for SHARE-seq, all stored at -80C.

The first part of the protocol describes tissue sectioning, lysis, and nuclei extraction using Miltenyi Biotec's gentleMACS Octo Dissociator with accessories. Due to the size of whole mouse liver, after weighing the tissue it is chopped into ~500 mg pieces while frozen. When nuclei are extracted and counted, we determine whether we have enough to fix for Split-seq and SHARE-seq and set aside 4 million and 1 million, respectively. Ideally, the second and third parts of this protocol are performed in parallel by at least two technicians to save time. The second part describes nuclei fixation using Parse Biosciences Evercode Nuclei Fixation Kit with v2 reagents (see attachment for original version). The third part describes nuclei fixation using a modified version of the SHARE-seq fixation protocol (see attachment for original version). Any remaining nuclei are flash-frozen as a dry pellet and stored at -80C.

ATTACHMENTS

[SO+10122022_Evercode+ 2022_07_15_GRO_nuclei_Fixation+v2.0.2+User+Manual.pdf](#) [2022_07_15_GRO_nuclei_prep_combo.docx](#)

Keywords: Fixation, Nuclei fixation, Split-seq, SHARE-seq, Evercode, snRNA-seq, Parse Biosciences, Nuclei isolation, UCI, Mortazavi, IGVF, Mouse, Liver

GUIDELINES

- We recommend using a 5 ml pipette for aspirations and resuspensions > 1 ml.
- We cut the frozen livers ahead of time on dry ice to make lysis easier.
- Tips to get frozen chopped tissues out of tubes: **Use forceps** or a 1 ml tip to break up frozen tissue pieces so they can fit on the sides of the C tube. Tap tubes firmly upside-down on bench. Flick/shake open tube over the C tube with lysis buffer. Flick/shake C tube with liver pieces so that they are at least touching the lysis buffer.
- Record everything in the [IGVF spreadsheet](#), "Samples into experiment" tab.
- When possible, after nuclei isolation and during the first round of counting, 2 technicians should continue on with Parse fixation after establishing the volume needed for roughly 4 million nuclei per sample (4 samples processed per technician). The remaining technician should determine the exact volume needed for 1 million cells and proceed with SHARE-seq fixation. Parallel fixation saves about an hour of time.

MATERIALS

Name	Manufacturer	Cat. #
Nuclei Fixation Kit v2	Parse Biosciences	ECF2003
Nuclei Extraction Buffer	Miltenyi Biotec	130-128-024
RNase Inhibitor, murine	New England Biolabs	M0314L
PBS	HyClone	SH30256.02
7.5% BSA	Life Technologies	15260037
1 M HEPES pH 7.3	Sigma	H0887-100ml
NaCl	Fisher	BP358-1
MgCl ₂	Fisher	AA12315A7
Tween-20	Fisher	BP337-500
5% digitonin	Promega	G944A
Enzymatics RI	Enzymatics	Y9240L
SUPERase RI	Invitrogen	AM2696
Yeast tRNA	Invitrogen	AM7119
Glycine	Fisher	BP381-500

Name	Manufacturer	Cat. #
1M Tris pH 8.0	Thermo	AM9855G
Formaldehyde (methanol-free)	EMS	15710
gentleMACS C Tube	Miltenyi Biotec	130-093-237
gentleMACS Octo Dissociator	Miltenyi Biotec	130-095-937
MACS SmartStrainers (70 µm)	Miltenyi Biotec	130-110-916
MACS SmartStrainers (30 µm)	Miltenyi Biotec	130-098-458
NucBlue Fixed Cell ReadyProbes	Thermo Fisher	R37606
Hemocytometer	Fisher Scientific	02-671-51B
Mr. Frosty	Sigma-Aldrich	635639

Reagents/equipment, manufacturer and catalog number

Name	Reagent	Volume (for 8 samples)	Final concentration
1% BSA-DEPC	BSA	1 g	1%
	DEPC water	100 ml	
Lysis buffer	Nuclei Extraction Buffer	50 ml	
	40 U/ul RNase inhibitor	250 ul	0.2 U/ul
NB-BSA + RNase inhibitor (Make 2 aliquots)	Nuclei Buffer (Parse Biosciences)	3.15 ml	NA
	7.5% BSA	350 ul	0.75%
	RNase inhibitor (Parse Biosciences)	44.1 ul	
RSB (Make 2 aliquots of 35	PBS	35 ml	
	7.5% BSA	467 ul	0.1%

Name	ml = 70 ml for 8 samples)	Reagent	Volume (for 8 samples)	Final concentration
		RNase inhibitor	175 ul	0.2 U/ul
		1 M HEPES pH 7.3	150 ul	10 mM
		5 M NaCl	30 ul	10 mM
		1 M MgCl ₂	45 ul	3 mM
		10% Tween-20	150 ul	0.1%
		H ₂ O	14.625 ml	
	SHARE-RSB	7.5% BSA	80.26 ul	0.04%
		5% digitonin	30 ul	0.01%
		Enzymatics RI	37.5 ul	0.1 U/ul
		SUPERase RI	18.75 ul	0.025 U/ul
		Yeast tRNA	150 ul	100 ug/ml

Buffers

Setup

- 1 Coat SHARE-seq nuclei prep tubes with BSA. Fill 8 1.5 ml tubes with **1.5 ml 1% BSA** in H₂O and incubate for **30 minutes**. After incubation, aspirate BSA solution and dry for **30 minutes**. Store at 4C.
- 2 Label tubes.
- 3 Pre-chill centrifuge to **4C**.
- 4 Prepare ice buckets.

- 5 Prepare **50 ml lysis buffer** in a 50 ml conical tube on ice. Distribute **3.5 ml** into 8 gentleMACS C Tubes on ice. For NZO mice, use **4 ml** lysis buffer. Add **250 ul RNase inhibitor** to the lysis buffer aliquot the day of the experiment.
- 6 Prepare **70 ml RSB** in 2 50 ml conical tubes on ice, **35 ml per tube**. Add **175 ul RNase inhibitor** to each tube the day of the experiment.
- 7 Prepare **2 aliquots of 3.5 ml NB + BSA**. Add **44.1 ul RNase inhibitor** included in Parse Biosciences fixation kit the day of the experiment to each aliquot.
- 8 Prepare **2.5 ml nuclei buffer + RNase inhibitor** for final resuspension. Add **31.5 ul RNase inhibitor** to **2.5 ml nuclei buffer**.
- 9 Prepare **15 ml SHARE-RSB** in a 50 ml conical tube at room temperature. To SHARE-RSB, add **30 ul digitonin**, **37.5 ul Enzymatics RI**, **18.75 ul SUPERase RI**, and **150 ul yeast tRNA** fresh the day of the experiment.
- 10 Thaw components of **2 Parse Biosciences Nuclei Fixation v2** kits at room temperature, then place on ice.
- 11 Distribute **20 ul** NucBlue Fixed Cell ReadyProbes into 16 PCR strip tubes for cell counting. Need 8 tubes for counting after nuclei extraction, and another 8 tubes for final fixed nuclei.

Tissue sectioning

- 12 Keep flash frozen tissue samples on dry ice.

- 13 Tilt frozen tissues into a plastic petri dish on dry ice.
- 14 Using a clean razor blade, roughly chop tissue into pieces **~500 mg**. Not necessary to chop tissues < 500 mg.
- 15 Using clean forceps, move chopped, frozen tissue back to the original 2 ml tube. Tissue should never be thawed during this process.

Tissue lysis and nuclei extraction

- 16 Keep chopped flash frozen tissue samples on dry ice until lysis.
- 17 Drop whole frozen tissue into a chilled gentleMACS C Tube with **3.5 ml** lysis buffer. For NZO mice, use **4 ml** lysis buffer. Close tubes firmly and invert immediately, ensuring tissue is not stuck to the bottom or side. Keep tubes on ice and proceed immediately to dissociation.
- 18 Run the gentleMACS Program **4C_nuclei_1** on the Octo Dissociator (**~5 minutes**).
- 19 Remove tubes, ensuring tissue did not get stuck on the sides, and spin down in a **4C** centrifuge for **~10 seconds** to bring liquid to the bottom, then place tubes back on ice.
- 20 Filter nuclei suspension through **70 um** MACS SmartStrainer into a 5 ml tube. Fit a tube rack in ice for extra stability while filtering.
- 21 Wash **70 um** MACS SmartStrainer with **2.5 ml additional lysis buffer**. Add 2.5 ml to C tubes, cap, and swish to recover any nuclei stuck to the sides and cap of the C tubes, then wash the strainer.

- 22 Discard strainer and centrifuge the 6 ml nuclei suspension at **4C, 350g** for **5 minutes**.
- 23 Discard supernatant and resuspend nuclei pellet in **3.5 ml RSB**.
- 24 Filter nuclei suspension through **30 um** MACS SmartStrainer into a 5 ml tube.
- 25 Dilute some nuclei **1:10** by adding **500 ul nuclei** to **4.5 ml RSB** in a new 5 mL tube. This should help the concentration reach around 4 million per ml.
- 26 Count 1:10 diluted nuclei. Use **1:5** dilution factor, **5 ul + 20 ul dye**. Final dilution factor = 1:50.

Parse nuclei fixation

- 27 Set aside **4 million nuclei in RSB** in a new 5 ml tube and spin down at **4C, 350g** for **5 minutes**.
- 28 Remove supernatant and resuspend nuclei in **750 ul NB-BSA + RNase inhibitor** and filter through a **40 um** strainer (provided in Parse Biosciences kit) into a new 5 ml tube.
- 29 Add **250 uL Nuclei Fixation Solution** and mix 3 times. Do not over-mix.

- 30** Incubate nuclei for **10 minutes** on ice. Set 1 P200 pipette to 80 ul and keep the P1000 at 250 ul.
- 31** Add **80 uL Nuclei Permeabilization Solution** and mix by pipetting 3 times with the P1000 still set to 250 uL. Do not over-mix.
- 32** Incubate **3 minutes** with nuclei on ice.
- 33** Add **4 ml Nuclei Neutralization Solution** and invert the tube once to mix.
- 34** Centrifuge at **4C, 750g** for **10 minutes**.
- 35** Aspirate and discard supernatant.
- 36** Resuspend the samples in **300 ul Nuclei Buffer with RNase inhibitor** without BSA and move through a 40 um filter into a labeled 1.5 ml tube.
- 37** Count nuclei. Use **1:11** dilution factor, e.g. **2 ul + 20 ul dye**.
- 38** Add Nuclei DMSO: **5 uL** and gently flick tubes to mix. One minute later, add another **5 uL** and flick to mix, then after another minute add a final **5 uL** for a total volume of **15 uL**. Mix by gently pipetting 5x with a P200 set to 150 ul.

39 Split nuclei suspension into 2 labeled tubes, **150 ul** per tube.

40 Place tubes in a Mr. Frosty at **-80C**. The next day, move tubes to boxes in -80C racks.



SHARE-seq nuclei fixation

41 Set aside **1 million nuclei** for each of the 8 samples in RSB and spin down at **4C, 750g** for **5 minutes**.

42 Remove supernatant and resuspend nuclei pellet in **1 ml room temperature SHARE-RSB**. Transfer tube to a room temperature rack.

43 At RT, add **13.34 ul of methanol-free formaldehyde** (16% stock solution). Final concentration for nuclei: 0.2%. Close tube and nutate cells at **RT** for **5 minutes**.

44 To quench fixation, per reaction, add **56.1 ul fresh 2.5M Glycine** (0.94g per 5 ml stock), **50 ul of 1M Tris pH 8.0**, **13.3ul of 7.5% BSA**, and mix using a pipette. Incubate on ice for **10 minutes**.

45 Spin **750g, 4C, 5 minutes**. Gently remove supernatant.

46 Add **200 ul of SHARE-RSB** and gently resuspend pellet. Store on ice until all samples are completed.

47 **Pool** 200 ul of resuspended nuclei from all 8 founders into 1 labeled 2 ml tube.

48

Spin **1,000g, 4C, 10 minutes**. Gently remove supernatant. Remove all fluid and freeze at **-80C** as a **dry pellet**.



Storage of leftover nuclei

49

Move remaining nuclei in RSB on ice to labeled 2 ml tubes.

50

Spin **750g, 4C, 5 minutes**.

51

Remove all supernatant and flash-freeze nuclei as a dry pellet in liquid nitrogen. Store at **-80C**.

