




Jun 28, 2021

# The method to evaluate antibacterial effect of moth-eye films

Kiyoshi Minoura<sup>1</sup><sup>1</sup>Next Generation Technology Development unit, Development group, Sharp Display Technology Corporation

1 Works for me

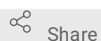

Share

[dx.doi.org/10.17504/protocols.io.ipkcdkw](https://dx.doi.org/10.17504/protocols.io.ipkcdkw)

Kiyoshi Minoura

## ABSTRACT

This method used a model designed to represent its practical application.

## DOI

[dx.doi.org/10.17504/protocols.io.ipkcdkw](https://dx.doi.org/10.17504/protocols.io.ipkcdkw)

## PROTOCOL CITATION

Kiyoshi Minoura 2021. The method to evaluate antibacterial effect of moth-eye films. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.ipkcdkw>

## LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Jun 28, 2017

## LAST MODIFIED

Jun 28, 2021

## PROTOCOL INTEGER ID

6604

### Preparation of microbial suspension

- 1 Staphylococcus aureus (S. aureus) NBRC12732 strain and Escherichia coli (E. coli) NBRC3972 strain were kindly supplied by the National Institute of Technology and Evaluation (Tokyo, Japan).
- 2 S. aureus and E. coli on nutrient agar medium (NA) plate (Nissui Pharmaceutical Co. Ltd., Tokyo, Japan) was cultured overnight at 37°C.
- 3 One colony was selected, inoculated on a new NA plate and cultured again overnight at 37°C
- 4 The bacteria from the second culture were then suspended in 1/500 nutrient broth (1/500 NB; Eiken Chemical Co. Ltd., Tochigi, Japan) at McFarland Standard 1.0, which contained around 10<sup>8</sup> cfu/ml bacterial cells.

The S. aureus or E. coli suspension was diluted 100 times with 1/500NB to a final bacterial density ranging from


5  $2.5 \times 10^6$  to  $10 \times 10^6$  cfu/ml.

#### Preparation of sample films


- 6 The moth-eye stamper was made from aluminum which was anodized and etched. The surface was fabricated into arrays of nano structures whose pitch and height were both approximately 200nm.
- 7 The resin that composed of urethane trifunctional acrylate having polyethylene glycol (PEG) chains as spacers was spread on base substrate film with 75  $\mu\text{m}$  poly ethylene terephthalate thickness.
- 8 The substrate film with the resin was laminated with a nip roll on the moth-eye stamper.
- 9 The resin was exposed to ultraviolet light, and was peeled off from the moth-eye stamper by hand after the completion of curing. The peeled film was used as sample A.
- 10 As for flat surface films (Sample B), the moth-eye stamper was substituted to the flat glass stamper. The film was fabricated with the same manner with step 9. The peeled film was used as sample B.
- 11 PET film was used as a control film.
- 12 All sample films were cut out in 5cm square.
- 13 All sample films were disinfected with ethanol.

#### Antibacterial assay of droplet model


- 14 The sample film was placed in the center of the petri dish.
- 15 Ten  $\mu\text{l}$  of the final microbial solution was dropped onto a sample film and kept in a safety cabinet at room temperature for the indicated time.

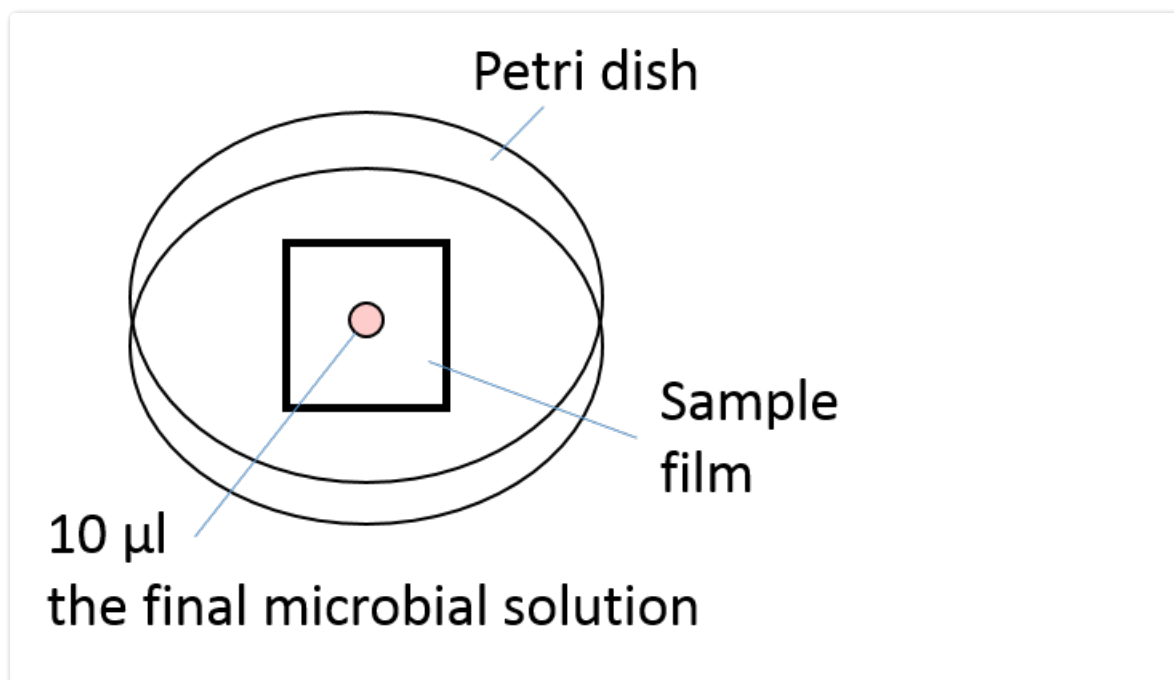

- 16 The medium consisting of 9.6 ml of fresh soybean-casein digest broth with lecithin and polysorbate 80 (SCDLP; Wako Pure Chemical Industries Ltd., Osaka, Japan) was poured directly onto the sample film in the petri dish, and the bacteria attached to the film surface were washed out by shaking the petri dish for 30 sec at 500 rpm/min using a vortex mixer.

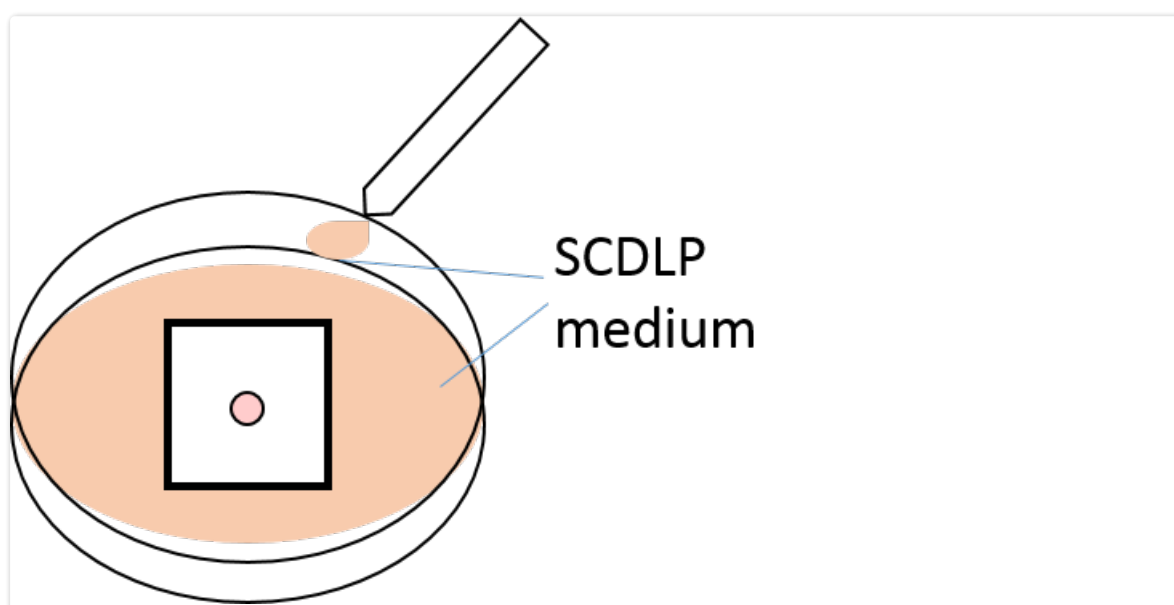

- 17 The sample film was washed out again with the SCDLP medium collected from the petri dish. This process was repeated three times.
- 18 All the SCDLP medium in petri dish, which was around 10ml, was collected in the 50mL conical centrifuge tube.

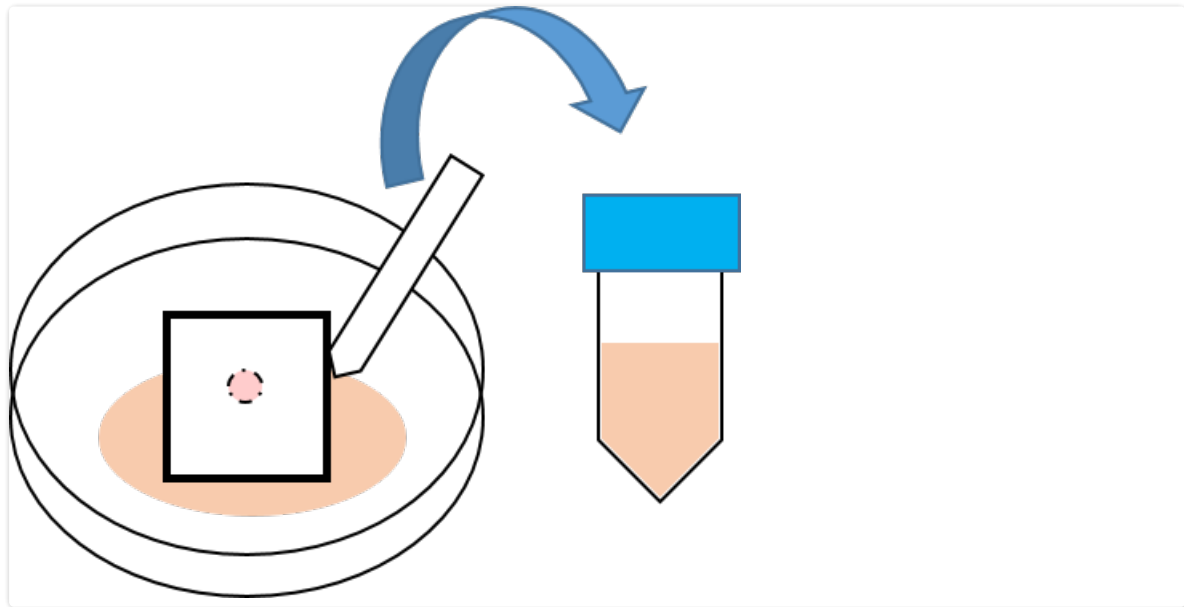

- 19 The bacterial suspension was then serially 10-fold diluted with PBS. One ml of each diluted sample was quickly mixed with 9 ml of melted Pearl Core plate count agar (PC; Eiken Chemical Co. Ltd., Tokyo, Japan) kept at 50°C in a water bath, and poured onto a petri dish and cooled down to room temperature.
- 20 The bacteria suspension in Pearl Core plate count agar was cultured 24 hours at 37°C and the number of visible colony in the dish was counted.