



AUG 25, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.n92ldpoe7l5b/v1](https://dx.doi.org/10.17504/protocols.io.n92ldpoe7l5b/v1)

**Protocol Citation:** Javier Espadas, Aurelien Roux 2023. Generation of membrane tubules by lipid-covered silica beads rolling. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.n92ldpoe7l5b/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jun 02, 2023

## 🌐 Generation of membrane tubules by lipid-covered silica beads rolling

Javier Espadas<sup>1</sup>, Aurelien Roux<sup>1</sup>

<sup>1</sup>Department of Biochemistry, University of Geneva, Geneva, 1211, Switzerland

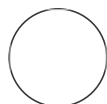

xinbo.wang

### ABSTRACT

This protocol explains the high throughput methodology to generate membrane tubules from lipid-covered silica beads.

### ATTACHMENTS

[733-1846.docx](#)

### MATERIALS

#### Materials:

##### ■ Lipids:

⊗ 1,2-dioleoyl-sn-glycero-3-phosphocholine (18:1, 18:1 PC) Avanti Polar Lipids, Inc. Catalog #850375P

⊗ 12-dioleoyl-sn-glycero-3-phospho-L-serine (sodium salt) Avanti Polar Lipids, Inc. Catalog #840035

⊗ 12-Dioleoyl-sn-glycero-3-phosphoethanolamine labeled with Atto 647N Merck MilliporeSigma (Sigma-Aldrich) Catalog #42247

■ Glass vials (2700 Supelco, Sigma-Aldrich).

■ ⊗ Silica Beads Microspheres-Nanospheres Catalog #140256-10

■ Parafilm.



■ Petri dish.

■ ⊗ Chloroform Merck MilliporeSigma (Sigma-Aldrich) Catalog #650498

**Last Modified:** Aug 25, 2023

**PROTOCOL integer ID:** 82811

**Keywords:** membrane tubules , lipid-covered silica beads

-  sticky-Slide VI  
0.4 Ibidi Catalog #80608
-  Bovine serum albumin 2 mg/ml Thermo Fisher Scientific Catalog #23209



#### Solutions:

- Lipid films hydration buffer A: 25 mM HEPES 7.4.
- Working buffer:

A
20mM HEPES 7.4
150mM NaCl
2.5mM MgCl <sub>2</sub>
5% Glycerol
2mM DTT

## Protocol


3h 15m

- 1 Mix of DOPC, DOPS and Atto 647N DOPE at 59.9:40:0.1 mol% respectively in a final volume of  200  $\mu$ L with chloroform and  0.5 g/L lipid final concentration in a glass vial.



- 2 Dry the lipid mixture in the glass vials for  02:00:00 in a vacuum chamber forming the dried lipid films on the bottom of the glass vials.

2h


- 3 Add  200  $\mu$ L of the lipid films hydration buffer A to the glass vial containing the dried lipid films.




- 4 Vortex the glass vials until visually seeing complete resuspension of the dried lipid films in the solution (seen by an increase in the turbidity of the lipid solution) forming the multilamellar vesicles (MLVs).

5 Mix  10  $\mu\text{L}$  of MLVs with  2  $\mu\text{L}$  of silica beads in an Eppendorf.






6 Deposit 6 drops of  2  $\mu\text{L}$  each containing the mixture of MLVs and silica beads on a parafilm slide placed in the bottom of a petri dish.

7 Dry the drops for  01:00:00 in the vacuum chamber until the liquid is completely dried.

1h



8 Stick a microfluidic device on a 1.5 borosilicate coverslip.

9 Passivate the microfluidic channels by adding  200  $\mu\text{L}$  solution of  2 g/L BSA for  00:15:00 .

15m

10 Clean the BSA solution by passing  200  $\mu\text{L}$  working buffer solution in each channel 5 times.



11 Add  200  $\mu\text{L}$  of working buffer to each channel with a final concentration of GFP-LRRK2 of  500 nanomolar (nM) .



12 Pick one of the dried drops and add it to the inlet of the microfluidic device.

Gently tilt the chamber towards you 60 degrees with the inlet in the upper part and the outlet in

- 13** the lower part, and wait until visually seeing the lipid-covered silica beads moving from the inlet to the outlet openings of the microfluidic device.
- 14** Wait until reaching the steady state of protein coverage on the lipid tubules.