

FEB 06, 2023

🌐 Inoue transformation buffer

Andreas Sagen¹

¹University of Oslo



Andreas Sagen

University of Oslo, The National Institute of Occupational H...

ABSTRACT

The Inoue method of transformation is a chemical transformation method of prokaryotes (bacteria), which add significant improvements to methods such as the classic Calcium chloride transformation method.

MATERIALS

LAF

Scale

Filter unit

OPEN ACCESS

Protocol Citation: Andreas Sagen 2023. Inoue transformation buffer.
protocols.io
<https://protocols.io/view/inoue-transformation-buffer-cnrbve2n>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working



Created: Feb 05, 2023



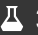

Last Modified: Feb 06, 2023

PROTOCOL integer ID:
76419


Keywords: Transformation, bacteria, chemical, Inoue

10x Inoue transformation buffer

1 Add  80 mL distilled water to a  100 mL tube

- 2 Measure and add  10.89 g Manganese chloride,  2.2 g Calcium chloride,  3.0 g PIPES and  18.65 g Potassium chloride


Materials:


 Manganese(II) chloride tetrahydrate **Sigma-aldrich Catalog #M3634**

 Calcium chloride dihydrate **Sigma-aldrich Catalog #C3881**

 PIPES **Sigma-aldrich Catalog #P1851**

 Potassium chloride **Sigma-aldrich Catalog #P3911**

- 3 Add distilled water to  100 mL

- 4 Filter sterilize with a 0.2 μm pore-size filter and store frozen at  -20 °C