



JAN 10, 2024

OPEN ACCESS



**Protocol Citation:** Abhinaya Srikanth, Brian Merritt, Peter Thielen 2024. Integrated Indexing One-Step RT-PCR with NEB LunaScript: Influenza A Virus Sequencing with Native Oxford Nanopore Software Compatibility. [protocols.io](https://protocols.io/view/integrated-indexing-one-step-rt-pcr-with-neb-lunas-c63ezgje) <https://protocols.io/view/integrated-indexing-one-step-rt-pcr-with-neb-lunas-c63ezgje>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**Created:** Jan 05, 2024

**Last Modified:** Jan 10, 2024

**PROTOCOL integer ID:** 92998

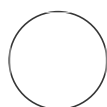
## Integrated Indexing One-Step RT-PCR with NEB LunaScript: Influenza A Virus Sequencing with Native Oxford Nanopore Software Compatibility

Abhinaya  
Srikanth<sup>1</sup>,

Brian  
Merritt<sup>1</sup>,

Peter  
Thielen<sup>1</sup>

<sup>1</sup>Johns Hopkins University Applied Physics Laboratory



Peter Thielen

### ABSTRACT

This integrated indexing RT-PCR protocol has been optimized for whole genome sequencing of influenza A virus (IAV) using Oxford Nanopore's Ligation Sequencing Kit V14. The updated protocol is an optimization of an [early approach](#) used across several laboratories for research purposes, and uses New England Biolabs reagent sets including LunaScript® Multiplex One-Step RT-PCR (#E1555), NEBNext dA-Tailing Module (E6053), and the NEBNext Quick Ligation Module (E6056). Previously demonstrated IAV whole genome sequencing approaches for Oxford Nanopore used Superscript III/Platinum Taq RT-PCR reagents input into the Native Indexing reagent set, which results in a complex, multi-step process that requires more than a full working day to carry out.

The improved integrated indexing protocol can be completed within a working day, and incorporation of the molecular index during first-strand cDNA synthesis in the RT-PCR step reduces the likelihood of sample cross-contamination during handling. Recent improvements to the approach additionally include Oxford Nanopore defined constant regions that are required for full compatibility with default informatics packages. This update now enables native demultiplexing within either Guppy or Dorado, effectively serving as a "molecular patch" where previous implementations required software modifications that were often broken by silent updates to MinKNOW.

As designed, this protocol has potential for extension to other viruses of interest by modifying primer sets. A simplified excel workbook is provided for design of new primer sets, for which the first worksheet requires input of gene-specific primer sets, the second worksheet provides fully constructed primers, and the third can be exported directly for submission to a primer synthesis vendor such as IDT DNA.

### ATTACHMENTS

[Process Overview.png](#)

### MATERIALS

New England Biolabs Components:

A	B	C	D	E	F
<b>NE B #</b>	<b>Component Name</b>	<b>Component #</b>	<b>Stored at (°C)</b>	<b>Amount</b>	<b>Concentration</b>
E6056	NEBNext Quick Ligation Reaction Buffer	E6058AVIAL	-20	1 x 0.2 ml	5 X
E6056	Quick T4 DNA Ligase	E6057AVIAL	-20	1 x 0.1 ml	Not Applicable
E6053	NEBNext dA-Tailing Reaction Buffer	E6055AVIAL	-20	1 x 0.1 ml	10 X
E6053	Klenow Fragment (3' → 5' exo-)	E6054AVIAL	-20	1 x 0.06 ml	Not Applicable
E1555	LunaScript <sup>®</sup> Multiplex One-Step RT-PCR Enzyme Mix	M1556SVIAL	-20	1 x 0.05 ml	25 X
E1555	LunaScript <sup>®</sup> Multiplex One-Step RT-PCR Reaction Mix	M1557SVIAL	-20	1 x 0.25 ml	5X
E1555	Nuclease-free Water	B1502AVIAL	-20	1 x 1.5 ml	Not Applicable

Oxford Nanopore Components:

A	B	C	D	E	F
<b>ONT #</b>	<b>Component Name</b>	<b>Component #</b>	<b>Stored at (°C)</b>	<b>Amount</b>	<b>Concentration</b>
SQK-LSK114	Ligation Sequencing Kit V14		-20	1 x 0.2 ml	6 Runs
	R10.4.1 Sequencing Flowcell		4		

Other Equipment:

A
---

A
<b><i>Component Name</i></b>
Qubit
Tapestation
Ethanol (molecular biology grade)
Custom Indices

#### BEFORE START INSTRUCTIONS

Note: High quality isolated Influenza A RNA must be used for this protocol as it is essential for successful complete genome amplification. There are many options for RNA purification, including solid phase reversible immobilization methods (e.g. spin columns, paramagnetic beads), and Trizol derivatives. Freshly isolated RNA with minimal or no freeze thaw cycles typically ensures the highest success rates for complete genome amplification.

## Protocol for Integrated Indexing of Influenza A Virus

1



Prepare the One-Step Reverse Transcriptase Multi-Segment PCR for Influenza A Virus as described below:

A	B
Reagent	Volume (uL)
Purified Viral RNA	5
5x Buffer (Reaction mix)	10
Nuclease Free Water (NFW)	31
Indexed Primer (MBTuni12-NB##)	1
Unindexed Primer (MBTuni13)	1
25x Enzyme	2
Total Volume	50

\*Resulting amplicons should range from 700 bp to 2.1kb

\*\* For no template control replace purified viral RNA with NFW

- 2 Flick the tube or pipet up and down to mix at minimum 10 times followed by a short spin



- 3 Incubate reaction in thermo-cycler following these steps:



A	B	C	D
Step	Temp (C)	Duration (min)	Cycles
RT Priming	55	2	1
Reverse Transcription	42	20	1
RT kill/ PCR Hot Start	98	2	1
PCR Round One	95	.5	5
	44	.5	
	68	3.5	
PCR Round Two	95	.5	26
	57	.5	
	68	3.5	
Final Extension	68	10	1
Hold	10	∞	∞

### 3. Library Pooling and Purification


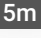
- 4 Quantify unpurified PCR product using Qubit High Sensitivity DNA reagents

- 5 Pool indexed PCR products together by equal mass into a single tube for purification

6 Add .6 volumes of AMPure XP to the pooled iiRT-PCR product and mix by finger flicking

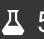
**Note**

Example: 60uL AMPure per 100uL pooled library

7 Incubate for  00:05:00 at room temperature or if available incubate on hula mixer at low speeds  IKA plate shaker at 500 RPM

8 Move to magnetic stand to separate beads, allowing them to sit until solution clears

9 Remove and discard supernatant

10 Add  500  $\mu$ L of freshly made 80% ethanol



11 Remove tube from magnet and finger flick

12 Repeat steps 8-11 for a total of two washes

13 Allow to dry for  00:05:00 or until bead pellet is no longer shiny 

#### Note

Do not over dry! Pellet will be matte, but will not have cracked.

- 14 Add  31  $\mu\text{L}$  nuclease-free water and resuspend with gentle pipetting. Incubate for  00:10:00 10m

#### Note

This 10 minute elution duration is longer than normal Ampure XP, and is intended to maximize the release of larger PCR products.

- 15 Move to magnetic stand and allow to clear for  00:05:00 5m

- 16 Transfer purified DNA to new tube




- 17 Quantify  1  $\mu\text{L}$  of final pooled iiRT-PCR product with high sensitivity qubit reagents




- 18 (Optional) Following purification of iiRT-PCR products, samples can be analyzed on an agarose gel, Tapestation, or similar platform.

## DNA A-Tailing with NEBNext dA-Tailing Module

42m 30s

- 19 Thaw all NEB consumables on ice and once thawed allow mixture to adjust to room temperature. The  30s  
 pipette the buffer up and down to. Break up the precipitate and then a  00:00:30 vortex.

- 20 Place 100-200 fmol of barcoded, pooled and purified DNA into a 1.5 Eppendorf tube and adjust volume to


 30 µL with nuclease free water

21 Mix in a 0.2 ml PCR plate and then flick and spin down:

A	B
Reagents	Volume (µL)
DNA	30
NEBNext dA-Tailing. Reaction Buffer (10x)	5
Klenow Fragment (3' --> 5' exo-)	3
Nuclease Free Water (NFW)	12
Total Volume	50

22 Incubate plate in thermal cycler at  37 °C for  00:30:00 and then  65 °C for  00:05:00 

23 Transfer DNA into 1.5 ml tubes

24 Add  60 µL of AMPure beads to the tube and flick to mix

25 Incubate on hula mixer for  00:05:00 

26 Spin down and place on magnet. Once solution is clear pipette off supernatant



- 27 Add  200  $\mu\text{L}$  of 70% ethanol and remove ethanol without disturbing beads
- 28 Repeat step 27
- 29 Take tube off of magnet and spin down, then place back on magnet and remove any residual ethanol
- 30 Remove tubes from magnet and resuspend pellets in  61  $\mu\text{L}$  of nuclease free water
- 31 Incubate for  00:02:00 at room temperature 2m
- 32 Place tubes on magnet and once solution is clear remove supernatant and place in clean 1.5 ml tubes
- 33 Quantify  1  $\mu\text{L}$  using qubit

## Adapter ligation and clean-up

37m

- 34 Thaw all NEB reagents at room temperature, spin down for 5 seconds, and mix by pipetting. DON'T VORTEX T4 DNA Ligase. Spin down the Ligation Adapter and Quick T4 DNA Ligase, pipette to mix then place on ice. Thaw the Elution Buffer and Short Fragment Buffer at room temp, mix by vortexing, spin down, and then place on ice.
- 35 Mix by pipetting and spin down the following reaction:



A	B
Reagents	Volume (uL)
Pooled Barcoded Sample	30
Ligation Adapter (LA)	5
NEBNext Quick Ligation Reaction Buffer (5X)	10
Quick T4 DNA Ligase	5
Total	50


36 Incubate the reaction for  00:20:00 at room temperature


20m



37 Transfer DNA into 1.5 ml tubes




38 Add  20  $\mu$ L of AMPure beads to each tube and flick to mix

39 Incubate on hula mixer for  00:05:00

5m







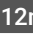
40 Spin down and place on magnet. Once solution is clear pipette off supernatant

41 Add  125  $\mu\text{L}$  of Short Fragment Buffer, flick to resuspend, then place back on tube and allow beads to pellet, then remove short fragment buffer without disturbing beads

42 Repeat step 41

43 Take tube off of magnet and spin down, then place back on magnet and remove any residual Buffer

44 Remove tubes from magnet and resuspend pellets in  15  $\mu\text{L}$  of NFW

45 Incubate for  00:10:00 at  37 °C and every  00:02:00 agitate by lightly flicking sample  12m




46 Place tubes on magnet and once solution is clear remove supernatant and place in a clean 1.5 ml tube

#### Note

This is the final clean library, and should be stored on ice or at 4C until sequencing

## Final Library Quantification

47 Quantify  1  $\mu\text{L}$  of the final clean library using a Qubit DNA fluorometer

48 Optional: Run final library on a Tapestation to evaluate average peak size

- 49** Calculate the final volume required for 200fmol to be sequenced on a MinION flowcell

**Note**

Tip: Use the NEBioCalculator <https://nebiocalculator.neb.com/#!/dsdnaamt>. Assuming 1kb product mean, a final loading concentration of 200fmol requires 125ng of library into the flowcell.

- 50** Proceed to load the Oxford Nanopore flowcell using the most up to date protocols from the vendor.