




SEP 30, 2023

Cryo-EM Grid preparation

Annan SI Cook^{1,2}

¹University of California, Berkeley; ²Aligning Science Across Parkinson's

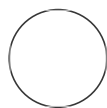

Annan SI Cook

ABSTRACT

This protocol details Cryo-EM grid preparation.

ATTACHMENTS

[852-2198.pdf](#)

MATERIALS

Materials

- PI3KC3-C1~RAB1A sample with n-Octyl- β -D-Glucopyranoside (OG) at a final concentration of 0.05%
- QUANTIFOIL R 2/1 mesh Cu 300 holey carbon grids
- Vitrobot cryo-plunger (Thermo Fisher Scientific)
- Cryo-EM microscope (300 kV Titan Krios with X-FEG and energy filter)
- K3 Summit direct electron detector (Gatan)
- SerialEM software
- Liquid nitrogen
- Humidity control system

OPEN  ACCESS



Protocol Citation: Annan SI Cook 2023. Cryo-EM Grid preparation. **protocols.io** <https://protocols.io/view/cryo-em-grid-preparation-c2pbydin>


License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Sep 22, 2023


Last Modified: Sep 30, 2023

PROTOCOL integer ID:
88515

- 1 Freshly glow-discharge the QUANTIFOIL R 2/1 mesh Cu 300 holey carbon grids using a glow-discharge system (PELCO easiGlow) at 25 mA current for  00:01:00 .

1m

Sample Application

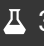- 2 Mix the PI3KC3-C1~RAB1A sample with 0.05% n-Octyl- β -D-Glucopyranoside (OG) to achieve the desired final concentration.
- 3 Apply a small volume (typically  3-4 μ L) of the sample onto the glow-discharged side of the grids.

Vitrification

3s

- 4 Load the grids with the sample into the Vitrobot cryo-plunger.

- 5 Set the Vitrobot to the following conditions:

- 5.1 Maintain 100% humidity inside the chamber.

- 5.2 Set the temperature to  4 °C .

- 5.3 Wait for  00:00:03 .

3s

- 5.4 Blot using a blotting force of -15 (may need calibration depending on your particular instrument) .

Data Acquisition Setup

- 6 Prepare the Cryo-EM microscope (300 kV Titan Krios) for data acquisition.
- 7 Used X-FEG and an energy filter set to 20 eV.

Data Collection Parameters

- 8 Use SerialEM software for automated data collection.
- 9 Set the following parameters for data collection:
- Magnification: 81,000x.
 - Super-resolution pixel size: 0.525 Å.
 - Defocus range: -0.8 to -2.2 micrometers.
 - Collect 50-frame image stacks.

Data Collection

- 10 Initiate the data collection process using the specified parameters.
- 11 Monitor the data acquisition progress and ensure that the images are being recorded.

Cumulative Dose Control

- 12 Each movie has a cumulative electron dose of $50 \text{ e}/\text{\AA}^2$.