



VERSION 2
JAN 11, 2024

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.kxygx3y6kg8j/v2

Protocol Citation: Verena Dederer, chatterjeedeep, Sebastian Mathea, Stefan Knapp 2024. LRRK2 thermal shift assay. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.kxygx3y6kg8j/v2>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Sep 07, 2023


🌐 LRRK2 thermal shift assay V.2

Verena

Dederer^{1,2}, chatterjeedeep^{1,2}, Sebastian Mathea^{1,2},
Stefan
Knapp^{1,2}

¹Institute of Pharmaceutical Chemistry, Goethe University Frankfurt, Max-von-Laue-Straße 9, Frankfurt 60438, Germany;

²Structural Genomics Consortium, Buchman Institute for Molecular Life Science (BMLS), Max-von-Laue-Straße 15, Frankfurt 60438, Germany

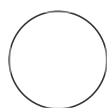

Verena Dederer

ABSTRACT

Thermal shift assay or differential scanning fluorimetry analyzes the effect of small molecules on the thermostability of a protein by gradual heat denaturation and monitoring absorption of the fluorescent dye SYPR Orange at 488 nm.

Last Modified: Jan 11, 2024

PROTOCOL integer ID:
87480

Funders

Acknowledgement:

Aligning Science Across

Parkinson's grant

Grant ID: ASAP-000519

Fluorescent-based thermal shift assay

- 1 Prepare 4 μ M master mix of protein in buffer (20 mM Hepes pH 7.4, 150 mM NaCl, 5% glycerol) and add 1:1000 dilution of SYPR Orange.
- 2 Aliquot 20 μ L of the master mix into a white 96 well plate.
- 3 Add DMSO or small molecule binder with a final concentration of 10 μ M.
- 4 Seal plate, mix well and centrifuge 30 sec at 500xg.
- 5 Place plate into MX3005P real-time PCR instrument.
- 6 Measure fluorescence with excitation and emission filters set to 465 and 590 nm while gradually increase temperature 3K/min during 71 cycles from 25 to 95°C.