-



NOV 09, 2022



WORKS FOR ME 1

ONA extraction from fecal samples


DO

## dx.doi.org/10.17504/protocols.io.5jyl8jnorg2w/v1

Anique Ahmad<sup>1</sup>, Arya Gautam<sup>1</sup>, Tsedenia Denekew<sup>1</sup>, Aashish Jha<sup>1</sup>

<sup>1</sup>New York University, Abu Dhabi

Aashish Jha: jhaar@nyu.edu

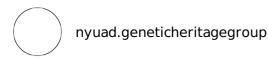

COMMENTS 0

**ABSTRACT** 

DNA Extraction from fecal samples

DOI

# dx.doi.org/10.17504/protocols.io.5jyl8jnorg2w/v1

PROTOCOL CITATION

Anique Ahmad, Arya Gautam, Tsedenia Denekew, Aashish Jha 2022. DNA extraction from fecal samples. **protocols.io** 

https://dx.doi.org/10.17504/protocols.io.5jyl8jnorg2w/v1

LICENSE

This is an open access protocol distributed under the terms of the <u>Creative</u> <u>Commons Attribution License</u>, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

protocols.io

1

**Citation**: Anique Ahmad, Arya Gautam, Tsedenia Denekew, Aashish Jha DNAÃÂ extraction from fecal samples <a href="https://dx.doi.org/10.17504/protocols.io.5jyl8jnorg2w/v1">https://dx.doi.org/10.17504/protocols.io.5jyl8jnorg2w/v1</a>

**CREATED** 

Sep 16, 2022

LAST MODIFIED

Nov 09, 2022

PROTOCOL INTEGER ID

70128

### **GUIDELINES**

Sample Type	Maximum Input
Feces	200 mg
Soil	250 mg
Liquid Samples <sup>1</sup> and Swab Collections <sup>2</sup>	250 µl
Cells (isotonic buffer, e.g. PBS)	50-100 mg (wet weight) (109 bacterial and 108 yeast cells)
Samples in DNA/RNA Shield™,3	≤ 1 ml

Table 1. Sample type and maximum input for DNA extraction using this protocol

#### **MATERIALS TEXT**

ZymoBIOMICS DNA Extraction Kit
Fecal samples
Homogenizing pestle
Liquid nitrogen or dry ice
Ice box with ice
1.5 ml microcentrifuge tubes
Vortex
Omni Bead Ruptor Elite bead beater
Centrifuge

#### **BEFORE STARTING**

Label the following:

ZR BashingBead<sup>TM</sup> Lysis Tubes (0.1 & 0.5 mm)

 ${\sf ZymoSpin^{TM}\,III\text{-}F\,Filter\,in\,Collection\,Tube}$ 

**Collection Tube** 

ZymoSpin<sup>TM</sup> II-CR Column in Collection Tube

ZymoSpin<sup>TM</sup> III-HRC Filter in Collection Tube

2 Sets of 1.5 ml microcentrifuge tubes (not provided with kit)

■ Spin the labelled ZR BashingBead<sup>TM</sup> Lysis Tubes (0.1 & 0.5 mm) for 10 seconds in a mini-centrifuge to



2

ensure that the beads have settled at the bottom

- Include 1 control per batch and assign its position randomly.
- Transfer fecal sample to ZR BashingBead™ Lysis Tube (0.1 & 0.5 mm) and add ZymoBIOMICS™ Lysis Solution:

10s

- Flash frozen fecal samples will be solid and an aliquot ( ♣ 15 mg ) of feces is expected in a 2mL container (Table 1). Add ♣ 250 μL of ZymoBIOMICS™ Lysis Solution to the fecal sample and homogenize evenly using the homogenizing pestle while submerged in liquid nitrogen (or on dry ice). Once the sample is homogenized, add another ♣ 500 μL of ZymoBIOMICS™ Lysis Solution such that the final volume is ♣ 750 μL . Mix well by vortexing.
  - Spin the sample container for 00:00:10 in a minicentrifuge.
  - Transfer the entire  $\bot$  750  $\mu$ L of the sample mix into the ZR BashingBead  $\bot$  Lysis Tube (0.1 & 0.5 mm) and proceed to Step 2.
- 1.2 Fecal samples collected using Zymo DNA/RNA Shield™ will exist as a solution (Table 1). Mix the sample by shaking well or vortexing. Add up to Δ 250 μL of such sample to a ZR BashingBead™ Lysis Tube (0.1 & 0.5 mm). Adjust final volume to Δ 1 mL with ZymoBIOMICS™ Lysis Solution. Cap the tube tightly and proceed to Step 2.


Secure in a Omni Bead Ruptor Elite bead beater fitted with a 2 ml tube holder assembly and process at max speed (30 m/s) for 00:05:00 . Rest for 00:05:00 .

3 Centrifuge the ZR BashingBead<sup>™</sup> Lysis Tubes (0.1 & 0.5 mm) in a microcentrifuge at  $\geq$  10000 x g, 00:01:00.

Transfer up to  $\Delta$  600  $\mu$ L supernatant to the Zymo- $Spin^{7}$  III-F Filter in the labelled Collection Tube and centrifuge at  $8.000 \times g$ , 00:01:00. Discard the Zymo- $Spin^{7}$  III-F Filter.

10m

5	Add $\bot$ 600 $\mu$ L of ZymoBIOMICS <sup>™</sup> DNA Binding Buffer to the filtrate in <i>the labelled Collection Tube</i> from Step 4. Mix well.	
5.1	Repeat so that the final volume of ZymoBIOMICS™ DNA Binding Buffer added is	
6	Transfer $\  \  \  \  \  \  \  \  \  \  \  \  \ $	11
6.1	Discard the flow through from the Collection Tube and repeat step 6.	
7	Add	11
8	Add   ZymoBIOMICS™ DNA Wash Buffer 2 to the <i>Zymo-Spin™ IICR Column</i> in a Collection Tube and centrifuge at 10.000 x g, 00:01:00 . Discard the flow-through.	11
9	Repeat with ∠ 200 µL ZymoBIOMICS™ DNA Wash Buffer 2 to the <i>Zymo-Spin™ IICR Column</i> in a Collection Tube and centrifuge at 10.000 x g, 00:01:00	11
10	Transfer the Zymo-Spin™ IICR Column to a clean 1.5 ml microcentrifuge tube and add (50 µl minimum) ZymoBIOMICS™ DNase/RNase Free Water directly to the column matrix and incubate for (50 00:05:00). Centrifuge at (10.000 x g, 00:01:00) to elute the DNA	61
11	Place a <i>Zymo-Spin™ III-HRC Filter</i> in <i>a new Collection Tube</i> and add ∠ 600 µL ZymoBIOMICS™	31

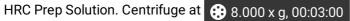

3m

- 12 Transfer the eluted DNA (Step 10) to a prepared Zymo-Spin™ III-HRC Filter in a clean 1.5 ml microcentrifuge tube and centrifuge at exactly ( 16.000 x g, 00:03:00
- 13 The filtered DNA is now suitable for PCR and other downstream applications. Eluted DNA should be frozen ( $-30 \text{ to } -15^{\circ}\text{C} \text{ or } -90 \text{ to } -65^{\circ}\text{C}$ )