




MAR 22, 2023

Sequential Microbial Biomass and Nitrogen Extraction

In 1 collection

maggie.bowman¹

¹Environmental Molecular Sciences Laboratory

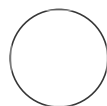

nicholas.sconzo

ABSTRACT

The method is modified from the Hofmockel lab at PNNL, modified from Suding Lab protocol, modified from S. E. Hobbie, 5 May 1998.

MATERIALS

- [M] 0.5 Molarity (M) K₂SO₄ (87.13 g or 8.71 g K₂SO₄ per 1 L extractant made in ultrapure water)
Always test background concentrations of new lots of K₂SO₄!
- Buchner Funnel
- Ethanol cleaned forceps
- Whatman No.42 filter paper, pre-leached 3x with K₂SO₄ ([M] 0.5 Molarity (M) or [M] 0.05 Molarity (M) depending on sample extractant) and 1x with ultrapure water; (Whatman Part#1442-047)
See Rinsing Filters section below
- Vacuum Pump
- 50 mL conical tubes with V-bottom (new sterile or acid washed)-label with Sample_ID, date, initials; 1 per sample + blanks
- 2-gallon glass jar with bottle-top repeat dispenser, calibrate to 24 mL
- Filter supports
- Acid washed funnels
- 20 mL glass scintillation vials (acid washed)- 2 per sample + blanks- label with Cat#, date, and either "background" or "fumigated"
- Ethanol-free chloroform (for analysis of carbon) (VWR # BJ049-1L)
- Viton/Butyl Gloves
- Phosphoric acid for preservation
- Dropper for acid
- Manual 2 mL pipet pump
- Data sheets

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.q26g7yr6qgwz/v1

Protocol Citation: maggie.bowman 2023. Sequential Microbial Biomass and Nitrogen Extraction. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.q26g7yr6qgwz/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


Protocol status: Working

Created: Sep 14, 2022

Last Modified: Mar 22, 2023

PROTOCOL integer ID:
70055

Subsample soil

- 1 Subsample (~  8 g) place in V- bottomed conical tube labeled with Cat# and weight. Label caps with Cat#
- 2 Gravimetric water content will also be essential for final calculations on a dry weight basis.

Data to Collect






- 3 Data to Collect
 - 3.1 Weight of soil and Cat# for each sample on the data sheet
 - 3.2 Time onto and off the shaker 1st extraction (by batch)
 - 3.3 Time of chloroform addition (by batch)
 - 3.4 Time of chloroform evacuation (by batch)
 - 3.5 Time onto and off the shaker 2nd extraction (by batch)


Rinsing Filters

- 4 Gather the vacuum pump, Erlenmeyer filter flask, number of filters needed, K₂SO₄, Nano pure water, Buchner funnel, and tweezers.
- 5 Set up a tray with bench paper. Wipe the plastic side of the bench paper with ethanol.
- 6 Set up vacuum, Erlenmeyer flask, and Buchner funnel.
- 7 Place a stack of 5-8 Filters in the funnel.
- 8 to cover the filters. Turn on the vacuum pump.
- 9 Once the K₂SO₄ has been pulled through, repeat with K₂SO₄ two more times.
- 10 Wait a few seconds to ensure all K₂SO₄ is pulled through (you will hear the vacuum noise change) Then add ultrapure water until filters are completely covered.
- 11 Wait a few seconds to vacuum off as much liquid as possible then turn off the vacuum and use tweezers to remove the filters.

- 12 Spread the filters out on the tray with the clean bench paper. Place the tray in the hood overnight to dry.

Extraction

- 13 Work in batches (when possible randomized experimental treatments among batches, but do not mix isotopically labelled and natural abundance samples in order to prevent cross contamination)
- 14 Include multiple (minimum 3 per day) empty  50 mL conical tubes. Carry these through the entire process as blanks.
- 15 Add  24 mL K₂SO₄ to non-fumigated subsample using repeat dispenser. Cap well.
- 16 Place on shaker for  200 rpm, 02:00:00 or low setting at room temp ( 20 °C)
- 17 Balance samples and centrifuge @  4000 rpm, 00:20:00 . 20m
- 18 Place funnels in support racks and insert a pre-leached filter in each funnel.
- 19 Remove sample from centrifuge and place scintillation vial with matching Cat# below funnel.

- 20 Pour supernatant into funnel, take care to check that filtrate is dripping directly into the vial.
- 21 Leave at least 1 cm head space in vials to prevent bursting (collect excess in waste container)
- 22 Using ethanol clean forceps transfer filter paper back into conical tube with soil pellet.
- 23 In the fume hood- add 2 drops of phosphoric acid to the extract (for preservation for TOC analysis).
- 24 Store extracts at  -20 °C .
- 25 Using ethanol clean forceps transfer filter paper back into conical tube with soil pellet.
- 26 Rinse and acid wash funnels for next extraction.




FUMIGATION


3d

27 FUMIGATIONS MUST BE DONE IN THE FUME HOOD

CHLOROFORM IS A KNOWN CARCINOGEN

Wear Viton/Butyl gloves with disposable nitrile gloves inside during steps 28-31. Always wear a lab coat and safety glasses with side guards when working with chloroform.

- 28 Add  2 mL chloroform to soil pellet, plus filter.
- 29 Add additional labeling to tube and cap if needed, chloroform will cause ink to run if it comes in contact with ink
- 30 Let sit in the dark in the hood for  24:00:00 : cover the tube racks with a black garbage bag (darkness prevents the chloroform from breaking down). 1d
- 31 Remove garbage bags and caps and allow to vent  48:00:00 . Be sure to keep track of which cap goes to each tube. 2d
- 31.1 Note: Be sure to place notice on the fume hood that chloroform containing samples are inside
- 32 Extract the sample as above: as per steps #13-24.
- 33 As an added precaution against any residual chloroform, perform the K₂SO₄ addition, filtration, and phosphoric acid addition steps inside the hood
- 34 The filter from the second extraction can be discarded.

- 34.1** For all experiments leveraging isotopic labeling the leftover soil pellets should be archived at  20 °C. This includes 1abeled samples and the unlabeled controls.
- 34.2** For experiments not requiring the pellet: At the end of the extraction open all tubes inside the fume hood and leave to dry completely to remove any trace of chloroform before disposing of tubes containing filter and soil.