



APR 28, 2023

🌐 Molecular Diagnosis of Viral Hepatitis B Infection

Frank Twum Aboagye¹, Maame Ekua Acquah²

¹Biomedical and Public Health Research Unit, CSIR - Water Research Institute;

²West African Centre for Cell Biology of Infectious Pathogens, College of Basic and Applied Sciences, University of Ghana, Legon, Accra, Ghana



Frank Twum Aboagye

Biomedical and Public Health Research Unit, CSIR - Water Res...

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.kqdg39ozeg25/v1

Protocol Citation: Frank Twum Aboagye, Maame Ekua Acquah 2023. Molecular Diagnosis of Viral Hepatitis B Infection. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.kqdg39ozeg25/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Apr 14, 2023

Last Modified: Apr 28, 2023

PROTOCOL integer ID:
80544

Keywords: Viral Hepatitis , Hepatitis B, RT-PCR, DNA Extraction, Infectious Diseases

DISCLAIMER

This is an optimized protocol for the quantitative detection of Hepatitis B using Zymo Quick-DNA Miniprep Kit (200 prep) and Bosphore HBV Quantitative Kit. The authors do not accept any liability for the collection and handling of both samples and reagents, results from the use of the protocol and its interpretation as well as any errors or omissions that may be made. The reader should make his/her own evaluation as to the appropriateness of the procedures described.

ABSTRACT

With over 400 million HBV infections, viral hepatitis B remains a global public health concern. Diagnosis is primarily based on an immunological assay approach, which utilizes the Hepatitis B surface antigen in detection amongst other markers. This method, however, has several limitations which include the inability to detect mutation in the viral genome resulting in diagnostic escape and low antigen titers in study samples. Using an alternative approach, which is the molecular diagnostic technique such as real-time polymerase chain reaction (RT-PCR) would circumvent the limitations of immune detection. This protocol, thus, provides a step-by-step process of HBV diagnosis using RT-PCR which is a sensitive tool for diagnosis. The steps involved include sample collection and preparation, nuclei acid isolation, HBV detection and quantification using RT-PCR as well as the interpretation of results.

This protocol combines the high nuclei acid yield from isolation using Zymo Quick DNA Mini-prep Kit and Bosphore HBV Quantitative Kit for amplification of HBV DNA. The Bosphore HBV Quantitative Kit has a low detection limit of 1×10^1 IU/ml and with a turn-around time of less than 4 hours when combined with Zymo Quick DNA Mini-prep Kit. The nucleic acid isolated in this protocol can be amplified using the Bosphore HBV Quantitative kit on several thermocyclers, which makes the protocol robust, cost-efficient, and cost-effective in resource-scarce areas.

IMAGE ATTRIBUTION

<https://www.anatoliagenetworks.com/en/diseases/hepatitis-b/>

GUIDELINES

Reagent Preparation for Nucleic Acid Isolation

Nucleic Acid Isolation

Detection and Quantification of HBV Nucleic Acid

Results and Interpretation

Troubleshooting

MATERIALS

Consumables

1. Microcentrifuge tube (1.5 ml)
2. 100 – 1000 µl filtered tips
3. 20 – 200 µl filtered tips
4. 5 – 20 µl filtered tips
5. 96-well PCR plate or 8-well PCR strip
6. Permanent marker

Reagents

1. ZYMO Quick DNA Miniprep Kit
2. Bosphore HBV Quantitative Kit (Includes internal control, standards (4), PCR master mix, positive control and nuclease-free water)
3. Absolute Ethanol (Molecular Grade)
4. Proteinase K

Equipment

1. Thermocycler (Real-Time PCR)
2. Incubator
3. Vortex
4. Microcentrifuge

SAFETY WARNINGS



1. Handle all blood specimens and reagents as a potential biohazard.
2. Discard all waste materials in the appropriate receptacles.
3. β -mercaptoethanol has a pungent smell, it should be opened in a well-aerated space.
4. Ensure the use of the appropriate PPEs at all times

BEFORE START INSTRUCTIONS

Sample Collection and Preparation

Aseptically collect venous blood from the antecubital vein of the forearm and dispense it into either a serum separator tube (SST), a plain tube (red top), or an EDTA tube (plasma).

Centrifuge the specimen at 3500 rpm for 10 minutes and transfer the serum or plasma into a sterile plain tube or a cryovial.

NB:

For samples collected into SST or plain tubes, allow the specimen to stand vertically undisturbed for at least 1 hour for the specimen to clot completely before centrifuging it to yield serum



Samples collected into EDTA tubes can be centrifuged immediately to yield plasma.

Care should be taken during sample collection and processing to avoid haemolysis of the specimen

Handling of PCR Reagents

Allow the PCR Master Mix and the other components of the Bosphore HBV Quantification kit to thaw completely at 4°C before use and avoid centrifuging to thaw.

Reagent Preparation for Nucleic Acid Isolation




- 1 Add  500 µL of β-mercaptoethanol to  100 mL of Genomic Lysis Buffer



Note

Cap the bottle tightly after each use















HBV Nucleic Acid Isolation

5h 2m




- 2 Transfer  200 µL of  serum/plasma to a sterile  1.5 mL microcentrifuge tube (pre-labeled).

- 3 Add  5 µL of Internal Control (provided with the amplification kit) to each specimen and vortex for  00:00:30

30s

- 4 Add  400 µL of Genomic Lysis Buffer and  10 µL of Proteinase K to the sample.
- 5 Vortex and incubate the sample at  56 °C for 3 - 5 hours or  Overnight . 5h
- 6 Vortex the sample at  3000 rpm for  00:00:30 30s
- 7 Transfer the entire content of the 1.5 ml microcentrifuge tube into a Zymo-Spin IIC Column in a collection tube.
- 8 Centrifuge at  10000 rpm for  00:01:00 . Discard the flow-through liquid. 1m
- 9 Transfer the Zymo-Spin IIC column into a new collection tube.
- 10 Add  200 µL of DNA Pre-Wash Buffer to the spin column and centrifuge  10000 rpm for  00:01:00 . Discard the flow-through liquid. 1m
- 11 10. Add  500 µL of gDNA Wash Buffer to the spin column and centrifuge at  10000 rpm for  00:01:00 . Discard the flow-through liquid. 1m


12 11. Transfer the Zymo-Spin IIC Column into a sterile 1.5 mL microcentrifuge tube

13 12. Add  100 µL (70 - 100) of DNA Elution Buffer to the spin column and incubate at  Room temperature for  00:30:00


30m

14 13. Centrifuge at  13000 rpm for  00:00:30 to elute the DNA.

30s

15 14. Store the Viral Nucleic acid at  -20 °C pending further analysis.

RT-PCR Detection and Quantification of HBV Nuclei Acid

16 The PCR is done in a  25 µL reaction as described below:



A	B
Reagent/Component	X1 (µL)
PCR Master mix	15
Test Nucleic acid/Standards/Negative Control/Positive Control	10
Total Reaction Volume	25 µL

Allow the PCR Master Mix and the other components of the Bosphore HBV Quantification kit to thaw completely at 4°C before use and avoid centrifuging to thaw.


17 The reaction is done under the following cycling conditions:



A	B	C	D
Step	Temperature (°C)	Time (minutes)	Cycles
Initial Denaturation	95	14:30	NA
Denaturation	97	00:30	50 cycles
Annealing and Synthesis (Data Collection)	54	01:30	
Hold	32	01:00	NA


- 18** Select the appropriate fluorophore depending on the analyzer been used. For Bio-Rad CFX 1000 Series, select the **FAM dye for HBV** and **HEX for Internal Control (IC)**.

A1 FAM and HEX	A2 FAM and HEX	A3 FAM and HEX	A4 FAM and HEX	A5 FAM and HEX	A6 FAM and HEX
A7 FAM and HEX	A8 FAM and HEX	A9 FAM and HEX	A10 FAM and HEX	A11 FAM and HEX	A12 FAM and HEX
B1 FAM and HEX	B2 FAM and HEX	B3 FAM and HEX	B4 FAM and HEX	B5 FAM and HEX	B6 FAM and HEX
B7 FAM and HEX	B8 FAM and HEX	B9 FAM and HEX	B10 FAM and HEX	B11 FAM and HEX	B12 FAM and HEX
C1 FAM and HEX	C2 FAM and HEX	C3 FAM and HEX	C4 FAM and HEX	C5 FAM and HEX	C6 FAM and HEX
C7 FAM and HEX	C8 FAM and HEX	C9 FAM and HEX	C10 FAM and HEX	C11 FAM and HEX	C12 FAM and HEX
D1 FAM and HEX	D2 FAM and HEX	D3 FAM and HEX	D4 FAM and HEX	D5 FAM and HEX	D6 FAM and HEX

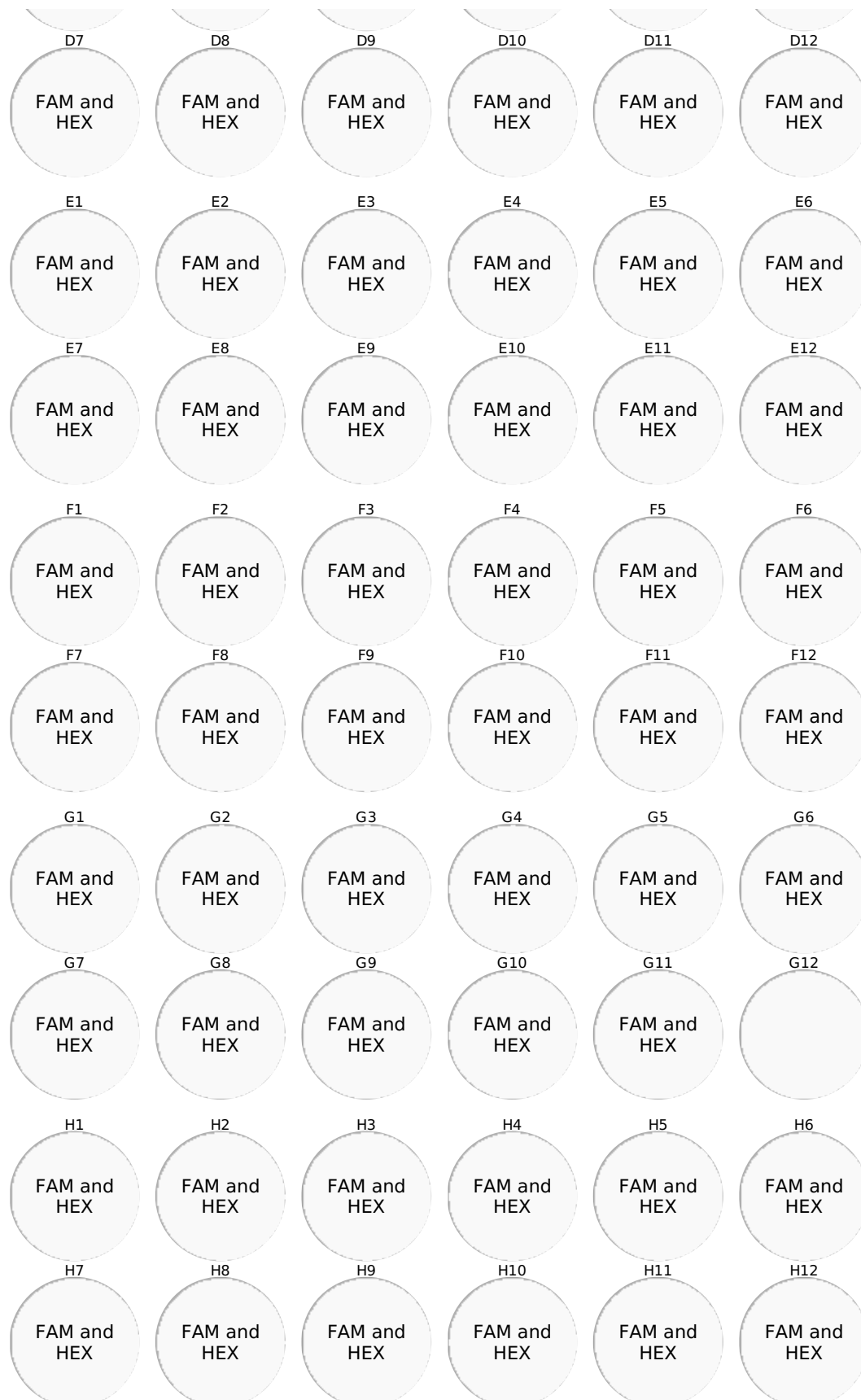

19 Identify unknown samples (test samples), standards, positive control, and negative control on the plate as described below. Assign the right quantitative values to the standards

A1 Unknown	A2 Unknown	A3 Unknown	A4 Unknown	A5 Unknown	A6 Unknown
A7 Unknown	A8 Unknown	A9 Unknown	A10 Unknown	A11 Unknown	A12 Negative Control
B1 Unknown	B2 Unknown	B3 Unknown	B4 Unknown	B5 Unknown	B6 Unknown
B7 Unknown	B8 Unknown	B9 Unknown	B10 Unknown	B11 Unknown	B12
C1 Unknown	C2 Unknown	C3 Unknown	C4 Unknown	C5 Unknown	C6 Unknown
C7 Unknown	C8 Unknown	C9 Unknown	C10 Unknown	C11 Unknown	C12 Standard 1
D1 Unknown	D2 Unknown	D3 Unknown	D4 Unknown	D5 Unknown	D6 Unknown
D7 Unknown	D8 Unknown	D9 Unknown	D10 Unknown	D11 Unknown	D12 Standard 2
E1 Unknown	E2 Unknown	E3 Unknown	E4 Unknown	E5 Unknown	E6 Unknown
E7 Unknown	E8 Unknown	E9 Unknown	E10 Unknown	E11 Unknown	E12 Standard 3

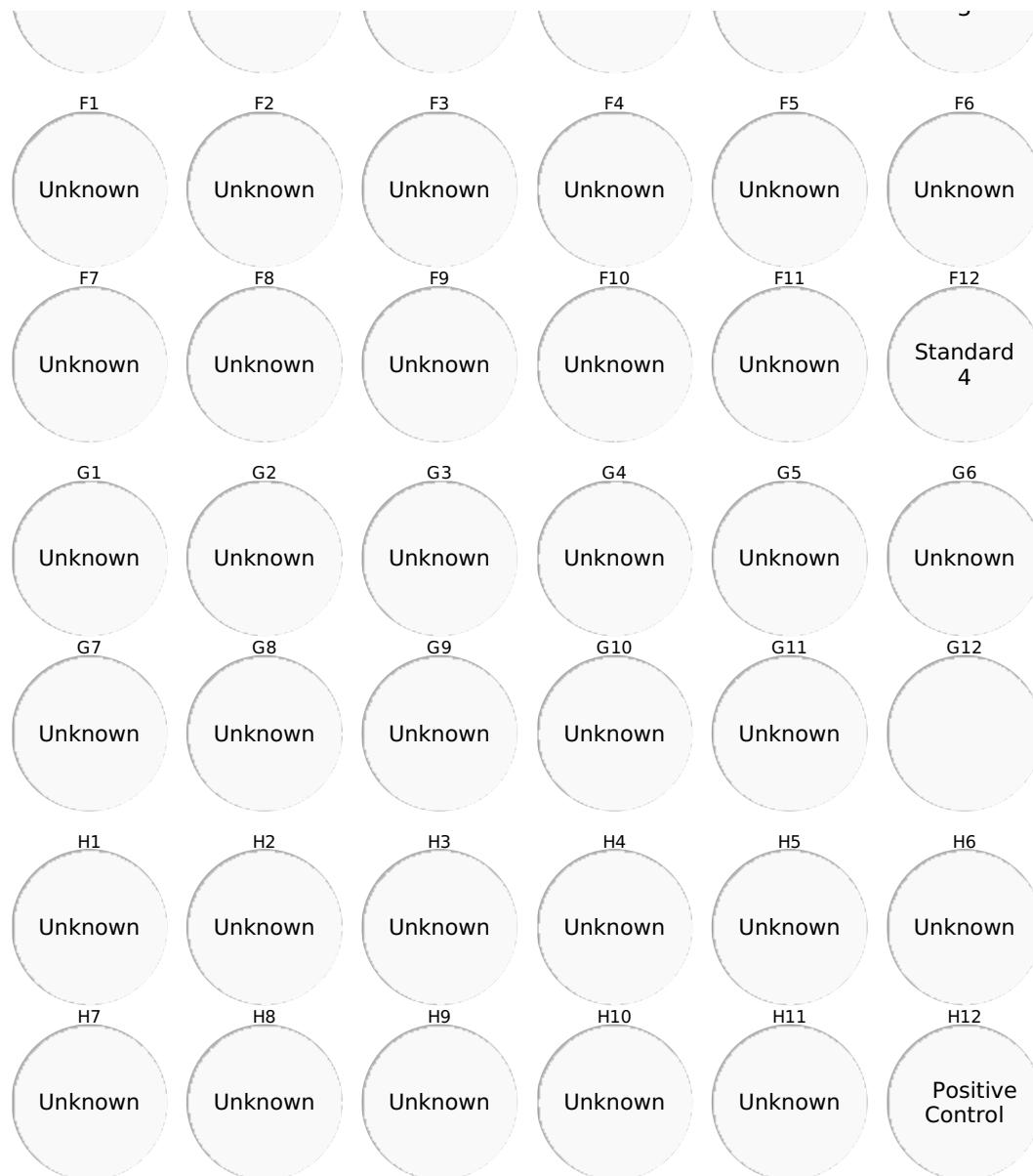

20 Initiate the protocol

Results and Interpretation

- 21** Samples that cross the threshold in the FAM channel are displayed with their starting quantities and the corresponding Ct-values. Samples that do not cross the threshold are displayed as "No Ct" or "–"

A	B	C	D
+	-	HBV Positive	The sample contains a high viral load of HBV and may suppress the amplification of the Internal Control. No need to check the HEX
+	+	HBV Positive	The sample contains HBV Nucleic acid
+	-	HBV Positive*	Caution: If the sample has a viral load <10 IU/ml and no Internal Control amplification, it indicates PCR inhibition which requires assay repetition (see troubleshooting)
-	+	HBV Negative	No HBV Nucleic acid was detected. Internal control amplified indicating good DNA isolation and PCR procedures.
-	-	Invalid**	Repeat the assay (see troubleshooting)

* and **: refer to the troubleshooting section

Troubleshooting

*If the sample has a viral load [M] <10 IU/ml and no Internal Control amplification, it requires that the assay is repeated once. Freeze-thaw[freeze the sample (temperature $\leq -10^{\circ}\text{C}$) for 20 minutes and defrost at room temperature for 10 minutes] the nucleic acid sample and dilute it with nuclease-free water in a ratio of 1:2. Use the diluted nucleic acid for repeating the assay

** Repeat the assay once taking into consideration possible pipetting errors in the first assay.

If both repeat as indicated above fail to yield the expected/better results, a new specimen should be collected from the patient and the process from nucleic acid isolation to amplification should be repeated with the new sample.