



NOV 14, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.bp2l6xm8klqe/v1

Protocol Citation: Michael X. Henderson 2023. Tissue Sectioning. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bp2l6xm8klqe/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working


Created: Nov 06, 2023

Last Modified: Nov 14, 2023

Tissue Sectioning

Michael X. Henderson¹

¹Van Andel Institute

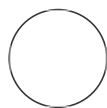

Maria Matos
Yale University

ABSTRACT

This protocol details formalin-fixed paraffin-embedded tissue sectioning on a microtome for subsequent tissue staining.

ATTACHMENTS

[n2fcbpqr.pdf](#)

Tissue Sectioning


2h 20m




- 1 Clean microtome blade with xylene (to remove paraffin) and then 100% ethanol (to dry).
- 2 Mount blade on microtome, and tighten screws to hold it in place.
- 3 Mount paraffin block in holder. Manually move the block near, but not touching the blade.
- 4 Align the block to the blade in both left-right and top-bottom directions by iterative adjustments followed by checking the alignment.
- 5 Once the block is directly apposed to the blade and aligned, switch cutting depth to **→|+ 6 μ m**.
- 6 Hydrate block if needed with ammonia. Dry block after hydration.
- 7 Begin sectioning the block. Grasp the first section lightly. Continue sectioning until you have a ribbon of sections (~6-10 sections).

- 8 Use camel hair brush to lightly remove ribbon by brushing from one side to the other. Alternatively, lightly separate the section attached to the block from the remainder of the sections by scoring with a razor blade.

Note

This will allow you to keep a non-wrinkled section attached to the front of the next ribbon.

- 9 Transfer the ribbon to a  37 °C water bath. Burst bubbles with a syringe needle. Score between sections with a razor blade.
- 10 Gently spring apart sections with forceps and mount them in order on glass slides. Transfer the slides to an upright rack to allow the water to wick off of the sections before transferring to a metal rack.
- 11 Continue sectioning (Steps 6-10) until the block is depleted of tissue.

- 12 Transfer slides to oven to dry at  42 °C for  02:00:00 or more or  60 °C for  00:20:00 2h 20m more.

