




Version 2

Sep 14, 2021

Immunofluorescence of Galectin-3 Puncta after lysosomal damage with LLoMe V.2

Vinay V. Eapen¹, Sharan Swarup¹, Melissa Hoyer¹, Harper JW²¹dx.doi.org/10.17504/protocols.io.bxgypjxw; ²wade_harper@hms.harvard.edu

1 Works for me

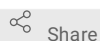

Share

dx.doi.org/10.17504/protocols.io.bx8pprvn

Harper JW

DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

ABSTRACT

Lysophagy-the selective elimination of damaged lysosomes by the autophagy pathway-is a critical housekeeping mechanism in cells. This pathway surveils lysosomes and selectively demarcates terminally damaged lysosomes for elimination. Among the most upstream signaling proteins in this pathway are the glycan binding proteins- Galectins-which recognize N and O linked glycan chains on the luminal side of transmembrane lysosomal proteins. These glycosyl modifications are only accessible to galectin proteins upon extensive lysosomal membrane rupture and serve as a sensitive measure of lysosomal damage and eventual clearance by selective autophagy. Indeed, prior work has shown that immunofluorescence of Galectin-3 serves as a convenient proxy for lysophagic flux in tissue culture cells (Aits et al., 2015; Maejima et al., 2013). Here we describe our method for monitoring galectin-3 puncta clearance as a proxy for turnover of damaged lysosomes via immunofluorescence and confocal imaging.

ATTACHMENTS

[d7fzbhjd.docx](#) Immunofluorescence of
Galectin-3 Puncta after
lysosomal damage with
LLoMe.cproj

DOI

dx.doi.org/10.17504/protocols.io.bx8pprvn

PROTOCOL CITATION

Vinay V. Eapen, Sharan Swarup, Melissa Hoyer, Harper JW 2021. Immunofluorescence of Galectin-3 Puncta after lysosomal damage with LLoMe. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bx8pprvn>
Version created by Harper JW

KEYWORDS

Immunofluorescence, Galectin-3, Lysosomal damage, LLoMe

LICENSE

— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 14, 2021

LAST MODIFIED

Sep 14, 2021

PROTOCOL INTEGER ID

53231

MATERIALS TEXT

Materials:

A	B	C
REAGENT or RESOURCE	SOURCE	IDENTIFIER
Chemicals		
LLoMe (L-Leucyl-L-Leucine methyl ester (hydrochloride))	Cayman Chemical	16008
Dulbecco's MEM (DMEM), high glucose, pyruvate	GIBCO / Invitrogen	11995
Phosphate Buffered Saline 1X	Corning	21-031-CV
Fetal Bovine Serum	Fisher	SH3008003
Bovine Serum Albumin	Gold biotechnology	A-420-250
paraformaldehyde	Electron Microscopy Sciences	15710
Triton-X	Sigma	T8787
Antibodies		
LAMP1 (D401S) Mouse mAb	Cell Signaling Technology	15665S
Galectin-3/ LGALS3 (M3/38) Rat mAb	Santa-Cruz	sc-23938
Software		
Cell Profiler	CellProfiler v4.0.6	https://cellprofiler.org/
Fiji	ImageJ V.2.0.0	https://imagej.net/software/fiji/
Metamorph	Metamorph v	https://www.moleculardevices.com/products/cellular-imaging-systems/acquisition-and-analysis-software/metamorph-microscopy#gref

Chemicals:

[L-Leucyl-L-Leucine methyl ester \(hydrochloride\)](#) **Cayman Chemical**

Company Catalog #16008

[DMEM, high glucose, pyruvate](#) **Thermo**

Fisher Catalog #11995065

[Phosphate Buffered Saline](#)

[\(PBS\)](#) **Corning Catalog #MT21-031-CV**

 **Paraformaldehyde Electron Microscopy**

Sciences Catalog #15710

 **Triton X-100 Sigma**

Aldrich Catalog #T8787

Antibodies:

 **LAMP1 (D401S) Mouse mAb #15665 Cell Signaling**

Technology Catalog #15665S

 **Anti-galectin-3 Antibody (M3/38): sc-23938 Santa Cruz**

Biotechnology Catalog #sc-23938

DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

Immunofluorescence of Galectin-3 Puncta after lysosomal damage with LLoMe

4h 12m

- 1 Plate the cells (selected by investigator) into 12 well glass bottom dishes (No. 1.5, 14 mm glass diameter, MatTek) are grown to 50-70% confluency in media.

For HeLa cells, we use Dulbecco's MEM (DMEM), high glucose (4500 mg/L), pyruvate (100 mg/L) supplemented with 10% fetal bovine serum.

- 2 Treat the cells with **[M]500 Micromolar (μM) – [M]1 Milimolar (mM)** of LLoMe for  **01:00:00** . 1h

Note: The exact dosage varies with cell line and should be determined empirically depending on the line used. This dose range has been tested extensively in Hela cells, and routinely generated lysophagic flux.

- 3 Remove LLoMe containing media from the cells and replace with fresh media not containing LLoMe.

- 4 

15m

After 10h, Wash the cells one time with phosphate buffered saline (PBS) and then fix with 4% paraformaldehyde in PBS for ⌚ 00:15:00 at 🌡 Room temperature .

5 


15m

Remove 4% paraformaldehyde in PBS, wash the cells once with PBS, and solubilize cells with 0.1% triton-X in PBS for ⌚ 00:15:00 at 🌡 Room temperature .

6 Block the cells are for ⌚ 00:30:00 at 🌡 Room temperature with sterile filtered blocking buffer (1% bovine serum albumin, 0.1% triton-X in PBS). ^{30m}

7 

1h 1m

Add primary antibodies to blocking buffer at 1:300 and then spun down for ⌚ 00:01:00 at 🌀 10000 x g . Remove the blocking buffer is completely and then apply the antibody in blocking buffer to the cells ( 100 µl applied to the center of the well where the glass coverslip is attached) for ⌚ 01:00:00 at 🌡 Room temperature .

8 

5m

Wash the cells 4 times with PBS (⌚ 00:05:00 for each wash).

9 

1h 1m

Add fluorescently conjugated secondary antibodies to blocking buffer at 1:300 and then spun down for ⌚ 00:01:00 at 🌀 10000 x g . Remove the blocking buffer is completely and then apply the antibody in blocking buffer to the cells (100µL applied to the center of the well where the glass coverslip is attached) for ⌚ 01:00:00 at 🌡 Room temperature .

10 

5m

Wash the cells 4 times with PBS (⌚ 00:05:00 for each wash) and left in PBS.

11 

Image the cells at 🌡 Room temperature using a Yokogawa CSU-X1 spinning disk confocal on a Nikon Ti-E inverted microscope at the Nikon Imaging Center in Harvard Medical School. Use Nikon Perfect Focus System to maintain cell focus over time. Equip the microscope with a Nikon Plan Apo 40x/1.30 N.A or 100x/1.40 N.A objective lens.

445nm (75mW), 488nm (100mW), 561nm (100mW) & 642nm (100mW) laser lines are controlled by AOTF. All

images are collected with a Hamamatsu ORCA-ER cooled CCD camera (6.45 μm^2 photodiode) with MetaMorph image acquisition software.

- 12 Detect Galectin-3 puncta using CellProfiler with the same pipeline applied for each condition (see attached CellProfiler pipeline).

Each cell area is first defined using a “identify primary objects” module that included objects 200 to 1000 pixels units, and each puncta is marked using a “identify primary objects” module that included objects 2 to 20 pixels units both with an optimized “robust background” threshold. Each cell for each condition is thresholded in the same way with a consistent pipeline. Object size and shape is measured, and each punctum is related to its respective cell to yield a puncta per cell readout. Each channel z series are brightness and contrast adjusted equally and then converted to rgb for publication using FIJI software.

- 13 Each channel z series are brightness and adjust contrast equally and then convert to rgb for publication using FIJI software.