




Protocol to isolate and fix nuclei from flash frozen mouse kidneys for IGVF

Elisabeth Rebboah¹

¹University of California, Irvine

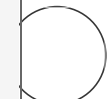

Elisabeth Rebboah
University of California, Irvine

AUG 17, 2023

ABSTRACT

This protocol describes isolation of nuclei from 10 week old mouse kidneys from 8 founder strains (B6J, AJ, 129S1J, NZOJ, WSBJ, NODJ, PWKJ, and CASTJ), preparation of a single nucleus suspension, and fixation for single nucleus RNA-seq using Parse Biosciences. We process 1 rep from each strain per day; e.g. female rep 1 across all 8 strains. The main products we use are Parse Biosciences Nuclei Fixation Kit (v2) and Miltenyi Biotec's gentleMACS Octo Dissociator with accessories. This protocol takes about 3.5 hours from start to finish.

The results are 2 aliquots of fixed single-nucleus suspensions for Parse per each of the 8 samples.

OPEN ACCESS



Protocol Citation: Elisabeth Rebboah 2023. Protocol to isolate and fix nuclei from flash frozen mouse kidneys for IGVF. [protocols.io](https://protocols.io/view/protocol-to-isolate-and-fix-nuclei-from-flash-froz-cyubxwsn) <https://protocols.io/view/protocol-to-isolate-and-fix-nuclei-from-flash-froz-cyubxwsn>

GUIDELINES

- We recommend using a 5 ml pipette for aspirations and resuspensions > 1 ml.
- Record everything in the [IGVF spreadsheet](#), "Samples into experiment" tab.

MATERIALS

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Aug 17, 2023

Last Modified: Aug 17, 2023

PROTOCOL integer ID: 86627

| Name | Manufacturer | Cat # |
|-----------------------------|---------------------|-------------|
| Nuclei Fixation Kit v2 | Parse Biosciences | ECF2003 |
| Nuclei Extraction Buffer | Miltenyi Biotec | 130-128-024 |
| RNase Inhibitor, murine | New England Biolabs | M0314L |
| PBS | HyClone | SH30256.02 |
| 7.5% BSA | Life Technologies | 15260037 |
| gentleMACS C Tube | Miltenyi Biotec | 130-093-237 |
| gentleMACS Octo Dissociator | Miltenyi Biotec | 130-095-937 |

| Name | Manufacturer | Cat # |
|--------------------------------|-------------------|-------------|
| MACS SmartStrainers (30 µm) | Miltenyi Biotec | 130-098-458 |
| NucBlue Fixed Cell ReadyProbes | Thermo Fisher | R37606 |
| Hemocytometer | Fisher Scientific | 02-671-51B |
| Mr. Frosty | Sigma-Aldrich | 635639 |

Reagents/equipment, manufacturer and catalog number

| Name | reagent | Volume for 8 samples | Final concentration |
|--|-------------------------------------|-------------------------------------|---------------------|
| Lysis buffer | Nuclei Extraction Buffer | 40 ml | NA |
| | 40 U/ul RNase inhibitor | 200 ul | 0.2 U/ul |
| NB + BSA + RNase inhibitor (make 2 aliquots) | Nuclei Buffer (Parse Biosciences) | 3.15 ml | NA |
| | 7.5% BSA | 350 ul | 0.75% |
| | RNase inhibitor (Parse Biosciences) | RNase inhibitor (Parse Biosciences) | 44.1 ul |
| RSB | PBS | 42 ml | |
| | 7.5% BSA | 560 ul | 0.1% |
| | RNase inhibitor | 210 ul | 0.2 U/ul |
| SHARE-RSB | 1 M HEPES pH 7.3 | 150 ul | 10 mM |
| | 5 M NaCl | 30 ul | 10 mM |
| | 1 M MgCl ₂ | 45 ul | 3 mM |
| | 10% Tween-20 | 150 ul | 0.1% |
| | H ₂ O | 14.625 ml | |
| | 7.5% BSA | 80.26 ul | 0.04% |
| | 5% digitonin | 30 ul | 0.01% |
| | Enzymatics RI | 37.5 ul | 0.1 U/ul |
| | SUPERase RI | 18.75 ul | 0.025 U/ul |

| | Name | reagent | Volume for 8 samples | Final concentration |
|--|------|------------|----------------------|---------------------|
| | | Yeast tRNA | 150 ul | 100 ug/ml |

Buffers

Setup

- 1 Label tubes.
- 2 Pre-chill centrifuge to 4C.
- 3 Prepare ice buckets.
- 4 Prepare lysis buffer in a 50 ml conical tube on ice. Distribute 2.5 ml into 8 gentleMACS C Tubes on ice. Add RNase inhibitor to the lysis buffer aliquot the day of the experiment.
- 5 Prepare RSB in a 50 ml conical tube on ice. Add NEB RNase inhibitor the day of the experiment.
- 6 Prepare NB + BSA + RNase inhibitor. Add RNase inhibitor included in Parse Biosciences fixation kit the day of the experiment.
- 7 Prepare 2.5 ml Nuclei Buffer + RNase inhibitor for final resuspension. Add 31.5 ul Parse RNase inhibitor to 2.5 ml nuclei buffer.

- 8 Thaw components of 2 Parse Biosciences Nuclei Fixation v2 kits at room temperature, then place on ice.
- 9 Distribute 20 ul NucBlue Fixed Cell ReadyProbes into 16 PCR strip tubes for cell counting. Need 8 tubes for counting after nuclei extraction, and another 8 tubes for final fixed nuclei.

Tissue lysis and nuclei extraction

- 10 Keep flash frozen tissue samples on dry ice until lysis.
- 11 Drop whole frozen tissue into a chilled gentleMACS C Tube with 2.5 ml lysis buffer. Close tubes firmly and invert immediately, ensuring tissue is not stuck to the bottom or side. Keep tubes on ice and proceed immediately to dissociation. There should be 2 kidneys.
- 12 Run the gentleMACS Program 4C_nuclei_1 on the Octo Dissociator (~5 minutes).
- 13 Remove tubes, ensuring tissue did not get stuck on the sides, and spin down in a 4C centrifuge for ~10 seconds to bring liquid to the bottom, then place tubes back on ice.
- 14 Filter nuclei suspension through 70 um MACS SmartStrainer into a 5 ml tube. Fit a tube rack in ice for extra stability while filtering.
- 15 Wash 70 um MACS SmartStrainer with 2 ml additional lysis buffer. Add 2 ml to C tubes, cap, and swish to recover any nuclei stuck to the sides and cap of the C tubes, then wash the strainer.

- 16 Discard strainer and centrifuge the 4.5 ml nuclei suspension at 4C, 350g for 5 minutes.
- 17 Discard supernatant and resuspend nuclei pellet in 3 ml RSB. Filter nuclei suspension through 30 um MACS SmartStrainer into a 5 ml tube.
- 18 Take 200 ul and add to a 5 ml tube that contains another 1.8 ml RSB for a 1:10 dilution (same as male gonads and similar to liver protocols).
- 19 Count nuclei. Use 1:11 dilution factor, 2 ul + 20 ul dye.
- 20 Centrifuge 4 million nuclei at 4C, 500g for 5 minutes and remove supernatant.
- 21 Resuspend pellet in 750 ul NB-BSA + RNase inhibitor.

Parse nuclei fixation

- 22 Filter the nuclei suspension (750 ul) NB-BSA + RNase inhibitor through a 40 um strainer (provided in Parse Biosciences kit) into a new 5 ml tube.
- 23 Add 250 uL Nuclei Fixation Solution and mix 3 times. Do not over-mix.
- 24 Incubate nuclei for 10 minutes on ice. Set 1 P200 pipette to 80 ul and keep the P1000 at 250 ul.

- 25** Add 80 uL Nuclei Permeabilization Solution and mix by pipetting 3 times with the P1000 still set to 250 uL. Do not over-mix.
- 26** Incubate 3 minutes with nuclei on ice.
- 27** Add 4 ml Nuclei Neutralization Solution and invert the tube once to mix.
- 28** Centrifuge at 4C, 750g for 10 minutes.
- 29** Aspirate and discard supernatant.
- 30** Resuspend the samples in 300 ul Nuclei Buffer with RNase inhibitor without BSA and filter through a 40 um filter into a labeled 1.5 ml tube.
- 31** Count nuclei. Use a 1:11 dilution factor, e.g. 2 ul + 20 ul dye.
- 32** Add Nuclei DMSO: For 300 ul samples: add 5 ul and gently flick tubes to mix. One minute later, add another 5 ul and flick to mix, then after another minute add a final 5 ul for a total volume of 15 ul. Mix by gently pipetting 5x with a P200 set to 150 ul.

- 33 Split nuclei suspension into 2 aliquots of equal volume in labeled 1.5 ml tubes.
- 34 Place tubes in a Mr. Frosty at -80C. The next day, move tubes to boxes in -80C racks.
- 35 Move leftover nuclei suspension to labeled 2 ml tubes and spin at 4C, 750g for 5 minutes. Remove supernatant and flash-freeze nuclei in liquid nitrogen as dry pellets. Store at -80C.

SHARE-seq nuclei fixation

- 36 Set aside 1 million nuclei for each of the 8 samples in RSB and spin down at 4C, 750g for 5 minutes.
- 37 Remove supernatant and resuspend nuclei pellet in 1 ml room temperature SHARE-RSB. Transfer tube to a room temperature rack.
- 38 At RT, add 13.34 ul of methanol-free formaldehyde (16% stock solution). Final concentration for nuclei: 0.2%. Close tube and nutate cells at RT for 5 minutes.
- 39 To quench fixation, per reaction, add 56.1 ul fresh 2.5M Glycine (0.94g per 5 ml stock), 50 ul of 1M Tris pH 8.0, 13.3ul of 7.5% BSA, and mix using a pipette. Incubate on ice for 10 minutes.
- 40 Spin 750g, 4C, 5 minutes. Gently remove supernatant.
- 41 Add 200 ul of SHARE-RSB and gently resuspend pellet. Store on ice until all samples are

completed.

- 42** Pool 200 ul of resuspended nuclei from all 8 founders into 1 labeled 2 ml tube.
- 43** Spin 1,000g, 4C, 10 minutes. Gently remove supernatant. Remove all fluid and freeze at -80C as a dry pellet.

Storage of leftover nuclei

- 44** Move remaining nuclei in RSB in 5 ml tubes on ice to new labeled 2 ml tubes.
- 45** Spin 750g, 4C, 5 minutes.
- 46** Remove all supernatant and flash-freeze nuclei as a dry pellet in liquid nitrogen. Store at -80C.