



NOV 21, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.ewov1yr6kvr2/v1

Protocol Citation: Monique Copeland 2023. 4% Paraformaldehyde in .1M PB preparation. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.ewov1yr6kvr2/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working


Created: Aug 16, 2019

Last Modified: Nov 21, 2023

4% Paraformaldehyde in .1M PB preparation

Monique Copeland¹

¹HHMI Janelia Research Campus

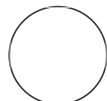

Monique Copeland
 HHMI Janelia Research Campus

ABSTRACT

For perfusion solution.

MATERIALS

MATERIALS

Paraformaldehyde Catalog #19200

PROTOCOL MATERIALS

Paraformaldehyde Catalog #19200 Materials, Step 4

Sodium Hydroxide, 10.0N Thomas Scientific Catalog #C871W32 Step 3

Corning™ Disposable Vacuum Filter/Storage Systems Fisher Scientific Catalog #09-761-5

Step 8

Falcon Tube (50 mL) Fischer Scientific Step 9

Sodium Phosphate Dibasic Dihydrate (USP/FCC/EP/BP), Fisher Chemical Fisher Scientific Catalog #S472-500

Step 5.1

Sodium Phosphate Monobasic Anhydrous (Colorless-to-White Crystals), Fisher BioReagent Fisher Scientific Catalog #BP329-500

Step 5.2

Paraformaldehyde prep


1 Add 400 mL of MilliQ water to a 1000 mL beaker

Equipment

Milli-Q® Integral Water Purification System for Ultrapure Water	NAME
filter	TYPE
Millipore Sigma	BRAND
ZRXQ010WW	SKU
http://www.emdmillipore.com/US/en/product/Milli-Q-Integral-10-Water-Purification-System,MM_NF-ZRXQ010WW	LINK

2 Add stir bar in beaker and set on hot plate at 50 °C

3 Add 5 drops (500ul) of 10N NaOH (from 60ml dropper)

 Sodium Hydroxide, 10.0N Sigma
 Aldrich Catalog #C871W32


(this makes the solution basic to help dissolve the paraformaldehyde)

4 Weigh 20 g of Paraformaldehyde and add to MilliQ water.
 (Stir about 01:00:00 or until all parafromaldehyde is dissolved and solution is clear)


 Paraformaldehyde Sigma
 Aldrich Catalog #19200

5 Make 1 Molarity (M) Phosphate buffer at PH 7.4

5.1 Prepare 1 Molarity (M) Sodium Phosphate Dibasic Dihydrate (Store stock at Room temperature .)

 Sodium Phosphate Dibasic Dihydrate (USP/FCC/EP/BP), Fisher Chemical Sigma Aldrich Catalog #S472-500

5.2 Prepare 1 Molarity (M) Sodium Phosphate Monobasic Anhydrous (store stock at Room temperature)

 Sodium Phosphate Monobasic Anhydrous (Colorless-to-White Crystals), Fisher BioReagent Sigma Aldrich Catalog #BP329-500


5.3 To make 100 mL 1M PB at PH 7.4 add 19 mL Monobasic plus 81 mL Dibasic (store as stock at Room temperature)


5.4 Add 50 mL 1M PB PH 7.4 to paraformaldehyde solution.


6 Test PH and adjust solution to final desired pH using HCl or NaOH.

7 Top off to 500ml with MilliQ water.

8 Filter solution through .2um CA member filter.

 Corning™ Disposable Vacuum Filter/Storage Systems Sigma
Aldrich Catalog #09-761-5

9 Lable container with '4% paraformaldehyed in .1M PB' and date. Store at  4 °C for up to 1 week.
(Can also be aliquouted into

 Falcon Tube (50 mL) Sigma
Aldrich

and stored at  -20 °C for 3 months.)