



SEP 23, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocol.s.io.eq2lyj33plx9/v1

Protocol Citation: Elias Adriaenssens 2023. SDS-PAGE and western blot analysis. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.eq2lyj33plx9/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Jun 28, 2023

Last Modified: Sep 23, 2023

PROTOCOL integer ID:
 84147

Keywords: SDS-PAGE, western blot analysis

SDS-PAGE and western blot analysis

Elias Adriaenssens¹

¹Sascha Martens lab, University of Vienna, Max Perutz Labs - Vienna

ASAP Collaborative Research Network

Sascha Martens lab, University of Vienna, Max Perutz Labs - Vienna



Elias Adriaenssens

Sascha Martens lab, University of Vienna, Max Perutz Labs - ...

ABSTRACT

This protocol describes SDS-PAGE and western blot analysis.

ATTACHMENTS

[758-1926.pdf](#)

MATERIALS

Materials

- 100 mM DTT
- PageRuler Prestained protein marker (Thermo Fisher)
- nitrocellulose membranes (RPN132D, GE Healthcare)
- Mini Trans-Blot Cell (Bio-Rad)
- 0.1% Tween 20
- horseradish peroxidase (HRP)-coupled antibodies

RIPA buffer

A	B
Tris-HCl pH 8.0	50 mM
NaCl	150 mM
sodium deoxycholate	0.50%
SDS	0.10%
NP-40	1%

⊗ cOmplete™, Mini, EDTA-free (Protease Inhibitor) Roche Catalog ##11836170001


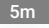
⊗ Roche PhosSTOP™ Merck MilliporeSigma (Sigma-Aldrich) Catalog #4906837001

⊗ Pierce™ Detergent Compatible Bradford Assay Kit Thermo Fisher Catalog #23246


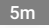
⊗ SuperSignal™ West Femto Maximum Sensitivity Substrate Thermo Fisher Catalog #34096

SDS-PAGE and western blot analysis

10m




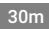
1 For SDS-PAGE and western blot analysis, collect cells by trypsinization and subsequent centrifugation at  300 x g, 4°C, 00:05:00 . 



2 Wash the cell pellets in PBS and centrifuge once more at  300 x g, 4°C, 00:05:00 . 






3 Remove the supernatant and lyse the cell pellets in RIPA buffer supplemented by cOmplete EDTA-free protease inhibitors (11836170001, Roche) and phosphatase inhibitors (Phospho-STOP, 4906837001, Roche).

4 After incubating in RIPA buffer for  00:20:00  On ice , clear the samples by centrifugation at  20000 x g, 4°C, 00:10:00 . 





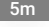



5 Collect the soluble supernatant fraction and measure the protein concentrations using the Pierce Detergent Compatible Bradford Assay Kit (23246, Thermo Fisher).

6 Then, adjust the samples for equal loading and mix it with 6x protein loading dye, supplemented with 100 mM DTT. Then, boil it for  00:05:00 at  95 °C . 






7 Load the samples on 4-12% SDS-PAGE gels (NP0321BOX, NP0322BOX, or NP0323BOX, Thermo Fisher) with PageRuler Prestained protein marker (Thermo Fisher).






8 Transfer the proteins onto nitrocellulose membranes (RPN132D, GE Healthcare) for  01:00:00 at  4 °C using the Mini Trans-Blot Cell (Bio-Rad). 

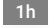9 After the transfer, block the membranes with 5% milk powder dissolved in PBS-Tween (0.1% Tween 20) for  01:00:00 at  Room temperature . 



10 Incubate the membranes  Overnight at  4 °C with primary antibodies dissolved in the blocking buffer and wash. 



- 10.1** Wash for  00:05:00 (1/3) 5m
- 10.2** Wash for  00:05:00 (2/3) 5m
- 10.3** Wash for  00:05:00 (3/3) 5m
- 10.4** Then, incubate with species-matched secondary horseradish peroxidase (HRP)-coupled antibodies diluted 1:10,000 in blocking buffer for  01:00:00 at  Room temperature 1h.

11 Afterwards wash the membranes three times with PBS-T and further process for western blot detection.



- 12** Incubate the membranes with SuperSignal West Femto Maximum Sensitivity Substrate (34096, Thermo Fisher) and image it with a ChemiDoc MP system (Bio-Rad).



- 13** Analyze the images with ImageJ [57].

