


PMN- 04 - Measurement of Intracellular Reactive Oxygen Species (ROS) levels in Human PMN

Marco Cosentino¹, Elisa Storelli¹, Alessandra Luini¹, Massimiliano LM Legnaro¹, Emanuela Rasini¹, Marco Ferrari¹, Franca Marino¹

¹Center for Research in Medical Pharmacology, University of Insubria (Varese, Italy)

1 Works for me dx.doi.org/10.17504/protocols.io.biywkfxe

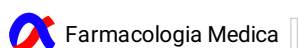

ABSTRACT

List of published work using this protocol:

- Immune Profile of Hyaluronic Acid Hydrogel Polyethylene Glycol Crosslinked: An in Vitro Evaluation in Human Polymorphonuclear Leukocytes. Franca Marino, Marco Cosentino, Massimiliano Legnaro, Alessandra Luini, Julia Sigova, Roberto Mocchi, Trello Lotti, Nicola Zerbinati. *Dermatol Ther.* 2020 Apr 10; e13388. doi: 10.1111/dth.13388. Online ahead of print.
- A Novel Standardized Cannabis sativa L. Extract and Its Constituent Cannabidiol Inhibit Human Polymorphonuclear Leukocyte Functions. Alex Mabou Tagne, Franca Marino, Massimiliano Legnaro, Alessandra Luini, Barbara Pacchetti, Marco Cosentino. *Int J Mol Sci.* 2019 Apr 13;20(8):1833. doi: 10.3390/ijms20081833.
- β 2 -Adrenoceptors Inhibit Neutrophil Extracellular Traps in Human Polymorphonuclear Leukocytes. Franca Marino, Angela Scanzano, Laura Pulze, Monica Pinoli, Emanuela Rasini, Alessandra Luini, Raffaella Bombelli, Massimiliano Legnaro, Magda de Eguileor, Marco Cosentino. *J Leukoc Biol.* 2018 Sep;104(3):603-614. doi: 10.1002/JLB.3A1017-398RR. Epub 2018 Apr 18.
- Adrenergic Modulation of Migration, CD11b and CD18 Expression, ROS and interleukin-8 Production by Human Polymorphonuclear Leukocytes. Angela Scanzano, Laura Schembri, Emanuela Rasini, Alessandra Luini, Jessica Dallatorre, Massimiliano Legnaro, Raffaella Bombelli, Terenzio Congiu, Marco Cosentino, Franca Marino. *Inflamm Res.* 2015 Feb;64(2):127-35. doi: 10.1007/s00011-014-0791-8. Epub 2015 Jan 6.
- The Essential Oil of Bergamot Stimulates Reactive Oxygen Species Production in Human Polymorphonuclear Leukocytes. Marco Cosentino, Alessandra Luini, Raffaella Bombelli, Maria T Corasaniti, Giacinto Bagetta, Franca Marino. *Phytother Res.* 2014 Aug;28(8):1232-9. doi: 10.1002/ptr.5121. Epub 2014 Jan 23.
- Immunomodulatory Properties of Achyrocline Satureioides (Lam.) D.C. Infusion: A Study on Human Leukocytes. Marco Cosentino, Raffaella Bombelli, Elena Carcano, Alessandra Luini, Franca Marino, Francesca Crema, Federico Dajas, Sergio Lecchini. *J Ethnopharmacol.* 2008 Mar 28;116(3):501-7. doi: 10.1016/j.jep.2007.12.014. Epub 2008 Jan 6.

DOI

dx.doi.org/10.17504/protocols.io.biywkfxe

PROTOCOL CITATION

Marco Cosentino, Elisa Storelli, Alessandra Luini, Massimiliano LM Legnaro, Emanuela Rasini, Marco Ferrari,

Franca Marino 2020. PMN- 04 - Measurement of Intracellular Reactive Oxygen Species (ROS) levels in Human PMN. [protocols.io](https://doi.org/10.17504/protocols.io.biywkfxe)
[dx.doi.org/10.17504/protocols.io.biywkfxe](https://doi.org/10.17504/protocols.io.biywkfxe)



LICENSE

———— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 24, 2020

LAST MODIFIED

Jul 24, 2020

PROTOCOL INTEGER ID

39670

MATERIALS

NAME	CATALOG #	VENDOR
NaCl	S9625	Sigma Aldrich
C-DCFH-DA	C2938	Invitrogen - Thermo Fisher
KCl	P9541	Sigma
MgSO4	1.05886	Sigma
Glucose	1.08337	Sigma
CaCl2	1.02382	Sigma
Hepes	H7006	Sigma

MATERIALS TEXT

Dubnoff bath

EQUIPMENT

NAME	CATALOG #	VENDOR
Spectrofluorimeter	LS50B	Perkin Elmer
Centrifuge	5804	Eppendorf Centrifuge

- 1 Isolate PMN using either protocol PMN- 01a or PMN 01b.
- 2 Suspend isolated PMN 1×10^6 cells/ml in **HBSS/HEPES medium**.





SOLUTION- 12 – HBSS-HEPES for ROS measurement in PMN
by [Farmacologia Medica](#)


- 3 Add 2×10^6 cells in 2 mL of **HBSS/HEPES** for each tube.

4 Incubate for ⌚ **01:00:00** with 2 Micromolar (μM) of **C-DCFH-DA** in the **dark** at **37 °C** (Dubnoff Bath).



5 Wash twice with **2 mL** of **HBSS/HEPES**.

Each time centrifuge at **400 x g, Room temperature 00:05:00**

  SOLUTION- 12 – HBSS-HEPES for ROS measurement in PMN
by **Farmacologia Medica**

 Centrifuge
Eppendorf 5804

6 Aspirate supernatant and resuspend with **2 mL** of **HBSS/HEPES**.

  SOLUTION- 12 – HBSS-HEPES for ROS measurement in PMN
by **Farmacologia Medica**

7 Recover PMN in cuvette of 2 ml in the **dark** at **37 °C** (Dubnoff Bath).

8 Measurements of fluorescence must be performed by means of a **Spectrofluorimeter with excitation wavelength of 488 nm**.



Spectrofluorimeter
Perkin Elmer LS50B

9

Cuvette must be maintained at 37°C by means of a thermostatic bath, that will be turned on before the starting ROS measurement.

- 10 Fluorescence emission must be collected at **525 nm** and recorded in a specific data table that will be signed for each specific experimental plan (time for recording and complete time analysis), that will be signed by the researcher.
- 11 Intracellular ROS levels must be expressed as fluorescence intensity, in arbitrary units (AU) or according to the specific experimental plan.