

MAR 08, 2024

# 🌐 DOH Workshop Protocol Part 1: Purification of High- Molecular Weight Genomic DNA from Gram-Negative Bacteria (MagAttract HMW DNA)

📁 In 2 collections

Vesa Qarkaxhija<sup>1</sup>, Bryan Wee<sup>1</sup>, Natalie Ring<sup>1</sup>

<sup>1</sup>The Roslin Institute, University of Edinburgh



Bryan Wee  
The University of Edinburgh

OPEN  ACCESS



## ABSTRACT

This protocol enables extraction of high-molecular-weight DNA from gram-negative bacterial cultures.

## ATTACHMENTS

[pyxtbpxx.pdf](#)

## DOI:

[dx.doi.org/10.17504/protocols.io.5jyl8p3p8g2w/v1](https://dx.doi.org/10.17504/protocols.io.5jyl8p3p8g2w/v1)

**Protocol Citation:** Vesa Qarkaxhija, Bryan Wee, Natalie Ring 2024. DOH Workshop Protocol Part 1: Purification of High- Molecular Weight Genomic DNA from Gram-Negative Bacteria (MagAttract HMW DNA).

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.5jyl8p3p8g2w/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

Created: Mar 05, 2024

## GUIDELINES


Last Modified: Mar 08, 2024

### Important points throughout protocol:

PROTOCOL integer ID: 96352



Be gentle when pipetting, too much fast pipetting/ extended vortexing can shear the DNA which will result in poorer sequencing outcomes.

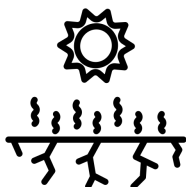

Do not let the beads dry after removal of supernatant. This can be avoided by keeping Eppendorf tubes closed if beads are not submerged.

**Note:** When bead pellet is moist, it appears shiny. As they start to dry, the shine reduces, and cracks start to form.



When working with beads ensure they are thoroughly mixed before using. This can be achieved by vortexing for at least 10 seconds before use and vortexing between use to prevent beads settling.

## MATERIALS

### Starting Materials

- 1.5 mL of overnight bacterial cell culture.

### Equipment

- P1000 pipette (Micropipette with 100 µL – 1000 µL range)
- P200 pipette (Micropipette with 20 µL – 200 µL range)
- P20 pipette (Micropipette with 2 µL – 20 µL range)
- P10 pipette (Micropipette with 0.5 µL – 10 µL range)
- Thermomixer (for heating, cooling and mixing 1.5ml tubes)
- Benchtop centrifuge for 1.5 µL (up to 5000 Gs) OR Bento Lab
- Minicentrifuge (for brief spin downs of 1.5ml eppendorfs and 0.2ml PCR tubes) OR Bento lab with improvised adaptor made from 0.5ml eppendorf within a 1.5ml eppendorf with caps cut off.
- PCR thermoblock ( 21 °C - 80 °C required) OR Bento Lab
- DNA fluorometer (Promega Quantus or Themofisher Qubit)
- Vortex mixer

### Consumables

- P1000 filter pipette tips (with 100 µL – 1000 µL range)
- P200 filter pipette tips (with 20 µL – 200 µL range)
- P20 filter pipette tips (with 2 µL – 20 µL range)
- P10 filter pipette tip (with 0.5 µL – 10 µL range)
- 1.5 ml Eppendorf DNA LoBind tubes
- 0.2 ml thin-walled PCR tubes
- Absolute ethanol (>96%)
- Distilled water ( 1400 mL per sample)
- HMW Magattract Kit
  1. Nuclease Free Water
  2. Buffer MW1 (Ethanol must be added)
  3. Buffer PE (Ethanol must be added)
  4. MagAttract Suspension G
  5. Proteinase K



6. RNAase A (  100 undetermined )


7. Buffer ATL

#### BEFORE START INSTRUCTIONS

- If precipitate is formed in Buffer ATL, dissolve by incubating at 37°C with occasional shaking. The presence of precipitate can result in incomplete lysis.
- Prepare (but do not start) the thermomixer to 56°C at 900rpm for 30min for the lysis step.

### Purification of High- Molecular Weight Genomic DNA from Gram-Negative ... 59m 12s

1 Pellet  1.5 mL of the Gram negative bacterial cell culture by spinning at  5000 x g for 3m

 00:03:00 .

#### Note

When removing tube from centrifuge, be careful not to shake the tube too much which can dislodge the pellet.

2

Remove and discard supernatant without disturbing the pellet. Use a P200/P20 pipette if there is a small amount of supernatant remaining.



#### Note

The cell pellet can be stored at  $-30\text{ }^{\circ}\text{C}$  to  $-15\text{ }^{\circ}\text{C}$  or  $-90\text{ }^{\circ}\text{C}$  to  $-60\text{ }^{\circ}\text{C}$  for future use, or can be used immediately.

3

Resuspend the bacterial pellet in  $180\text{ }\mu\text{L}$  Buffer ATL, add  $20\text{ }\mu\text{L}$  Proteinase K and mix by flicking/stirring the tube.



#### Note



Pipette gently, be careful not to shear DNA

4

Incubate for  $00:30:00$  on thermomixer at  $56\text{ }^{\circ}\text{C}$  shaking at 900rpm.




30m



### Things to do while waiting:

- Ensure that Buffers MW1 and PE were prepared according to instructions on them (i.e. adding appropriate amount of Ethanol).

5

Add  4  $\mu\text{L}$  RNase A to the sample, mix by pulse vortexing or by tapping the tube several times, and incubate for  00:02:00 at  Room temperature .

2m




#### Note




Pipette gently, be careful not to shear DNA

6

Set the thermomixer to  21  $^{\circ}\text{C}$  to give it time to cool for next steps.

7

DOH Workshop Protocol Part 1 Fully mix MagAttract Suspension G. Vortex the MagAttract Suspension G vigorously until thoroughly mixed ~  00:00:10 .

10s








#### Note

It is important that beads are fully mixed and re-mixed between every few samples.



Shake well before use

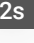8

Add  15  $\mu$ L MagAttract Suspension G and  280  $\mu$ L Buffer MB to the sample. Mix by pulse vortexing  2s  
( 00:00:01 -  00:00:02 quick presses x 3).






#### Note




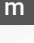
Pipette gently, be careful not to shear DNA

9

Place the tube onto the thermomixer and incubate at  21  $^{\circ}$ C for  00:03:00 at 1400 rpm..  3m



10


Place tubes onto magnetic rack and wait until beads pellet (~  00:01:00 ).  1m

- 11 While leaving tubes on the magnet, carefully remove supernatant without disturbing the beads using a P1000.

#### Note





Pipette gently, be careful not to shear DNA

- 12 Remove tube from magnetic rack and add  700 µL Buffer MW1 to the tube, flick the tube until bead pellet is unstuck from tube wall.




#### Note

MW1 Wash 1

- 13 Place the tube onto the thermomixer and incubate at  21 °C for  00:03:00 at 1400 rpm.



3m

- 14 Place tubes onto magnetic rack and wait until beads pellet (~  00:01:00 ).

1m




- 15 While leaving tubes on the magnet, carefully remove supernatant without disturbing the beads using a P1000.

Note





Pipette gently, be careful not to shear DNA

- 16 For the second wash, remove tube from magnetic rack and add  700  $\mu$ L Buffer MW1 to the tube, flick the tube until bead pellet is unstuck from tube wall.




Note

MW1 Wash 2

- 17 Place the tube onto the thermomixer and incubate at  21 °C for  00:03:00 at 1400 rpm.

3m



- 18 Place tubes onto magnetic rack and wait until beads pellet (~  00:01:00 ).


1m

- 19 While leaving tubes on the magnet, carefully remove supernatant without disturbing the beads using a P1000.



Note



Pipette gently, be careful not to shear DNA


- 20 Remove tube from magnetic rack and add  700  $\mu$ L Buffer PE to the tube, flick the tube until bead pellet is unstuck from tube wall.

PE Wash 1

- 21 Place the tube onto the thermomixer and incubate at  21 °C for  00:03:00 at 1400 rpm..

3m



- 22 Place tubes onto magnetic rack and wait until beads pellet (~  00:01:00 ).


1m

- 23 While leaving tubes on the magnet, carefully remove supernatant without disturbing the beads using a P1000.



#### Note




Pipette gently, be careful not to shear DNA

- 24** For the second wash, remove tube from magnetic rack and add  700  $\mu$ L Buffer PE to the tube, flick the tube until bead pellet is unstuck from tube wall.

PE Wash 2

- 25** Place the tube onto the thermomixer and incubate at  21 °C for  00:03:00 at 1400 rpm. 3m



- 26** Place tubes onto magnetic rack and wait until beads pellet (~  00:01:00 ). 1m

- 27** While leaving tubes on the magnet, carefully remove supernatant without disturbing the beads using a P1000. Then use a P200 to remove any traces of Buffer PE.

#### Note



Pipette gently, be careful not to shear DNA

28



While keeping the tube on the magnetic rack, slowly pipette  $700\ \mu\text{L}$  distilled water by pipetting down the side of the tube opposite the pellet, without disturbing any of the beads/pellet. Leave for 00:00:30 - 00:01:00 and then remove the supernatant without disturbing any of the beads/pellet.

#### Note



Pipette gently, be careful not to shear DNA

29

Repeat step 28.

#### Note


If necessary, use a P20 to remove any remaining distilled water.



Pipette gently, be careful not to shear DNA

30



Remove the tube from the magnetic rack and add  50  $\mu$ L of nuclease free water. Resuspend the beads by slowly pipetting or stirring with the pipette tip.



#### Note

Be as gentle as possible while ensuring that pellet is resuspended.



Pipette gently, be careful not to shear DNA

31

Place the tube onto the thermomixer and incubate at  21  $^{\circ}$ C for  00:03:00 at 1400 rpm.

3m



**32** Pellet beads on magnet and slowly pipette DNA eluate into a new lo-bind tube.

#### Note



Pipette gently, be careful not to shear DNA

Save the tube with magnetic beads in case of incomplete DNA elution, so you can repeat.