



# General Fungal DNA Extraction V.2

Version 1 is forked from [Kasson Lab DNA Extraction](#)

Matthew T  
Angie Macias<sup>1</sup>, Kasson<sup>1</sup>,  
<sup>1</sup>West Virginia University

Brian  
Lovett<sup>1</sup>

 Brian Lovett  
USDA-ARS

VERSION 2

AUG 24, 2023


## ABSTRACT

This is a routine protocol for extracting DNA from various fungi. This extraction method is suitable for follow-up molecular work such as PCR amplification.


## MATERIALS

Sterile micropestles, isopropyl alcohol, ethyl alcohol, cell lysis buffer, protein precipitation buffer, elution buffer, metal scraper.

## PROTOCOL MATERIALS

 Cell Lysis Solution, 1000ml (for Wizard Genomic) Promega Catalog #A7933


In [2 steps](#)

 Nuclei Lysis Solution, 1000ml Promega Catalog #A7943


Step 3

 Ethyl Alcohol Merck MilliporeSigma (Sigma-Aldrich) Catalog #E7023


In [2 steps](#)

 Protein Precipitation Solution 350ml Promega Catalog #A7953

Step 13

 Elution buffer pH 8.0 (250 mL) Alfa Aesar Catalog #J61558

In [3 steps](#)

 isopropyl alcohol Merck MilliporeSigma (Sigma-Aldrich) Catalog #W292907

In [2 steps](#)

OPEN  ACCESS



**Protocol Citation:** Angie Macias, Matthew T Kasson, Brian Lovett 2023. General Fungal DNA Extraction. [protocols.io](https://protocols.io/view/general-fungal-dna-extraction-cy5gxy3w) <https://protocols.io/view/general-fungal-dna-extraction-cy5gxy3w> Version created by [Brian Lovett](#)

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Aug 24, 2023

**Last Modified:** Aug 24, 2023

## Before you begin

- 1 Turn on hot water bath, set to  $65^{\circ}\text{C}$ .
- 2 Pull two Eppendorf  $1.5\text{ mL}$  centrifuge tubes per sample.

2.1 Label both sets of tubes with (short) sample names.



2.2 Label one tube set for each sample with an "I" for



isopropyl alcohol Sigma  
Aldrich Catalog #W292907



Sketch of "I"-labeled tubes (drawing from Angie Macias).

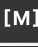
- 3 Add  $200\text{ }\mu\text{L}$  of  
Cell Lysis Solution, 1000ml (for Wizard Genomic) Promega Catalog #A7933 (or  
Nuclei Lysis Solution, 1000ml Promega Catalog #A7943) to **tubes without "I"**.


4 Add  600  $\mu$ L of  isopropyl alcohol Merck MilliporeSigma (Sigma-Aldrich) Catalog #W292907 to tubes labeled with "I".

5 Place tube with  Elution buffer pH 8.0 (250 mL) Sigma Aldrich Catalog #J61558 into  65  $^{\circ}$ C water bath.


## Extraction Protocol

1h 10m 3s

6 Sterilize some metal scrapers with flame and  95 % (v/v)

 Ethyl Alcohol Sigma Aldrich Catalog #E7023




7 Add 1/2 pea-sized amount of fungal tissue (young hyphae) to each tube containing



 Cell Lysis Solution, 1000ml (for Wizard Genomic) Sigma Aldrich Catalog #A7933




7.1 Flame-sterilize and cool scrapers between samples.


7.2 Alternatively, pellet a pea-sized amount of mycelium grown in liquid culture and transfer to each tube.


8 Macerate each sample with a new, sterile micropestle until tissue is homogenous.



9 Add  400  $\mu$ L of  Cell Lysis Solution, 1000ml (for Wizard Genomic) Promega Catalog #A7933 (to  600  $\mu$ L total volume added).



10 Add tubes to a floating rack to allow samples to incubate directly in  65 °C water bath for  00:30:00 . 30m

11 Remove samples and vortex for  00:00:03 before returning to  65 °C water bath for  00:30:00 . 30m 3s

11.1 Place a sufficient aliquot of  Elution buffer pH 8.0 (250 mL) Sigma Aldrich Catalog #J61558 in water bath to warm for Step 21.


12 Remove samples and allow them to cool on the bench for  00:05:00 . 5m

13 Add  200 µL of  Protein Precipitation Solution 350ml Sigma Aldrich Catalog #A7953 to each tube and vortex for 10 seconds.

14 Centrifuge samples for  00:03:00 at  14.000 rpm . 3m

#### Note

Proteins will form a large pellet: unload samples carefully into rack.

15 Using a P1000 micropipette, transfer supernatant to each tube containing  isopropyl alcohol Sigma Aldrich Catalog #W292907 and gently mix by inversion.

#### Note




It's better to leave some liquid than to carry bits of the protein pellet into the next step.

16 Centrifuge for  00:01:00 at  14.000 rpm . 1m

17 Carefully pour off the supernatant into waste container.

**Note**

Be careful to not lose your white DNA pellet!

18 Add  600  $\mu\text{L}$  of  70 % (v/v)  Ethyl Alcohol Sigma  
Aldrich Catalog #E7023 to each  
tube and mix gently by inversion.



19 Centrifuge for  00:01:00 at  14.000 rpm . 1m


20 Repeat Step 16.

21 Open and invert tubes onto a clean paper towel.

**Note**

A tube rack can be placed on the tube lids to secure inverted tubes onto the paper towel.

22 Add  100  $\mu\text{L}$  of warmed  
 Elution buffer pH 8.0 (250 mL) Sigma  
Aldrich Catalog #J61558 to each tube.

23 Store fully-labeled tubes in a box (not a tube rack) in the  -20  $^{\circ}\text{C}$  freezer.

