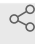




Oct 28, 2022

# 🌐 Calcium imaging in astrocytes

[gustavo.parfitt](#)<sup>1</sup><sup>1</sup>Icahn School of Medicine1 *Works for me* Share

This protocol is published without a DOI.

[Ahfeldt Lab](#)[gustavo.parfitt](#)

## ABSTRACT

Calcium imaging in astrocytes

## PROTOCOL CITATION

[gustavo.parfitt](#) 2022. Calcium imaging in astrocytes. **protocols.io**  
<https://protocols.io/view/calcium-imaging-in-astrocytes-ciheub3e>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED


Oct 27, 2022

## LAST MODIFIED

Oct 28, 2022

## PROTOCOL INTEGER ID

71942

- 1 Coat a 10 cm plate with 0.1% gelatin for  **00:20:00** 20m
- 2 From a 80% confluent 10 cm plate of astrocytes obtained from protocol <https://www.protocols.io/view/astrocyte-extraction-from-brain-organoids-261ge364wl47/v2>

- 3 Add 3 mL of [TrypLE™ Select Enzyme \(1X\), no phenol red Thermo](#)  
[Fisher Catalog #12563029](#) for  
🕒 00:05:00 5m
- 4 🌀 300 rcf, 25°C, 00:03:00 3m
- 5 Add pTALV-fUBIGW-GCAMP8SIREs-PUR0 lentivirus 1/100 in  
[Astrocyte](#)  
[Medium ScienCell Catalog # #1801](#) to the pellet and resuspend the cells
- 6 Incubate 🕒 00:05:00 5m
- 7 Add to a gelatin coated 10 cm plate with a 9 ml  
[Astrocyte](#)  
[Medium ScienCell Catalog # #1801](#) warmed media.
- 8 🕒 48:00:00 observe fluorescence levels to confirm infection 2d
- 9 Wait 🕒 168:00:00 for the GCAMPs expression to reach stable levels 1w
- 10 Coat a 2 cm plate with 0.1% gelatin for 🕒 00:20:00 20m
- 11 Plate 50k astrocytes of infected astrocytes
- 12 after 🕒 48:00:00 change to mature astrocyte media 2d

13 Before the start of imaging aim for 80% confluency.

14 Place the plate on a

Nikon Eclipse TE2000-U microscope	
Microscope	
Nikon	N/A

plate holder.

15 Continuously perfused with ACSF with the following composition (in mM): NaCl 125, KCl 5, D-Glucose 10, HEPES-Na 10, CaCl<sub>2</sub> 3.1, MgCl<sub>2</sub> 1.3. using a ValveBank8 II controller.

16 Using a 10X objective. GCaMP8s was excited using a 480 nm (Mic-LED-480A, Prizmatix), an HQ480/40x excitation filter, a Q505LP dichroic mirror, and an HQ535/50m emission filter (Semrock).

17 Sample at a rate of 4.7 fps with a frame exposure of 200 ms at 160x120 pixels (4x4 binning).

18 ROI segmentation of GCaMP8s, raw fluorescence extraction, and background correction can be performed with Nikon Elements software.